WORLD LEADER FOR FOOT PROFESSIONALS since 1975



INTERNATIONAL OFFER

podiatech.com

The Sidas company







Experience & loyalty since 1975

Since it was established in 1975, the Sidas group has been working with foot professionals and helping to improve the techniques used throughout the industry to ensure the advancement of the profession.

From analysis to training, Sidas equips foot professionals with analysis and measurement tools, materials and concepts, techniques and equipment, and even offers training in material applications and the making of plantar orthotics.



happy planet. 🛈

AUSTRALIA

happy p

Innovative, certified solutions that guarantee quality and reliability

The result of years of international experience in the medical and sports industries, the products developed by our R&D department are tailored to your technology and the needs of your patients.

Through its product traceability system and control procedures, Sidas guarantees quality and reliability. All Sidas solutions comply with European legislation and bear the "CE" mark.



AUSTRIA BELGIUM BRAZIL **CANADA CZECH REPUBLIC** DENMARK FINLAND FRANCE GERMANY INDIA ITALY ISRAEL JAPAN MOROCCO NEW ZEALAND **NETHERLANDS** POLAND RUSSIA SOUTH KOREA **SPAIN** SWEDEN SWITZERLAND THAILAND THE UNITED STATES TUNISIA **UK & IRELAND**

2 www.podiatech.com



Eco-designed headquarter

In 2009, the Sidas group moved to new headquarters in Voiron, Southeast France. Built with environmentally friendly materials according to a simple design, the Sidas offices reflect the group's commitment to the environment.

Maintaining firm ties in the country where it all started, with a production site located close to its headquarters, Sidas also has several production sites around Europe.





Podiatech Partners



UFSP



de la Podologie



INSEP medical







Association pour le Développement

EFOM school - Paris



SERVICES Training Product services: materials and modules OPCT® Factory service	Pages 8 - 9 Pages 10 - 11 Pages 12 - 15	
ANALYSE Platforms: Win-Track / Presscam / Presscam WIFI / Pressmat INTELISOLE®-EASY onboard technology LED Podoscope® – Posturoscope Goniometers - Accessories	Pages 17 - 21 Page 19 Page 20 Page 21	
EQUIPMENT History of PODIATECH Heat-sealing / thermoforming equipment Podiatech Premium station, Podium, Plateforme V-LAB XL [™] , Podolab [®] 2, V-LAB, Mobilab [®] 2 Lab options	Page 22 Pages 23 - 27	VACUUM BAGS WITHOUT HOUSINGS
 Foam base, Vacuum bags, Printlab 2, Sock kit V-vac, V-control, workshop accessories, Cooling system, Cable Positiv'lab, Press-moulding plate, PU positive casts, Oven Workshop 	Page 28 Page 29 Page 30	
 Hoods, Glues, Solvents and Workshop accessories Grinders - Grinder accessories 	Page 31 Pages 32 - 35	
 RAVV MATERIALS Eco-responsibility charter and eco range Top layers Thin materials Decorative materials Antiperspirant materials Antiperspirant materials Durable materials Microperforated materials Bactericidal/fungicidal materials Special use materials Multicoloured materials Top layers in rolls Resin and foam bases Foam bases Comfort and correction elements Technical materials and accessories Combination sheets 	Pages 38 - 39 Pages 38 - 51 Pages 40 - 41 Pages 42 - 43 Page 44 Pages 45 - 46 Pages 46 - 48 Pages 48 - 49 Page 50 Pages 50 - 51 Page 51 Pages 52 - 61 Pages 52 - 61 Pages 52 - 53 Pages 54 - 55 Pages 56 - 57 Pages 58 - 60 Page 61	
ORTHOTICS & ELEMENTS Modulo concept Intelligent top cover OPCT® - OPCT® concept - OPCT® City - OPCT® Multi, Multi sport, Specific Sport - OPCT® Confort, Patho-logic™ - OPCT® Kids - OPCT® specific biomechanichs - OPCT® OC Foam base solutions Resin base modules, latex sheet & PEVA adhesive-coated bases Elements 3Feet® & OTC - 3Feet® Patho-Logic® - 3Feet® Slim	Pages 62 - 63 Pages 64 Page 65 Pages 66 - 68 Pages 68 & 69 - 71 & 72 - 73 Pages 74 - 77 Pages 78 - 79 Page 79 Pages 80 - 83 Page 84 Pages 85 - 86 Pages 87 - 93 Pages 94 - 95 Page 94 Page 95	

3Feet[®] & OTC

FOOT & SHOE PRODUCTS

Foot and shoes accessories Resins, mould release agents and PU accessories Prestoform Ultimate accessories and consumables Materials for fitting and reproduction Shoe manufacturing materials Toe-ends

Pages **96-99** Pages 100 - 101 Pages **102 - 103** Page **104** Page **105** Page **105**



CHAIRS AND SEATS

Patient chairs Options and accessories Chair dimensions Operator's seats Pages **108 - 119** Page **120** Page **121** Pages **122 - 123**



TREATMENT UNITS

Treatment units Integrated unit Kart UK Nano Dimensions of treatment units Pages **126 - 141** Pages **142 - 143** Pages **144 - 145** Pages **146 - 147**



PERIPHERALS

Rotary instruments: micromotors - handpieces - turbines Lighting and accessories Compressors - Aspirators Pages **152 - 155** Pages **156 - 157** Pages **158 - 159**



STERILISATION

Autoclave Tests - printer - ion exchanger - water distiller Soaking tanks, heat sealer, sleeves, bags

MOBILE EQUIPMENT

Micromotors: portable - table / cases / accessories

Page **162** Pages **162 - 163** Pages **164 - 165**

Pages 166 - 167

Pages 170-171

Pages 172 - 175

Pages **176 - 178** Page **179**

Pages 180 - 182

Pages 184 - 185

Pages 188 - 190

Pages 192 - 193

Pages 198 - 199

Page 183

Page 191

Page 196

Page 197

CONSUMABLES

Gouge blades, scalpel blades, handles Nail clippers - scissors Handpiece burs Turbine burs Bandages/band-aids, protection, rasps and accessories Protection and hygiene / Masks, gloves, paper liners and social distancing kit Hygiene and disinfectants

SMALL APPLIANCE

Orthonyxia - Onychoplasty Night splints for hallux valgus - DOKI TAPE Orthoplasty - Ultra-fine mesh orthoplasty - injection and repairs

TREATMENTS


Warts - hyperhidrosis Onychomycosis - hyperhidrosis QENOA products

COMMUNICATION AND INFORMATION TOOLS

Marketing tools

Page 200



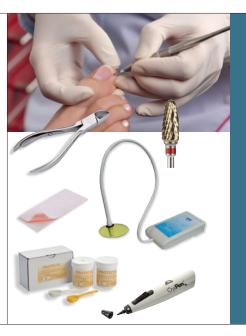



TRAINING, WORKSHOPS AND DEMOS

PODIATECH SIDAS MEDICAL

PODIATECH

SIDAS (1)

 \bigcirc


PRODUCTS: - MATERIALS: combinations, adhesive

coatings, etc. - MODULES: pre-cutting, pre-shaping, etc



SERVICES

COME TO OUR SHOWROOMS TO MEET OUR SALES

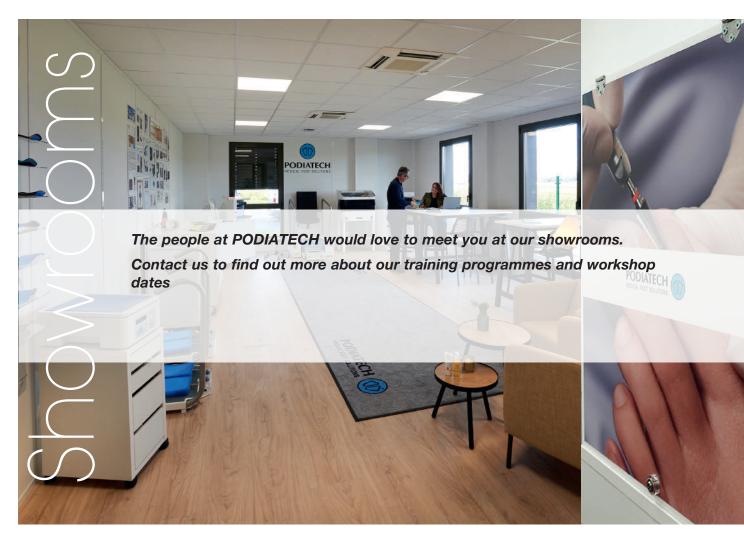





TRAINING Individual and group training sessions available: contact us






TECHNICIANS AND TRAINERS (by appointment)



Available at PODIATECH:


- Practical workshops and demonstrations by appointment for individuals and groups.
- Personalised training on thermoforming techniques and choosing the right materials for you.
- Training sessions with specialists on a variety of topics: running, skiing, analysis tools, clinical examinations, etc.

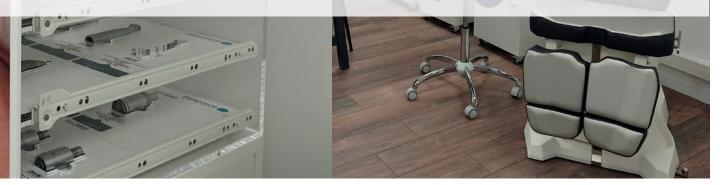






PODIATECH EXTRANET www.extranetsidas.com



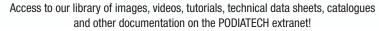

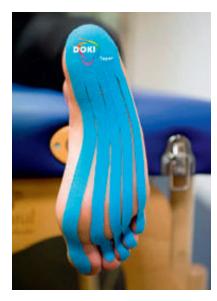

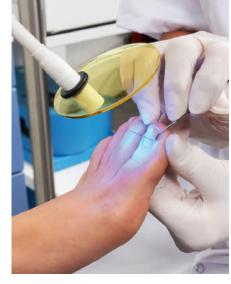

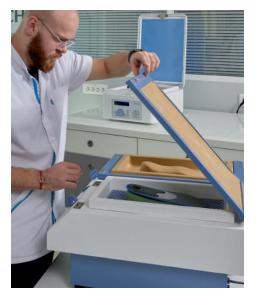



Materials

services +

Custom cut-outs using your shapes or our shapes of elements, top layers and bases in any material or thickness required. On request and subject to technical validation, the sheets can be supplied with an adhesive coating, pre-glued or in combination. Test samples can be sent on request.



Foam and resin base cut-outs



Pre-gluing sheets



Materials in sheets



Cutting out elements



Sticking together sheets



Materials in rolls



Applying adhesive coating to sheets



Combining sheets

services +


Shaping OPCT® S and SM



OPCT® S+

All parts of the OPCT[®] module are **assembled and heat-shaped** (base/middle layer/top layer).



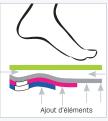

OPCT® SM+

The top layer is separate from the rest (base + middle layer), which is **heat-shaped**, making it possible to include supports and comfort elements right next to the foot.



Option to order certain* OPCT[®] products heatshaped with a glued top layer, and to order products separately <u>with no minimum quantity</u>. From sizes 35 to 45. Delivery within 3 weeks. Contact us.





Custom cut-outs

services

CUSTOM CUT-OUTS OF MODULES, BASES, ELEMENTS, SHEETS AND STRIPS

PODIATECH can create special OPCT® modules and elements for practitioners using its high-precision cutting equipment. This facility making use of the very latest technology is managed by the Methods Department at the PODIATECH production site near Grenoble, France.

Cut-outs from our shapes

All the materials available in sheet form from PODIATECH can be supplied in the form of pre-cut modules based on our patterns.

Cut-outs from your shapes

We can create any form of special shape you need. Contact us for details of minimum volumes and the process we use (technical documentation).


Customising your OPCT®



Choose the patterns and materials you want to create your special OPCT®s. We can supply them to you, subject to a minimum order of 30 pairs (mixed sizes possible). We can supply the models pre-moulded on a positive cast with met bar. For orders of 50 or more pairs, we can personalise your OPCT® modules with your four-colour logo.

Manufacture and delivery within 2 weeks guaranteed. For all your special requirements, contact us.




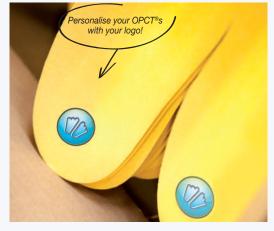







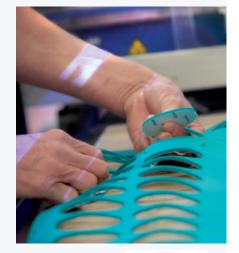

Customising your bases

We can supply Men's and Ladies' archless bases:

- in all the resins of the PODIATECH range (Podiaflex[®],
- Podiaflux[®], Podiaflex[®] on fabric 0.8 mm, 1.3 mm or 1.7 mm); - in all thicknesses: 0.8 mm, 1.2 mm, 1.9 mm (for Podiaflex[®])
- resin) and 1 mm or 1.2 mm (for Podiaflux® resin);
- for all shoe sizes from 34 to 50;
- with the adhesive-coated top layer of your choice* (microfibres, woven textiles, nonwoven textiles);
- posture-stimulation elements cut from various resins, of different thicknesses, for use in postural examinations (pack of testers) and for making posture insoles.

Customising your elements*

We can cut out your own elements to your required shapes from an accurate technical drawing on graph paper.

Your sales technician and trainer can provide details of the special procedure and minimum order number.

*Materials in the PODIATECH materials range.







PODIATECH app



Create 3D scans of your patients' feet.



Use the Structure Sensor and an iPad to scan the vacuum bag imprints



or use the Structure Sensor and an iPad to scan your patient's feet (off load)



or use the Structure Sensor and an iPad to scan the imprint box prints



or use a flatbed scanner.

Use the app to digitally design and personalise your OPCT®.



- Enter your patient information (patient data and clinical examination data).

0	tion that the most of	-0
ñ		
· · · · ·		
0		
A		
A		
·		
ē		
	B	
	·····	
	e •	
0	2	
	2.	

- Choose the d'OPCT[®] model you want from our full range.

(1992) (1			
	land of a		
	100000	0	
	0+ = =	0 = =	
	0	0	
	() = () = () =		

 Indicate which correction elements you want to add to your OPCT[®]. (Advanced and Elite levels only).

		8
		•
U		

- Position and draw the corrective elements wherever you like. (Elite level).
- Simple, intuitive online software.
- Available for iPad and PC.
- Free to download (subject to conditions).
- Compatible with most existing scanning technology.
- Several personalisation levels available.

THE DIGITAL SOLUTION FOR YOUR THERMOFORMED CUSTOM FOOT ORTHOTICS ----



SIDAS factory

We will make your $\mathsf{OPCT}^{\circledast}$ to your specifications at our SIDAS Factory.



- Machining a positive cast* from a scan validated in the Podiatech app.



- Thermoforming the selected $\mathsf{OPCT}^{\circledast}$ on the positive cast*.



- Adding the required correction elements.



Your OPCT® are made by a foot-health professional.

from a recyclable material.

Watch our OPCT Factory tutorials on our YouTube channel: PODIATECH medical foot solutions



Delivery

We send your OPCT® directly to your clinic.





- Option to order the positive cast* used in the thermoforming process.



You receive your custom foot orthotics moulded, made to measure and with the required corrective elements.



You can modify and adapt the OPCT® received.



USE THE STRUCTURE SENSOR AND AN IPAD TO SCAN THE VACUUM BAG IMPRINTS



OPCT® Factory: what our customers say



Used by

Groupama · FDJ

By Marc Retali, specialist in sports podiatry, and Thomas Arnaud, foot health professional.

For me, OPCT FACTORY provides all the advantages of digital technology with none of the disadvantages!"



Additional business and time optimisation

« Beyond the limited investment in machinery required, one of the main advantages is optimising my square meterage and time. No longer having to keep a stock also means we save time, money and space! » Thomas ARNAUD

Reproductibilité et base de données

« This service enables me to make my plantar orthotics using a fully featured app and means I can outsource production. » Marc RETALI



« It has also enabled me to optimise and improve my services because it means I am now able to reproduce identical orthotics at a later date, making repeat orders child's play! »

Thomas ARNAUD

Fully featured app and patient satisfaction

« The plantar orthotics are made at the PODIATECH plant by professionals. So I can focus on clinical examination, saving athletes precious recovery time.

I keep a database of scans, which means I can create instant reproductions and send the orthotics all over the world and at any time of the season. $\!$

Marc RETALI

« My long-standing patients say they can't tell any difference either. [...] And my patients like the fun, modern experience of using a tablet for their orthotics. » Thomas ARNAUD



Limited investment and stock

« [...] a module of excellence is nothing without faultless moulding, which imposes particular logistical constraints: the constant transporting of equipment to make imprints and plantar orthotics. » Marc RETALI

« With such a broad range of OPCT[®]s available, I can obtain an orthotic suitable for any type of activity without having to manage a vast stock. » Thomas ARNAUD





The healthier option

The advantages of OPCT®

« I've been a fan of the products in the Podiatech range for a few years now, and I know I can count on the quality of my finished products. » Thomas ARNAUD

« PODIATECH offers a wide range of technical materials and pre-designed plantar orthotics tailored to all fitting rooms and activities. [...] Luckily for me, the brand works hard on continually adapting its products. »

Marc RETALI



Pricing structure for an OPCT[®] made using OPCT[®] Factory



Choose your OPCT[®] from our full range. Price: refer to the currently valid PODIATECH catalogue (the reference price is the price of the SM version of the OPCT[®].



Different personalisation levels:

PRO

- Incorporates the scan,
- Select the OPCT[®] you want,
- Manufactured by the SIDAS Factory.

ADVANCED

- Incorporates the scan,
- Select the OPCT[®] you want,
 Select the corrective elements you
- want, • Manufactured by the SIDAS
- Manufactured by the SIDAS Factory.

ELITE

- Incorporates the scan,
- Select the OPCT[®] you want,
- You choose the number, design and position of corrective elements,
- Manufactured by the SIDAS Factory.



· Glued top layers

- Fully finished
- (Grinding of corrective elements + template cut-out) FREE OF CHARGE IN 2022
- Addition of corrective element(s) (All corrections to the structure of the positive cast are free of charge)
- · Return of positive cast

Subject to a minimum order number, shipment fees as stated in terms and conditions. Contact us for a copy of the applicable version of our special terms and conditions of sale.

• OPCT® Factory template PDVGABFACTORY

The 3D scanner solution recommended by PODIATECH




Protective case

Sensor structure

Sensor structure stand

iPad tablet

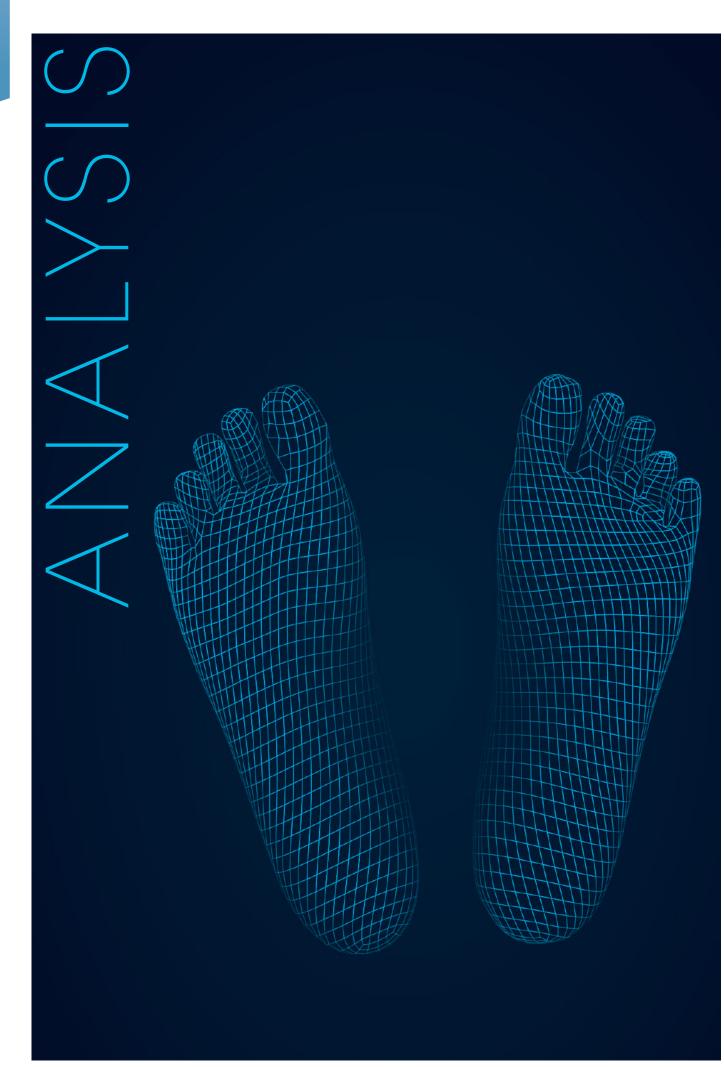




WIN-TRACK


Win-Track: a foot pressure sensitive walkway dedicated to static, posturographic and gait analysis with total freedom in stance and motion acquisition.

WIN-TRACK - AMWINTRACK


Dimensions	1610x652x30mm
(Length x Width x Height)	
Thickness	9mm
Sensing area	1500x500mm
Weight	11.8Kgs
Sensor technology	Resistive HD
Sensor's size	0.78x0.78cm
Number of sensors	12,288
Pressure range	From 0.4 to 100N
Temperature range	From 0°C to +85°C
PC/Platform interface	USB 2.0
Sampling rate	200Hz
Required operating system	Windows XP, Vista, Seven

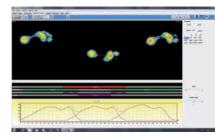

Dynamic analysis

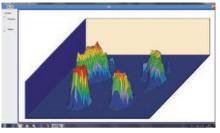

Static analysis

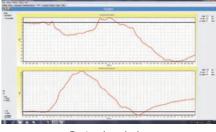

Postural analysis



WIN-TRACK MINI

WIN-TRACK Mini is a compact walkway designed for dynamic podometric analysis that can also be used for static posturographic data acquisition with the subject in any stance.

WIN-TRACK MINI - AMWINTRACKM

Dimensions (Length x Width x Height)	1115 x 652 x 30 mm
Thickness	9mm
Sensing area	1000x500mm
Weight	8Kgs
Sensor technology	Resistive HD
Number of sensors	8192
Pressure range	From 0.4 to 100N
Temperature range	From 0°C to +85°C
PC/Platform interface	USB
Required operating system	Windows 7, 8, 10





Video analysis with measurement of angles



You Tube

Watch our Presscam video



Static analysis



Dynamic analysis



Download the free Podocapture app to give you greater mobility in your podometric analyses! (For Presscam WiFi only.)

PRESSCAM

Presscam is a hybrid pressure-posturology platform. It is powered via a USB cable connected to a computer.



ADVANTAGES

- Clear, higher-resolution images and videos: 100 colour levels (compared to 20 in previous versions).
- Video connection with option to measure angles and lengths.

PRESSCAM PLATFORM - AMPRESSCAM5-0

- Dimensions (width/length): 550 x 661 mm
- Active surface area: 400 x 400 mm
- Base: 9 mm thickness
- Weight: 3 kg

UPDATES AND MAINTENANCE CONTRACT, 1 YEAR AMPRESSUP5-0

MAINTENANCE CONTRACT, 1 YEAR P52MAINTENAN

PRESSCAM CARRY CASE AMHOUSSE

PRESSCAM WIFI



New design featuring wireless connectivity

ADVANTAGES

- With wireless connectivity: can be used any time, anywhere
- Thin base: only 8 mm
- Mobility.

Technical information

The Presscam Wifi has the same technical characteristics and uses the same software as the Presscam. - Presscam Wifi 150 Hz / Presscam 100 Hz

- The Presscam WIFI package comprises the following:
- 1 Presscam Wifi platform
- 2 lithium-ion batteries
- 2 micro-USB cables for charging the batteries
- 1 x 220V / 5V mains adapter for recharging
- the batteries from the mains power supply
- 1 x 5-metre back-up USB cable
- 1 x Presscam software package on USB stick.

PRESSCAM WIFI PLATFORM AMPRESSCAMWIF





TRAINING Individual and group training sessions available: contact us



The new analysis tool that communicates with patients & medical professionals.

PRESSMATTM

Instant, accurate analyses for STATIC, STATIC POSTURAL and DYNAMIC assessments.

- 1 intuitive dedicated app
- 1 Bluetooth connection to your Android tablet.
- Size of sensors: 2 mm x 2 mm
- Number of sensors: $64 \times 64 = 4096$
- Type of sensor: resistive
- Frequency: 30 Hz (Bluetooth) 100 Hz (USB)
- Active surface area: 445 mm x 445 mm
- Device footprint: 600 mm x 630 mm
- Weight: 1.7 kg
- Measurement from 0 to 40 N/cm2
- 5-hour battery life in operation

ADVANTAGES

- Flexible, light, easy to transport and washable
- Option to add notes to each analysis, which are automatically appended to the patient report
- Non-slip pads

PRESSMAT™ PLATFORM - AMPRESSMAT







Pressmat[™] vidéo





INTELISOLE EASY

INTELISOLE Easy is used for investigating foot pressure distribution. It is your perfect partner for facilitating diagnosis, correcting and improving the gait of all kinds of patients, optimising athletic performance and measuring the effect of insoles in real conditions (at your practice and in the real world).

- Compact.
- 4 pairs of insoles and 8 separate sensors positioned in relevant zones (covering shoe sizes 35 to 46)
- CE-certified medical device
- Maximum pressure measured: 8 kg/cm2
- Bluetooth communication with PC and Android smartphone

Light, efficient control unit

- Acquisition at 100 Hz
- Automatic synchronisation of the 2 control units
- Fast USB charging (1 hr 30) - Simple clip fastening to laces or strip

High-precision pressure sensors

- Insole with 8 ultra-hardwearing calibrated sensors
- Guaranteed for 1 million activations (validated on test bench)

INTELISOLE EASY - AMWINSHOE

MACHINES



PLEXIGAS[®] **PODOSCOPE WITH LEDS**

- Automatic or manual illumination

- LED liahts
- Choice of 3 light levels
- Universal remote control
- Height 125 mm
- Capacity 150 kg

PLEXIGLAS® PODOSCOPE WITH LEDS - APODOSIDAS 20

Warranty duration 2 years (use of alcohol-based products prohibited)

POSTURSCOPE

Posture visualisation and analysis tool.

For use with a Podoscope or pressure platform for full analysis.

- Tilting mirror with return mark
- Plumb line
- Red background grid
- 4 independently adjustable feet

POSTUROSCOPE ONLY - AMPOSTURSCOP

Supplied in kit form

GRID BACKDROP

Backdrop featuring a grid for postural analyses.

Two eyelets for easy installation.

GRID BACKDROP - P80POSTPOSTURO



EXAMINATION COUCHES



 (ϵ)

FOLDING EXAMINATION COUCH IN BLACK OR WHITE - PFP021G

- Height: 80 cm
- Width: 56 cm
- Load capacity: 110 kg dynamic load
- Easy to transport
- Warranty: 2 years

ECOMAX EXAMINATION COUCH - PFP040

- Height: 80 cm
- Width: 60 cm
- Adjustable back
- Upholstered with 5cm-thick foam
- Load capacity 150 kg dynamic load
- Integrated paper roll dispenser
- Warranty: 5 years

NON-SLIP STEP, ABS PLASTIC (BLACK) - PFP044 - Warranty: 2 years

ELECTRIC PHYSIO EXAMINATION COUCH - PFP047

- Adjustable height
- Warranty: 2 years

New colours available: contact us



Colours available for ECOMAX

Heiaht 2.25 m Mirror amplitude +/- 10 cm



GONIOMETERS





Bunions STAINLESS-STEEL GONIOMETER AMGONIOHVAL Ideal for use with bunions



PROTRACTOR - AMGONIOANGLE High-precision semi-circular protractor, diameter 16 cm





ROM GAUGE - AMGONIOARTI Stainless-steel goniometer



Height 13 cm - Width 11.5 cm Depth 1.2 cm





ACCESSORIES

DERMATOGRAPH PENCIL

Red P46D002 Blue P46D003

STAINLESS-STEEL RULE, 30 CM





PLUMB LINE - P46D001

SKELETON FOOT - CMKPVPS



PACK OF 3 ANATOMICAL FEET: FLAT, HOLLOW AND NORMAL CMKPVPA







BABINSK REFLEX HAMMER P65M001



MADDOX KIT WITH LENS + EYE COVER - AMMADDOXKIT

PODOMETER - AMPODOMETRE



MONO FILAMENT -P81MONOFILAM

FEET REVELATOR™ - ABFEETREVE12



FEET MEASURER

ABMES-SPORT-DBLE12

-

P46D004







PODIATECH, the originator of thermoforming and sports podiatry

1984		Thermopodoform:
		The first use of thermoforming in podiatry.
1987		Podiaflex®: The first glue-free thermoformable material for insoles on the market.
1993		Podolab[®]: Launch of a new thermoforming machine.
1997		OPCT®: Innovation of the first module of pre-assembled and pre-cut materials.
2004	UFSP UNION FRANÇAISE POUR LA SANTÉ DU PIED	UFSP: PODIATECH has been actively supporting the French Union for Foot Health (UFSP) since 2004.
2013	CONNAISSANCE & EVOLUTION	Connaissance et Évolution: PODIATECH became a partner of this organisation, which promotes training, knowledge-sharing and development in podiatry.
2017		Astana pro cycling team: PODIATECH began supporting these elite athletes, supplying them with plantar orthotics.
2020	EVEN OPCT®	OPCT® Factory : Launch of a digital design and manufacture service for bespoke thermoformed plantar orthotics.



940 mm

630 mm

1100 mm

Step dimensions: H. 170 mm - D. 600 mm - L. 240 mm

PODIATECH PREMIUM STATION



ERGONOMIC AND COMPACT WORKSTATION

FOR WORKING UNDER PATIENT'S FULL WEIGHT OR HALF WEIGHT, IN FRONT OF OR BEHIND THE PATIENT

CLEVER SLIDE-RAIL SYSTEM

2 STANDARD VACUUM BAGS INSTALLED

MIRROR WITH MARKINGS TO CONFIRM CORRECT POSITION

COMFORTABLE ADJUSTABLE SEAT

COMPLETE PODIATECH PREMIUM STATION WITH STANDARD VACUUM BAGS

OPTIONAL EXTRA: 2 ADDED-VALUE HARD VACUUM BAGS

OPTIONAL EXTRA: ADDED-VALUE PREMIUM SEAT

3-year warranty* (does not include membrane)

BLUE PODIATECH SEAT - PSIE011_BL

PODIATECH PODIUM

PODIUM WITH STEP AND NO EQUIPMENT CMAVAPACKSTAT2_1

- Height 840 mm
- Depth 900 mm
- Height of base 350 mmHeight with bars 1210 mm



PODIATECH PLATFORM

PODIATECH PLATFORM - POPPLATEF2

- Height 940 mm
- Width 630 mm - Depth 600 mm
- Weight 17.6 kg
- Weight 17.0 Kg

600 mm

V-LAB XL™

The fast way to make several heat-sealed orthotics




BUILT-IN FANS FOR INTENSIVE USE

25% SAVING IN COOLING TIME


QUIET, HIGH-PERFORMANCE VACUUM PUMP

SOFT AND HARDWEARING MEMBRANE

BUILT-IN THERMOSTAT AND TIMER

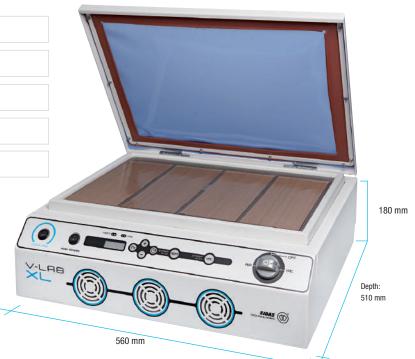

Technique for making insoles with the V-LAB XL™

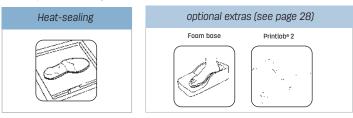

V-LAB XLTM- PGLVLABXL

Includes packaging and shipment (by carrier) 3-year warranty* (does not include membrane)

*Provided use is in accordance with user instructions, machine guarantee excludes damage caused by the use of plaster positive casts.

Delivered with the following accessories:

- 8 Teflon protective sheets
- 1 Large Teflon protective sheet
- 1 Bottle of talc
- 2 OPCT[®] protectors

V-LAB XL MEMBRANE- PZMAMEMB6

V-LAB XL FRAME SEAL- CZMAJOINTMO12

V-LAB XL[™] works best when combined with the PRINTLAB 2






PRINTLAB[®] 2 WITH STANDARD SILICON BAGS PGNPRINTLB2 + POPESV17_24PBL

Includes packaging and shipment (by carrier)

PRINTLAB[®] 2 WITH HARD SILICON BAGS PGNPRINTLB2 + POPESV17_24PBL

Includes packaging and shipment (by carrier)

More information on page 28 of the catalogue.

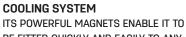

BE FITTED QUICKLY AND EASILY TO ANY PODIATECH MACHINE. More information on page 29 of the catalogue



640 mm

For making orthotics directly on the foot, in a corrected position.

580 mm

24 www.podiatech.com

80 mm

PODOLAB® 2 Design and Innovation



SUPPLIED WITH A PEDAL

TILTING MULTI-FUNCTION TOUCH SCREEN CONTROL PANEL

2-SETTING MEMORY FOR HEATING TIME AND TEMPERATURE


SOFT AND HARDWEARING COLOURED SILICONE MEMBRANE

TWO SEPARATE HEATING ZONES

BUILT-IN CARRY HANDLES



Thermoforming techniques:

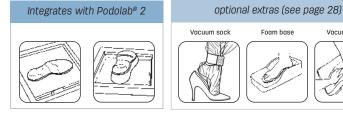

Podolab[®] 2 - PGLPODOLAB2

Includes packaging and shipment (by carrier)

- Delivered with the following accessories: - 1 Pack of 8 Teflon protective sheets
- 1 Large Teflon protective sheet
- 2 OPCT forefoot protectors
- 1 Bottle of talc
- I Bottle of tal
- 1 OK Net cleaning spray
- - 1 Pedal

Vacuum bag

PODOLAB II MEMBRANE - PZMAMEMB5


PODOLAB II FOAM SEAL - CZMAJOINTMO10

3-year warranty* (does not include membrane)

*Provided use is in accordance with user instructions, machine guarantee excludes damage caused by the use of plaster positive casts.



Vacuum socks kit: see page 27

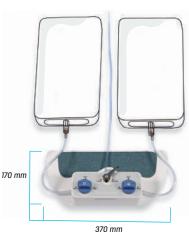

370 mm

V-CONTRÔL



Ergonomic control unit for use with the vacuum bags. Also compatible with Grey vacuum bags.

V-CONTRÔL + 2 COUPLINGS for Grey vacuum bags POPVCONTROL + POPTUBESV08

V-CONTRÔL (sold separately without vacuum bags) POPVCONTROL

More information on page 29 of the catalogue.

V-LAB


HEAT SEALER WITH QUIET, HIGH-PERFORMANCE VACUUM PUMP

SOFT AND HARDWEARING COLOURED SILICONE MEMBRANE

EASY-CHANGE MEMBRANE



Thermoforming technique:

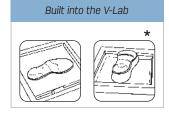


Heating system: Built-in thermostat and timer

Digital temperature and time display Vacuum start Temperature adjustment Cycle start (count-down + heating) Time adjustment

Heating surface:

350 mm x 275 mm

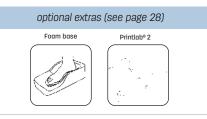



V-LAB - PGLVLABFUS

Includes packaging and shipment (by carrier)

Delivered with the following accessories:

- 1 Large Teflon protective sheet
- 1 OPCT® forefoot protector
- 1 Bottle of talc



V-LAB MEMBRANE - PZMAMEMB7

V-LAB COVER SEAL - CZMAJOINTMO11

V-LAB works best when connected to the V-VAC

V-VAC

Separate control unit. More information on page 29 of the catalogue.

VLAB V-VAC - PGNVVAC





480 mm

MOBILAB® 2 Design and Innovation

2 EXTERNAL VACUUM OUTLETS

SUPPLIED WITH A PEDAL FOR OPERATING THE VACUUM PUMP

MULTI-FUNCTION TOUCH SCREEN CONTROL PANEL

3-SETTING MEMORY FOR HEATING TIME AND TEMPERATURE

Thermoforming technique:



VACUUM SOCKS KIT 1 roll of jersey cloth + 1 pack of 50 pairs of vacuum sock bags + 1 filter attachment + 1 power strap P99KITCHAUSSETTESV



1 PACK OF 50 PAIRS **OF VACUUM SOCK BAGS** P99SACSVID



1 FILTER ATTACHMENT -PPREMBOUT



EQUIPMENT

CE

ands (f)

360 mm

Includes packaging and shipment (by carrier)

cover closed

180 mm

cover open 660 mm

505 mm

3-year warranty* (does not include membrane)

*Provided use is in accordance with user instructions, machine guarantee excludes damage caused by the use of plaster positive

MOBILAB II BLUE MEMBRANE - PZMAMEMB4

MOBILAB® 2 UNIT - PGLMEUBLMOB2

Unit sold without Mobilab® 2 and without accessories

1 POWER STRAP- P99GUETRE

THERMOFORMING TECHNICAL ACCESSORIES FOAM BASES

 \frown For making orthotics directly under the foot in the uncorrected position



- For thermoforming plantar orthotics of any size
- Removable bottom and side inserts
- Two types of inserts can be used (hard or soft) depending on the patient's foot shape
- Integrated positioning markings

Silicon bag full of microballs.

1 HARD 5 KG - POPESV17_50VBL

CONNECTING TUBE - POPTUBESV08

WEDGE - POPCALESV17

For setting the foot in a corrected position 1 STANDARD 2.4 KG - POPESV17_24PBL

PAIR OF BLACK FOAM BASES - POPTMN20_NR_PACK

PAIR OF PROTECTIVE MEMBRANES - POPTMBKITMEN

For setting the footprint without causing the bag to sag.

2-year warranty (does not include membrane)

CE

VACUUM BAGS





You Tube

View the video of our vacuum bag moulding techniques

PRINTLAB[®] 2



For making orthotics directly on the foot, in a corrected position.



PEDAL (option) - CMAVAHD09PED Vacuum-control pedal for hands-free use.

VACUUM SOCKS KIT





Two separate silicon vacuum bags filled with microballs connected to a built-in vacuum pump with filter. It can be used to set each foot imprint separately and obtain better stability during insole moulding.

ADVANTAGES

- Ergonomic
- Portable, ideal for working on the road
- Bags made of soft silicone, a material very resistant to heat and ageing
- Extremely malleable

See dimensions page 24

PRINTLAB® 2 WITH STANDARD SILICON BAGS PGNPRINTLB2 + POPESV17_24PBL

PRINTLAB® 2 WITH HARD SILICON BAGS PGNPRINTLB2 + POPESV17_24PBL

Includes packaging and shipment (by carrier)


Delivered with the following accessories:


- 1 OK Net cleaning spray
- 1 Dismantling spanner
- 1 Bottle of talc

VACUUM SOCKS KIT - P99KITCHAUSSETTESV

1 roll of jersey cloth + 1 pack of 50 pairs of vacuum sock bags

+ 1 filter attachment + 1 power strap

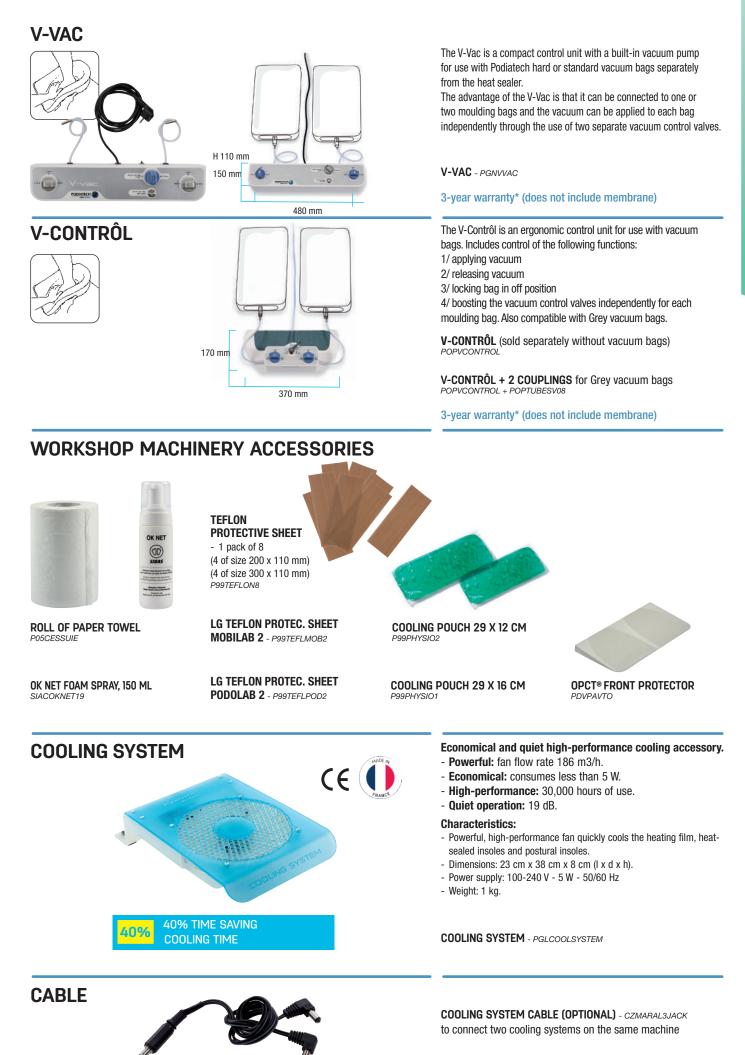

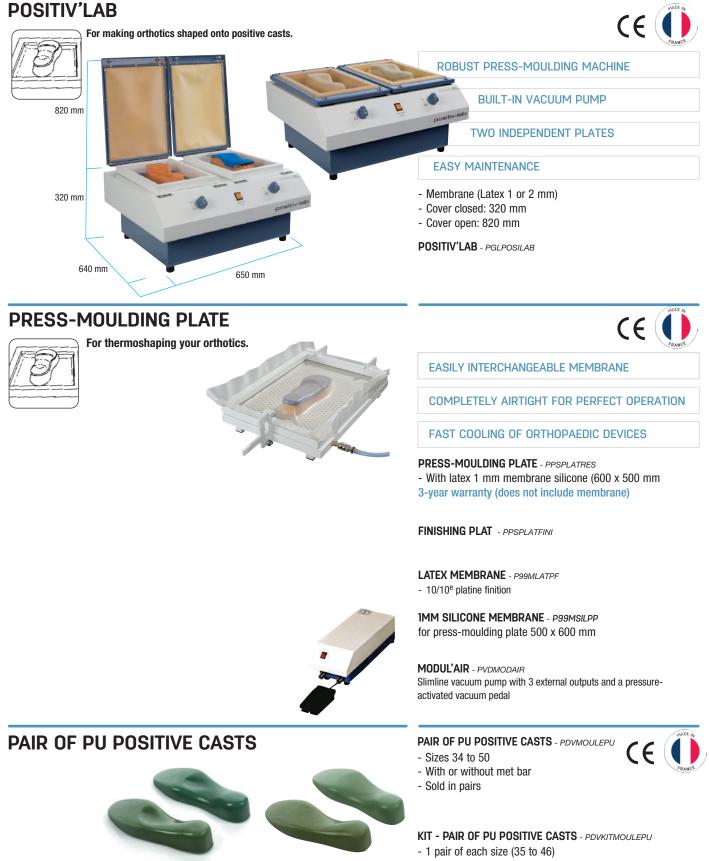

- With or without met bar
- Compact and easy to use.
- CE - Maximum temperature 105° C (variation of + or - 5° C).
- Height: 450 mm
- Width: 400 mm
- Depth: 580 mm
- Weight: 15 kg

OVEN - COUCHFOU2

30 www.podiatech.com

MOBILE FILTERING FUME HOOD

BUILT-IN LIGHT



750 mm

630 mm

Technical characteristics

Power output / 80 W
 Fold-back door:

Weight: 38 kg

- 230 V - 50 Hz

- Variable flow rate / 50 to 450 m3/h

H. 750 mm x D. 630 mm x W. 700 mm

The air expelled inside the workshop is completely safe to breathe.



- Built-in turbine operating at 50 to 450 m3/h.
- For use on work surface.
- Dim. H. 470 mm x D. 500 mm x W. 350 mm
- Weight: 7 kg
- 230 V, 50 Hz, 80 W

MOBILE FILTERING FUME HOOD - PEXPACKHOTTE Supplied with one carbon filter

Includes packaging and shipment (by carrier)

ONE SPARE FILTER (active charcoal) - P99FILCHARAC



- Fume hood suitable for gluing operations. Fold-back door for hand access, 250 mm clearance.
- Easy maintenance: periodic filter change

Using the fume hood

Recycling mode: active-charcoal filtering or dust filter Exterior discharge mode: exterior discharge filter + evacuation hose Ø 125 mm

Fume hoods of other dimensions available on request. Contact us for a quotation

FUME HOOD PACK WITH FOLD-BACK DOOR - PEXPACKHOTTE2 with active-charcoal filter included

ACTIVE-CHARCOAL FILTER - P99FILCHARAC



700 mm

www.podiatech.com 31

G5 V3 GRINDER Evolution





Streamlined grinder unit with built-in aspiration. Quiet operation 68 db (A), very powerful motor.

- 2 grinding stations:
- > Station 1 on the right (when facing the grinder):

manual self-locking chuck to which different grinder accessories can be fitted.

> Station 2 on the left (when facing the grinder):

rubber expanding wheel, Ø 90 mm, 100 mm wide, designed to hold abrasive belts.

Advantages:

- grinding belts can be changed quickly and easily, no tools required.
- generous guard clearance enables operator to grind large sizes easily (removable side door).
- fitted with a professional aspiration turbine which directs dust into the collection drawer.
- easy maintenance: bagless design (collection drawer).
- mounted on self-locking castors.
- generous guard clearance.
- additional side power socket
- emergency stop button for safety.

G5 V3 GRINDER - PPGTOURG5-V3

G5 GRINDER ACCESSORIES KIT - PXACC_TOURET

OPTIONAL G5 12CM HEIGHT EXTENDER - PPGTOURG5REH






4 Dust collection drawer

G4 GRINDER Evolution 2



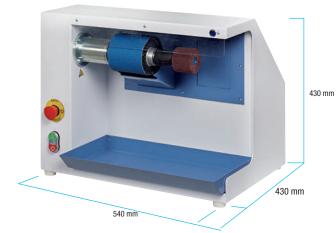

BELT DRIVE

2 ASPIRATION DUCTS

VERY HIGH TORQUE

EXPANDING WHEEL

STANDARD GRINDER ACCESSORIES KIT PXACC_PONSEU

COMPLETE GRINDER ACCESSORIES KIT PXACC_TOURET

Voir descriptifs des packs page 35

VACUUM CLEANER ONLY PPGASPIG4EVOLUTION Flow rate: 850 m³/hour



- Abrasive belts can be changed quickly and easily without tools.

- Vacuum cleaner attachment point at rear of grinder 2-year warranty

ÉVOLUTION 2 G4 GRINDER ONLY - PPGTOURG4EVOLUTION

ÉVOLUTION 2 G4 PACKAGE EVOLUTION 2 G4 GRINDER + ASPIRATOR PPGTOURG4EVOLUTION + PPGASPIG4EVOLUTION

G3 EVO GRINDER



TURBINE VACUUM CLEANER (600W)

EXPANDING WHEEL

- Abrasive belts can be changed quickly and easily without tools.

2-year warranty

G3 GRINDER - PPGTOURG3EVO_230

Not including delivery and installation

STANDARD GRINDER ACCESSORIES KIT - PXACC_PONSEU

COMPLETE GRINDER ACCESSORIES KIT - PXACC_TOURET

See package descriptions on page 35

(6

GRAVITY EVO GRINDER

(1) Height of user:

700 mm

650 mm min.

800 mm max

- < 165 cm: extensions not required
- 165 to 180 cm: 1 set of extensions required
- > 180 cm: 2 sets of extensions required

 $\langle 1 \rangle$



WIDE GRINDING AREA

- Enclosed grinder without vacuum extractor
- Grinding belts can be changed quickly no tool required.

CE

2-year warranty

NEW FEATURES:

530 mm

265 mm

GRAVITY

- Aluminium LED control switch
- Safety shutdown switch
- New shelf (height adjustable on installation)
- Extensions available to adjust grinder height to user's height (1)
- Improved ergonomic design: more space underneath expanding wheel (2)
- Suction inlet grid designed to facilitate cleaning (3)

GRAVITY EVO WITHOUT STAND - PPGTOURGRAVEVO

GRAVITY EVO WITH STAND - PPGTOURGRAVEVO

OPTIONAL KIT: SHELF + 4 LEGS + 8 EXTENSIONS *PZPGTRGRAV600MEUBLE*

STANDARD GRINDER ACCESSORIES KIT PXACC_PONSEU

COMPLETE GRINDER ACCESSORIES KIT PXACC_TOURE

Includes packaging and shipment (by carrier)



(2) Improved ergonomic design: more space underneath expanding wheel



800 mm

(3) Suction inlet grid designed to facilitate cleaning

See package descriptions on page 35





2-year warranty

- Workshop height extender/mount included

COMPLETE GRINDER ACCESSORIES KIT - PXACC_TOURET

Kit contents:

- 1 dust brush
- 3 coarse-grit abrasive belts
- 3 fine-grit abrasive belts
- 1 grinding tip holder
- 3 grinding tips
- 1 hard polishing roller
- 1 soft polishing roller
- 1 grinding roller
- 1 cleaning rubber
- 1 pair of safety goggles
- 1 pack of 2 masks
- 1 apron

GI EVO DEBRIS BUCKET

RIGHT-HAND DEBRIS BUCKET - PPGTOURG1EVO_AD

LEFT-HAND DEBRIS BUCKET - PPGTOURG1EVO_AG



MAGNETIC LED GRINDER LIGHT - 5 W COUCHLEDLIGHT

2 MAGNETIC HALOGEN GRINDER LIGHT - 20 W PPGLAMPAIMAN

TOURET GI LAMP BRACKET - PPGTOURG1EVO_SUPLAMP

STANDARD GRINDER ACCESSORIES KIT - PXACC_PONSEU

Kit contents:

- 3 coarse-grit abrasive belts
- 1 cleaning rubber
- 1 hard polishing roller
- 1 pair of safety goggles
- 1 pack of 2 masks

GRINDER ACCESSORIES



PROTECTIVE MASK -5 MASKS P99MASQUE



GRINDING BELTS Ø 90 mm - largeur 100 mm FINE GRAIN POTABRTSP_90_GF

COARSE GRAIN P07ABRTSP_90_GG



PAIR OF SAFETY GOGGLES



AND SUPPORT - POTSUPCAP

HOOD SHAPE GRINDER X 10 P07CAPUCHON



SANDING WHEEL [Ø 80 MM] P07FEUTRE2



CLEANING BLOCK FOR BELTS P07CREPE21_01

SET OF 2 CLEANING BLOCKS FOR BELTS - P07CREPE21_02



POLISHING WHEELS 180 SOFT POTFEUTRED

POLISHING WHEELS 280 HARD - PO7FEUTREF





MATERIALS IN SHEETS AND ROLLS

Specialised materials, at the forefront of your profession. The fruits of our research, these are pioneering products that will shape the future. They feature innovations and technological advances that are at the forefront of progress in your profession. They are available exclusively in the Podiatechrange.

Tried and tested materials that you can choose with complete confidence.

These materials have been around in modern podiatry for several years. Approved by thousands of users, they have become indisputable standards that you choose for their consistency and quality.

Access materials, accessible quality.

Like access drugs, these are materials that have recognised qualities. They comply with every aspect of the Podiatech specifications but are still available at very keen prices for optimum value for money.

Microfibre

S

ACCESS

A very fine, elastic and supple material. Very hardwearing. Presents no exposed threads during finishing grinding.

Nonwoven

Highly resistant to mechanical stress. Low elasticity.

EVA Very hardwearing and elastic. Heat-sensitive.

Woven

Hardwearing.

PE A comfortable foam. High temperature resistance.

Very good value for money.

Can reveal exposed fibres

during finishing grinding.

PE-EVA Comfortable and hardwearing. Suitable for most uses.



MAXIMUM RECOMMENDED TEMPERATURE for thermoforming the material. Above this temperature, the material may deteriorate (both in appearance and properties).



POSTURE : Material suitable for making postural orthotics.

SHORE A:

Shore is a measurement of the hardness of a foam. This hardness is directly related to comfort. Our technical and comfort materials have a shore A of 10 to 20. Our dynamic foams are 20 to 40, and for our correction and stabilising elements we have to use harder materials with a shore A of 40 to 70.

DENSITY:

Stated as kg/m3. The denser a material, the more compact it is, so density has an effect on the lifespan of a material. Low density is between 40 and 100, medium density is between 100 and 250, and high density is between 250 and 800.

www.podiatech.com 37

PODIATECH MATERIALS Eco-responsibility charter

For transparency, PODIATECH has classified its materials to indicate their level of eco-responsibility.

LOCAL: place where the material is made.



GREEN PRODUCT: analysis of all of the manufacturing processes and components of the material. Provides an overview of the product's environmental credentials (biosourced, recyclable, biodegradable, biocomponents, etc.).



SPOTLIGHT: OUR ECO MATERIALS

With the support of our research and development team, PODIATECH is today committed to developing the cleanest possible range of materials for making orthotics.

TOP LAYERS



REINFORCEMENTS Resin

	CONSTRUCTOR CONST	PODIAFLEX® SOFT 0.8 mm 1.2 mm Grey	- Sheet: 1000 mm x 750 mm - 0,8 mm Grey - 1.2 mm Grey
	S S S S S S S S S S S S S S S S S S S	PODIAFLEX® N 0.9 mm 1.3 mm 1.7 mm Plain	 Sheet: 1000 mm x 750 mm 0.9 mm Nature 1.3 mm Nature 1.7 mm Nature
FOAMS			
VEW	S C C C C C C C C C C C C C C C C C C C	POSTEVA ECO 55 1.5 mm Black	- Sheet: 1000 x 850 mm - 1.5 mm Black
VEN	Posture Posture Posture Posture T. 80°	POSTEVA ECO 70 1 mm Noir 2 mm Black 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0 0	 Sheet: 860 mm x 770 mm 1 mm Black 2 mm Black
CORRECTION ELEMENTS			
poployed for Biomechanics	S S S S S S S S S S S S S S	PODIEVA ECO 50 5 mm - 6 mm Carded Black	 Sheet: 1000 x 850 mm 5 mm Carded Black 6 mm Carded Black
Ropto had for Biomechanics	S CONTRACTOR CONTRACTOR S CT. 80°	PODIEVA ECO 60 8 mm - 10 mm Carded Black	 Sheet: 900 x 790 mm 8 mm Carded Black 10 mm Carded Black
TECHNICAL MATERIALS			
STOCK		LIGHTEN CONFORT 130 2 mm Green	- Sheet: 490 mm x 550 mm - 2 mm Green

TOP LAYERS Thin materials

Recommended for:





TOP LAYERS Thin materials

Recommended for:






TOP LAYERS Decorative materials

42 www.podiatech.com




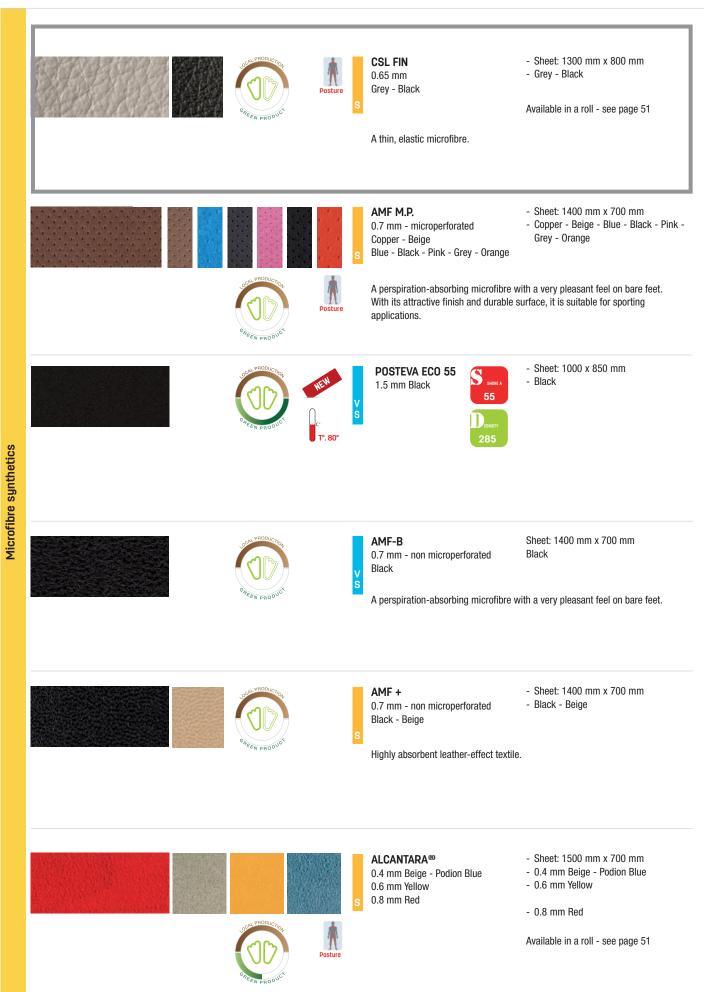

TOP LAYERS Decorative materials

Recommended for:




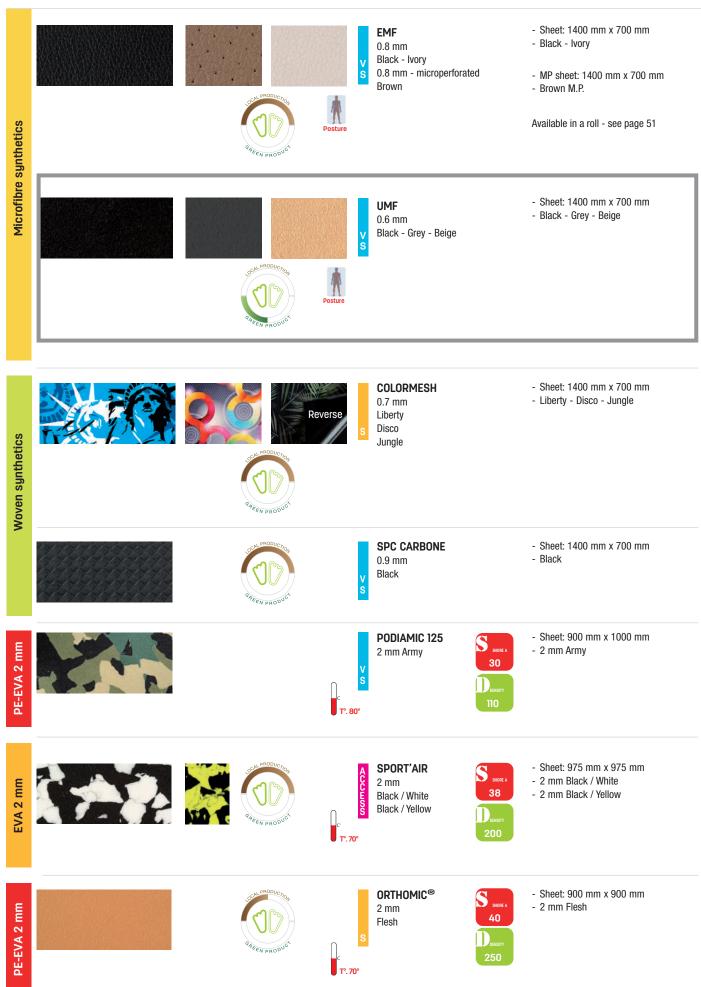

TOP LAYERS Antiperspirant materials

Recommended for:



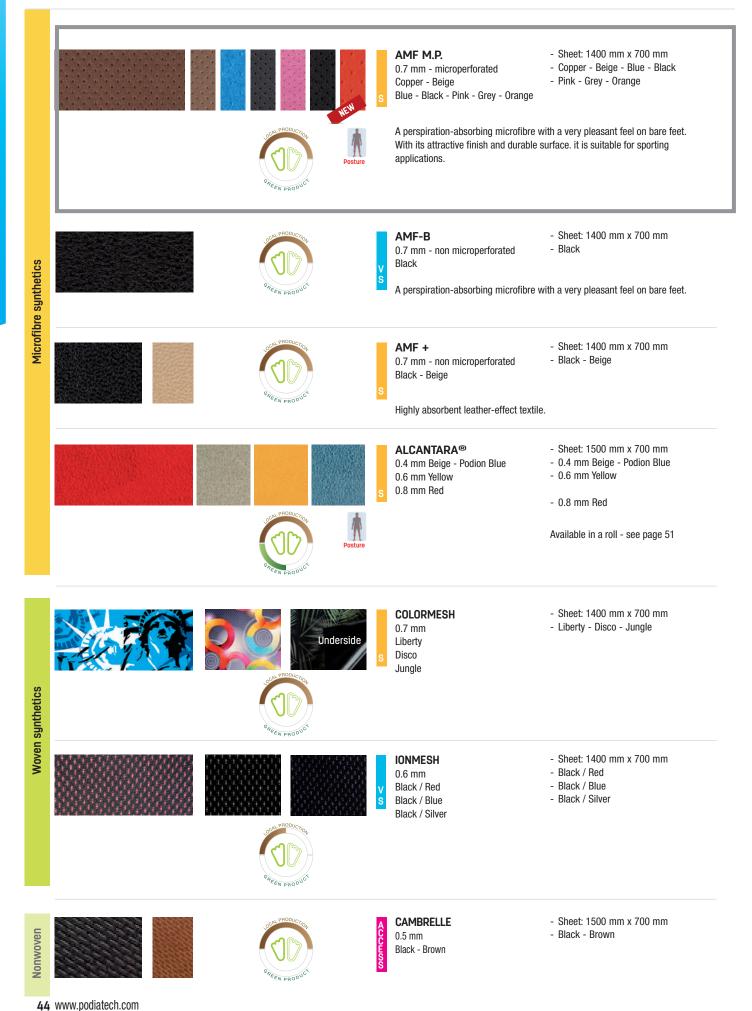

TOP LAYERS Durable materials

Recommended for:




				GIEN	mare	WOLK	aport
		ctus PRODucy of	S	AMF M.P. 0.7 mm - microperforated Copper - Beige Blue - Black - Pink - Grey - Orange	- C - P	heet: 1400 mm x 70 opper - Beige - Blu ink - Grey - Orange	e - Black
Microfibre synthetics		Posture		A perspiration-absorbing microfibre With its attractive finish and durable applications.			
		C. T. 80°	V S	POSTEVA ECO 55 1.5 mm Black		heet: 1000 x 850 m 5 mm Black	m
		SCOP PRODUCT OF	S	ALCANTARA® 0.4 mm Beige - Podion Blue 0.6 mm Yellow 0.8 mm Red	- 0 - 0	heet: 1500 mm x 70 .4 mm Beige - Podi .6 mm Yellow .8 mm Red	
crofibr		CAREEN PRODUCT Posture			Ava	ilable in a roll - see	page 51
Mic			V S	COSTO 0.8 mm Brown - Blue		heet: 1450 mm x 70 rown - Blue	00 mm
		SHEEN PRODUCTOR	V S	SPORTENE 3.2 mm Black Neoprene on lycra		heet: 550 mm x 600 .2 mm Black) mm
		PRODUCTION PRODUCTION PREEN PRODUCT	V S	RESISTENE 2.5 mm Black Elastomer DENSITY 200		heet: 1000 mm x 10 5 mm Black	000 mm
		CCAL PRODUCTION	V	SPC M.P. 0.8 mm - microperforated Black - Brown - Grey		heet: 1400 mm x 70 lack - Brown - Grey	
		GREEN PRODUCT	Š		Ava	ilable in a roll - see	page 51
				Recommended by chiropodists. Unl durability. Bactericidal and fungicid			sticity and
		CCM PRODUCTION		IONMESH 0.6 mm		heet: 1400 mm x 70 lack / Red	00 mm
Woven synthetics		SHEEN PRODUCT	V S	Black / Red Black / Blue Black / Silver	- B	lack / Blue lack / Silver	
			V S	SPC CARBONE 0.9 mm Black		heet: 1400 mm x 70 lack	00 mm
			V S	SPORTCOVER 0.8 mm Raspberry Beige - Black Grey - Brown Blue	- R - B	heet: 1400 mm x 7/ aspberry - Beige lack - Grey - Brown ilable in a roll - see	- Blue
		GREEN PRODUCT				boa.www	iatech.com 45

www.podiatech.com 45

TOP LAYERS Durable materials

Recommended for:

RAW MATERIALS

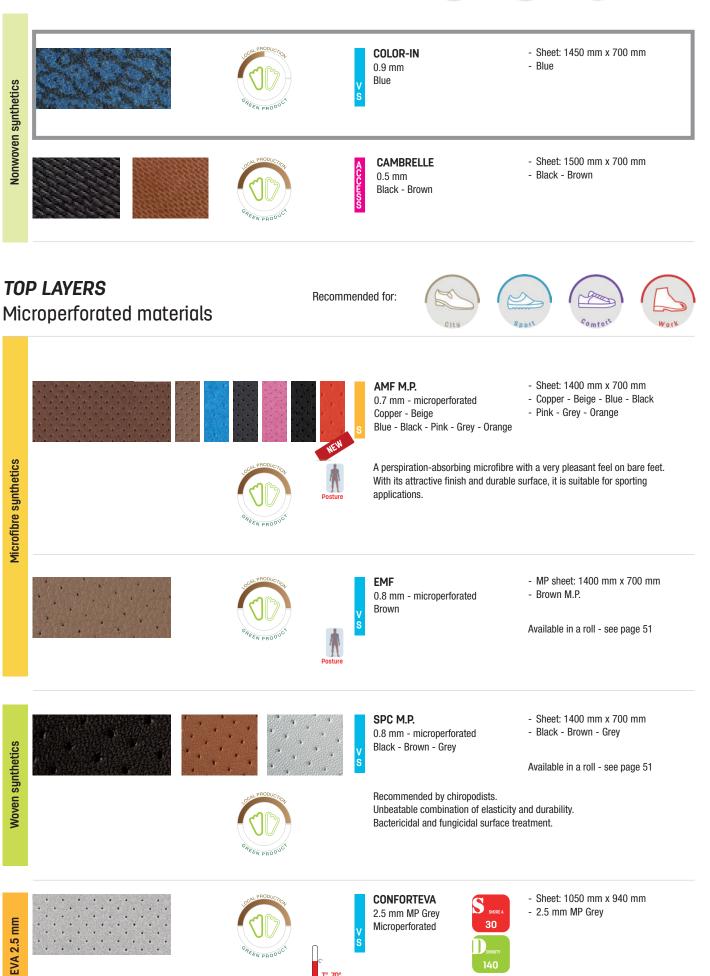


ENDROD

T°. 70

TOP LAYERS Microperforated materials




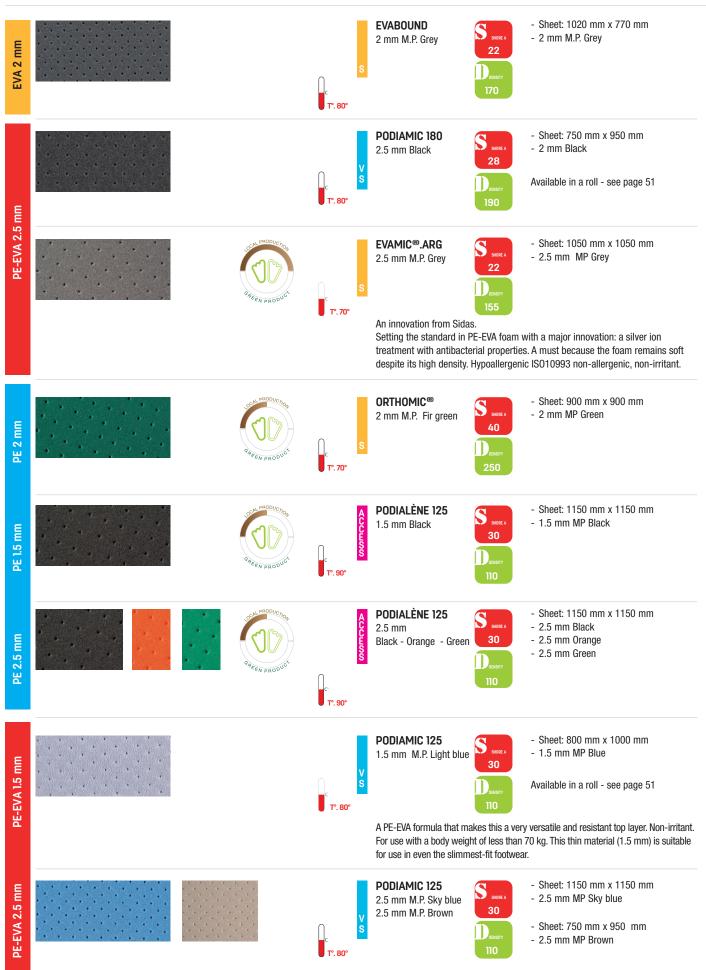


TOP LAYERS Microperforated materials


Recommended for:



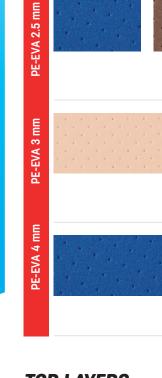



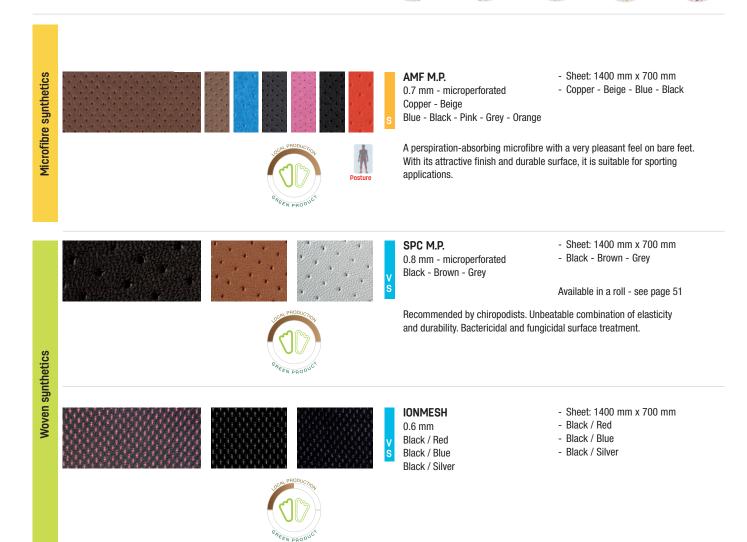


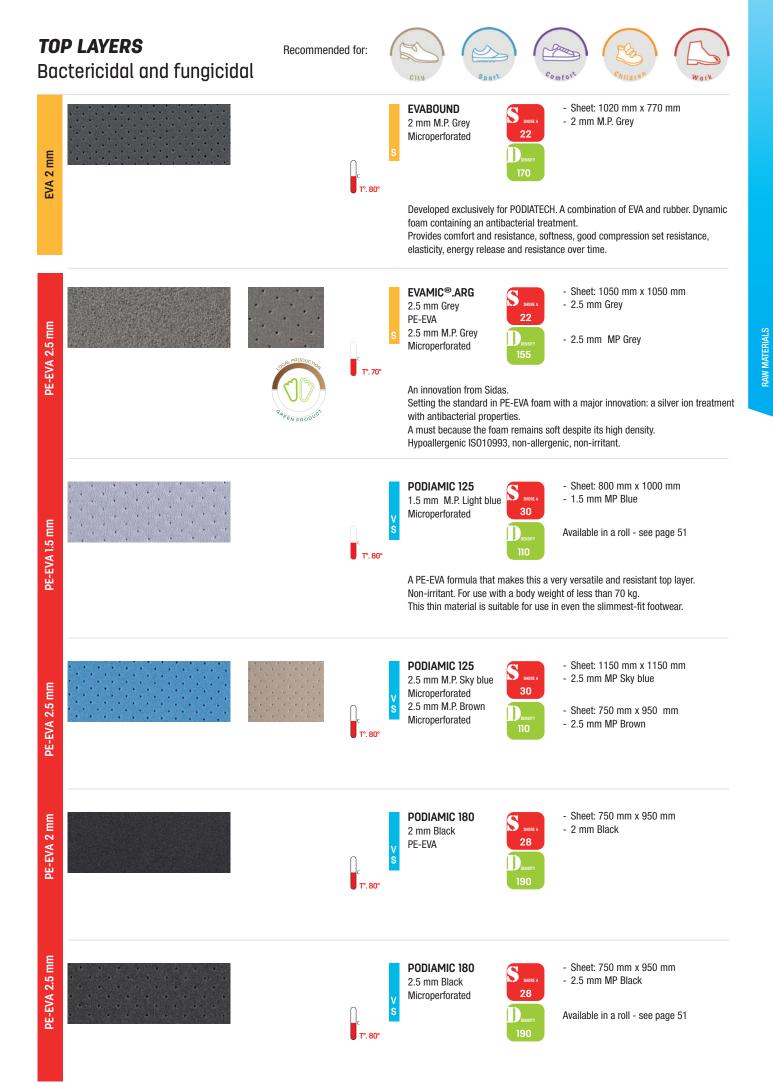

TOP LAYERS for highly specific uses

50 www.podiatech.com



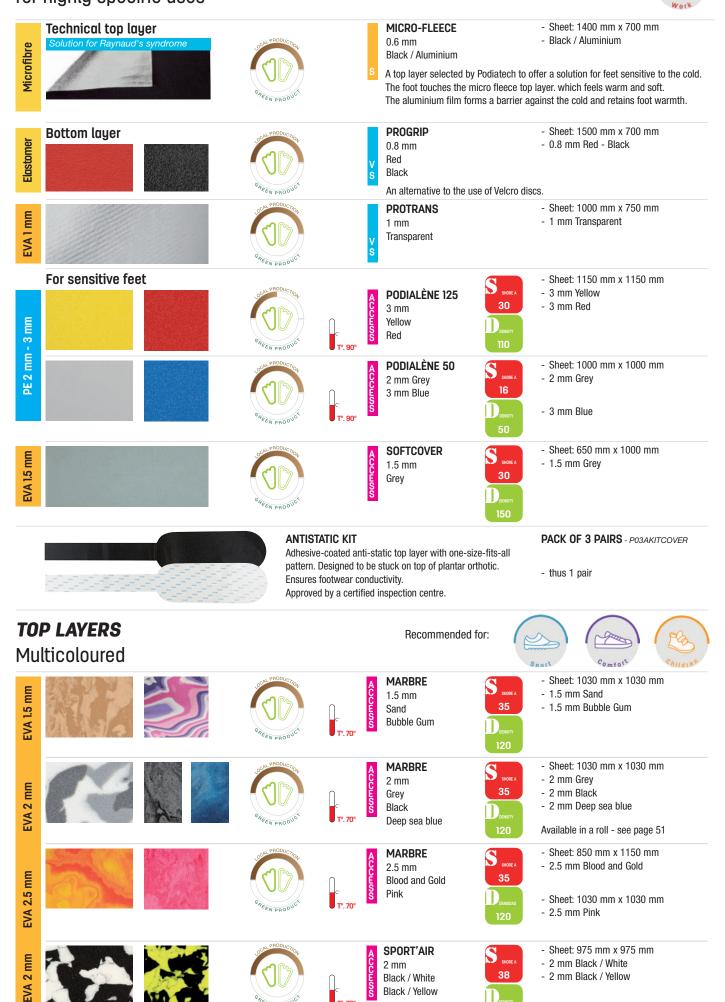

Black / Yellow

TOP LAYERS Multicoloured			Recommended for:	
EVA 1 mm	OF PRODUCTION	V S	PROTRANS 1 mm Transparent	- Sheet: 1000 mm x 750 mm - 1 mm Transparent
be-eva 2 mm	C.	V S 80°	PODIAMIC 125 2 mm Army 30 Deser 110	- 2 IIII Alliy
	PRODUCTION BREEN PRODUCT	ess	COLORMESH 0.7 mm Liberty Disco Jungle	- Sheet: 1400 mm x 700 mm - Liberty - Disco - Jungle
TOP LAYERS in rolls	Recommended for:			
	Sheen product	s	CSL FIN 0.65 mm Grey - Black A thin, elastic microfibre.	- Rolls: 12 cm x 10 m - Grey - Black
		v S	SPC M.P. 0.8 mm - microperforated Black - Brown - Grey	- Rolls: 12 cm x 10 m - Black - Brown - Grey
12 cm	CON PRODUCTION OFFICEN PRODUCTION OFFICEN PRODUCT	S S	ALCANTARA® 0.6 mm Podion Blue	- Roll: 12 cm x 10 m - Podion Blue
	PRODUCTION Preen provides	V S Ire	EMF 0.8 mm non microperforated Black non microperforated Brown microperforated Brown	Roll: 12 cm x 10 m - Black Roll: 12 cm x 10 m - Brown Roll: 12 cm x 10 m - Brown MP
	PRODUCTION PREEN PRODUCT	V S	AMF MP 0.7 mm Black Roll12 cm x 10 m	Roll: 12 cm x 10 m - Black
		s st V s	SPORTCOVER 0.8 mm Raspberry - Beige - Black- Grey Brown - Blue	- Rolls: 12 cm x 12.5 m - Raspberry - Beige y Black - Grey - Brown - Blue
30 cm	C T. 8	V S ^{BO°}	PODIAMIC 125 1.5 mm M.P. Light blue Microperforated	- 1.3 min wir blue
	OFFEEN PRODUCTION	A CCLLOS 70°	MARBRE 2 mm Black Deep sea blue	
C	C° T ". 6	V S 80°	PODIAMIC 180 2.5 mm Black Microperforated	

RESIN BASES From most flexible to most rigid



RAW MATERIALS

FLUENTE


Resins that fuse together to create more room. Do not require gluing or grinding.



A firm favourite for 25 years thanks to its qualities of robustness and flexibility.

RESIN BASES

From most flexible to most rigid

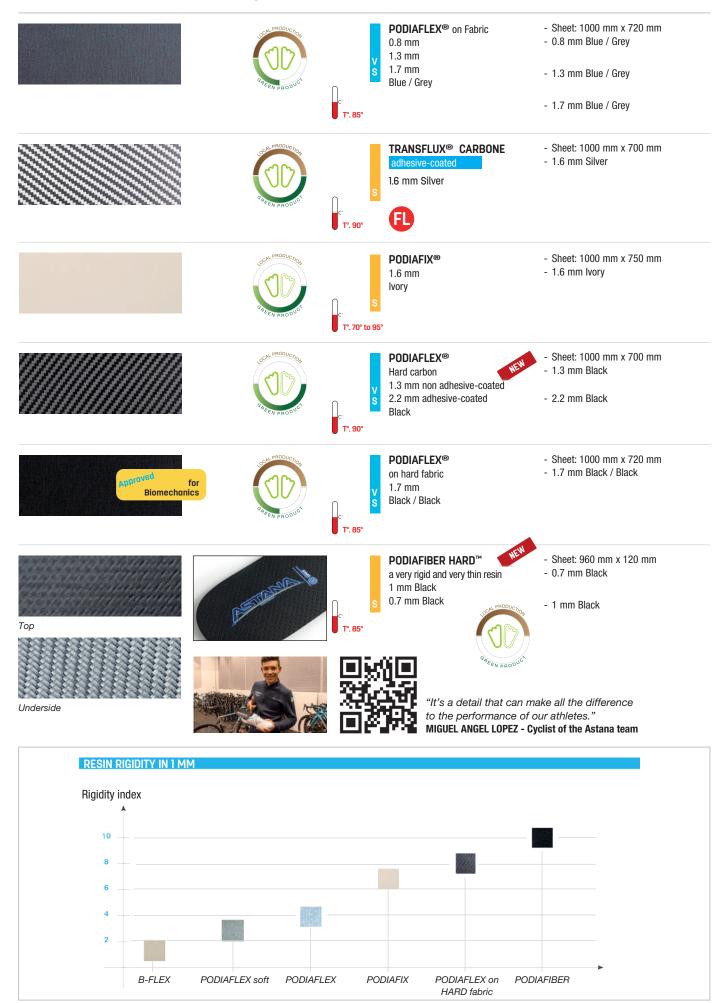

FOAM BASES thicknesses of 1 to 4 mm




FOAM BASES thicknesses of 5 to 20 mm

		CREEN PRODUCTION	C C T [°] . 100°	PODIALENE 160 B 5 mm 10 mm 15 mm White	 Sheet: 1250 mm x 1000 mm 5 mm White Sheet: 1150 mm x 1000 mm 10 mm White 1150 mm x 1000 mm 15 mm White
ΒE		CALL PRODUCTOR	C T*, 140°	Excellent thermoshaping. Can be PODIALENE 160 S 5 mm 10 mm 20 mm Sand	heat sealed to itself and to other materials. - Sheet: 1250 mm x 1000 mm - 5 mm Sand - Sheet: 1150 mm x 1000 mm - 10 mm Sand - Sheet: 1150 mm x 1000 mm - 15 mm Sand - Sheet: 1150 mm x 1000 mm - 20 mm Sand
		Original Production	V S	PODIALENE 200 6 mm Black 7 mm Blue 6 mm Black 0 mm Black 2 00 7 mm Blue 7 mm Blue 7 mm Blue	 Sheet: 950 mm x 950 mm 6 mm Black 7 mm Blue
		COL PRODUCTION GREEN PRODUCT	V C T°. 70°	CONFORTEVA 8 mm Light grey 30	- Sheet: 1050 mm x 940 mm - 8 mm Light grey
	Ap ^{proved} for Biomechanics	OCAL PRODUCTOR	۷ S ۳. 80°	PODIEVA 60 6 mm - 8 mm White Roughened	 Sheet: 920 mm x 810 mm 6 mm White 8 mm White
EVA	Roptoved for Biomechanics	OCAL PRODUCTOR	V C' T°. 80°	PODIEVA 65 6 mm 8 mm Red Roughened	- Sheet: 900 mm x 790 mm - 6 mm Red - 8 mm Red
	Approved for Biomechanics	COLUPRODUCTION COLORADO	Posture C' T'. 80°	PODIEVA 70 JEW 5 mm 6 mm 8 mm Black Roughened	 Sheet: 860 mm x 770 mm 5 mm Black 6 mm Black 8 mm Black
Ы			c S T°. 90°	PE SHORE 70 Roughened Non-Roughened 6 mm - Black Roughened 10 mm - Black Non-Roughened	 Sheet: 1100 mm x 840 mm 6 mm Black Roughened 10 mm Black Non-Roughened www.podiatech.com 55

COMFORT / CORRECTION ELEMENTS


Classified by Shore hardness

RAW MATERIALS

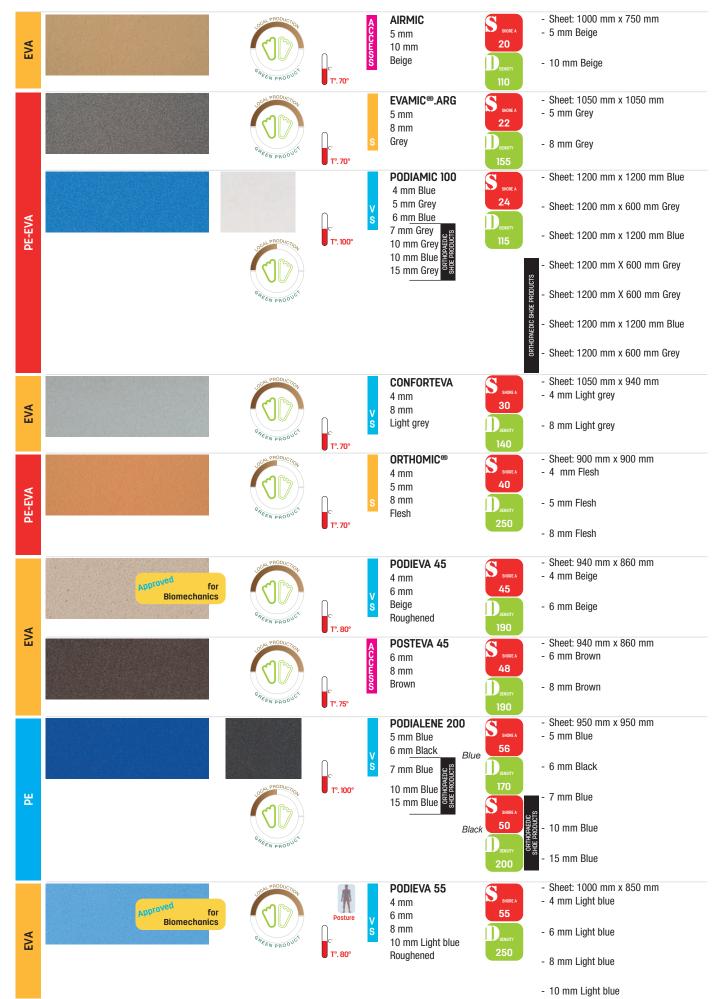

COMFORT / CORRECTION ELEMENTS

Classified by Shore hardness

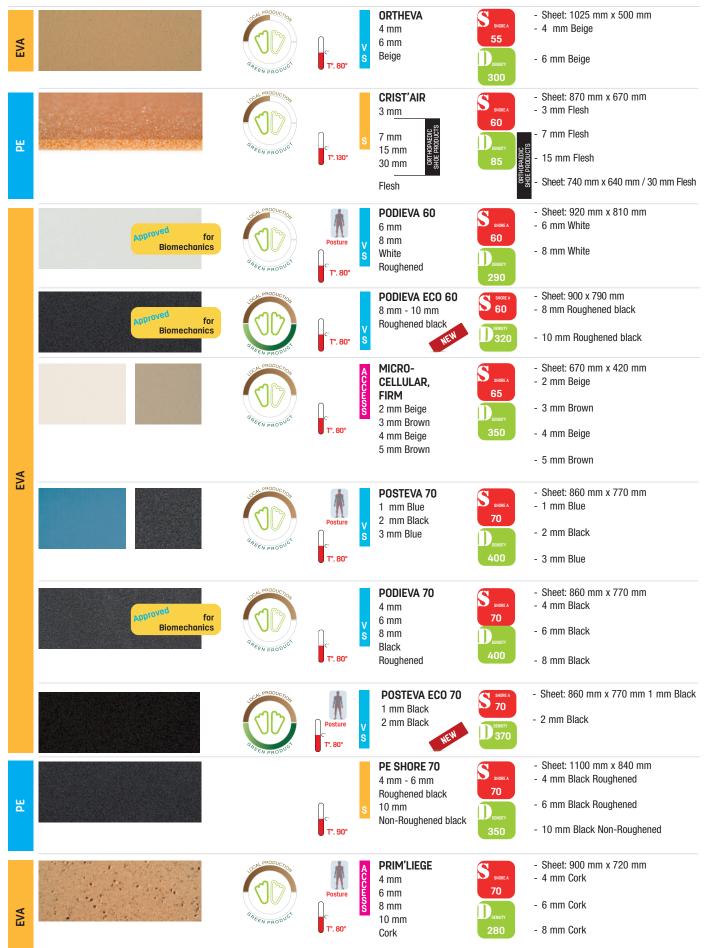

- 10 mm Cork

TECHNICAL MATERIALS

RAW MATERIALS



TECHNICAL MATERIALS

	OCAL PRODUCTION OFFICEN PRODUCT	PU	S	EXPERTENE 300 2.5 mm Green 4 mm Green	S INDE A 10 Density 300	- Sheet: 750 mm x 500 mm - 2.5 mm Green - 4 mm Green
	OCH PRODUCTOR	PU	S	EXPERTENE 400 1.5 mm Brown 4 mm Brown	SIGREA 10 Density 400	- Sheet: 750 mm x 500 mm - 1.5 mm Brown - 4 mm Brown
		PU	y s	PORON® DIABÈTE CARDE 2.5 mm Beige	S 5102E A 15 0653TY 270	- Sheet: 680 mm x 500 mm - 2.5 mm Beige
	 OCAL PRODUCTOR	E	S	DYNATENE 3 mm Orange	Sindle A 18 Density 420	- Sheet: 700 mm x 610 mm - 3 mm Orange
S		PU	v s	LIGHTEN DYNAMIQUE 130 4 mm Brown	SINDREA 15 Density 130	- Sheet: 490 mm x 550 mm - 4 mm Brown
All paintul areas		PU	v s	LIGHTEN DYNAMIQUE 300 2 mm Purple	SIGRE A 15 Density 300	- Sheet: 550 mm x 500 mm - 2 mm Purple
		PU	y s	LIGHTEN PROTECT 220 2 mm Green	S SHORE A 10 Density 220	- Sheet: 550 mm x 500 mm - 2 mm Green
		PU	v s	LIGHTEN DYNAMIQUE 160 3 mm Blue	Sinde A 15 Density 160	- Sheet: 490 mm x 550 mm - 3 mm Blue
		PU	v s	PORON® MS 3 mm Blue 6 mm Blue	S SHORE A 17 Density 330	 Sheet: 680 mm x 500 mm 3 mm Blue 6 mm Blue
	OCUL PRODUCTOR	E		PODITECH 3 mm Black 5 mm Black	S SHORE A 17 Density 500	- Sheet: 500 mm x 500 mm - 3 mm Black - 5 mm Black
	SCHL PRODUCTOR	PU		PRIM'TECH 3 mm Brown 5 mm Brown	S SHORE A 15 Denisty 350	 Sheet: 700 mm x 610 mm 3 mm Brown 5 mm Brown

TECHNICAL MATERIALS



Pack of 5 spare blades P99SCALPLAME

Fiskars® titanium scissors Thick materials, 21 cm P99CISRAZOREDGE **Fiskars**® Softouch® scissors Multi-purpose, 26 cm P99CISSOFTOUCH Template kit **CITY** lady

with small arch

Set of templates for Modulo range (Sizes 35/36 to 47/48) PDVCITY_LAD_AV_S R

(Sizes from 35 to 48) PDVGABMODULO



Insole templates 1 pack (Sizes 20 to 48) PDVGABSIDAS

60 www.podiatech.com

COMBINATION SHEETS





Bag of 50 ultra-thin black adhesive-coated "loop" Velcro® discs Ø 3 cm P99PAVELOURUF



Bag of 50 ultra-thin black adhesive-coated "hook" Velcro® discs Ø 3 cm P99PAVELCROUF



1 pack of 50 Velcro[®] discs P99VELCRO



Kit of 10 strips of ultra-thin Velcro® hooks and loops. 10 of each type (25 cm x 2.5 cm, adhesive-coated) P99VELCROKITBAND



RAW MATERIALS

Roll holder

PDVPRESENL12

Modulo Concept


This new way of presenting materials allows you to make your foot orthotics with a minimum of time spent in your workshop while keeping control of your patterns.

- Less storage and gluing
- Easy to use (rolls or strips)







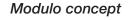



View our tutorials on our YouTube channel PODIATECH medical foot solutions

SET OF TEMPLATES FOR MODULO RANGE PDVGABMODULO

Sizes from 35 to 48

Bonus: your template This universal template allows you to cut the modulo plates in three patterns by size:

- Narrow template
- Standard template
- Wide template



Modulo - base/reinforcement RESIN/FOAM TAPE

Bases and reinforcements are cut out of your resin or foam strips with or without a pre-cut pad at the back of the foot.



12 cm

PODIAFLEX®

PODIAFLEX® on tissue with inclusion adhesive 3 sizes : P12MODRES01_35 (0.8 mm) / P12MODRES02_35 (1.3 mm) 35-39: 280 x 990mm P12MODRES01_40 / P12MODRES02_40 40-44: 310 x 875 mm P12MODRES01_44 / P12MODRES02_45 45-48: 335 x 980 mm

P12MODRES01

P12MODRES02

Description on page 53 of the catalogue.

PODIAFLEX®

in sheets/tapes 0.8 mm Blue / 1.2 mm Blue / 1.9 mm Blue Sheet: 1000 mm x 750 mm

Tapes: 12 cm x 1 m

Description on page 52 of the catalogue.

PODIALENE 200

in strip 6 mm Black *P12MODMOUSSE07_35* : 280 x 950 mm *P12MODMOUSSE07_40* : 310 x 950mm *P12MODMOUSSE07_45* : 335 x 950 mm

Description on page 53 of the catalogue.

62 www.podiatech.com

MODULO - INTERLAYER

MIDDLE FOAM TAPES

The intermediate strips in rolls respect the different needs of the rearfoot, midfoot and forefoot.

A real innovation that allows you to make your orthotics independently of the reason of consultation or the volume of the patient's footwear. Available in 3 different heights depending on the size:

- 35/39: 280 mm
- 40/44: 310 mm
- 45/48: 335 mm



ROLL COVERS

Customize your orthotics with our roll covers p.51 of the catalogue.

Discover with modulos a simplified protocol for the realization of your plantar orthotics.

Standard process



Thermo-galvanic process











ORTHOTICS & ELEMENTS



Intelligent top cover

30 pairs minimum

Integrate the intelligence of OPCT® into your milling, polypro covers, 3D printed covers, etc.

PODIATECH has developed for you a range of 8 intelligent topcover made of a combination of selected materials and adapted to all pathologies and all shoe volumes.

CITY RANGE



CITY LOW NARROW - P25CITY_LOW_NA

CITY LOW WIDE - P25CITY_LOW_WI

- AMF (Top cover)

- PE 180 2 mm

- PODIANE+ (primary layer)

PATHO RANGE




CITY STANDARD NARROW - P25CITY_STD_NA

CITY STANDARD WIDE - P25CITY_STD_WI

- EMF 0.8 mm

- ORTHOMIC 3 mm
- PODIANE+ 4 mm

Sizes: All our models are available in a narrow or wide in the following sizes 35/36 - 37/38 - 39/40 - 41/42 - 43/44 -45/46 - 47/48

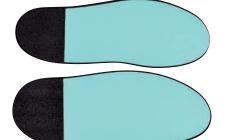

PATHO LOW NARROW - P25PATH_LOW_NA

PATHO LOW WIDE - P25PATH LOW WI

- PODIAMIC 150 1 mm Black

- PORON MEDICAL 2 mm

- PODIANE+ 2 mm

SPORT RANGE



PATHO STANDARD NARROW - P25PATH_STD_NA

PATHO STANDARD WIDE - P25PATH STD WI

- CONFORTEVA 2.5 mm
- LIGHTEN PROTECT 220
- PODIANE+ (primary layer)



PATHO VOLUME NARROW - P25PATH_VOL_NA

PATHO VOLUME WIDE - P25PATH VOL WI

- CONFORTEVA (top cover) - LIGHTEN DIABETE 160 4 mm (primary layer)



SPORT EXTRA LOW NARROW - P25SPT XLOW NA

SPORT EXTRA LOW WIDE - P25SPT_XLOW_WI

- PODIAMIC 150 1 mm Noir


- VISCOTENE 1.5 mm
- PE SHORE 1 mm
- NOENE 1 mm



SPORT LOW NARROW - P25SPT LOW NA

SPORT LOW WIDE - P25SPT_LOW_WI

- PODIALENE 125 1.5 mm
- PROPULSE 2 mm
- PE 200 2 mm
- PROTECT 2 mm

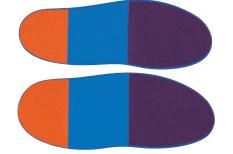

SPORT STANDARD NARROW - P25SPT STD NA

SPORT STANDARD WIDE - P25SPT_STD_WI

PODIAMIC 160 Blue

- LIGHTEN DYNAMIQUE 300 2 mm
- PODIALENE 200 2 mm
- PROTECH 2 mm



OPCT® – an innovative concept created by PODIATECH



PersonaliSand

with your logo (minimum

A constantly

evolving range

from €8/pair


incl VAT

Materials and shapes

customiSand to

your own specifications

Options to suit the way you work. Each $OPCT^{\otimes}$ comprises several layers: top layer + primary layer + base + reinforcements. PODIATECH makes different versions of the $OPCT^{\otimes} - S$, SM, O and OM – to match the level of customisability the practitioner wants to offer his patients and the amount of time he is prepared to spend in the workshop.

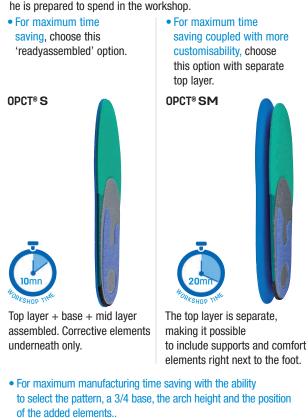



Top layer + mid layer assembled. Base separate from the assembly.




The top layer + mid layer + base are all separate from each other, making it possible to include supports and comfort elements right next to.

A thermowelding stage is required. The pattern used for the base is wider than for the S versions, enabling the professional to adapt the pattern.

Traceability





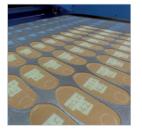







A solution

for every patient

PODIATECH OPCT®

Orthotics tested and

approved by athletes

and elite user groups

(GIGN special ops

forces)

Time-saving, enabling

the practitioner

to focus on clinical

Watch our video explaining

the OPCT® concept


Developed

in France b

PODIATECH's R&D department

www.podiatech.com 65

OPCT® CITY RANGE

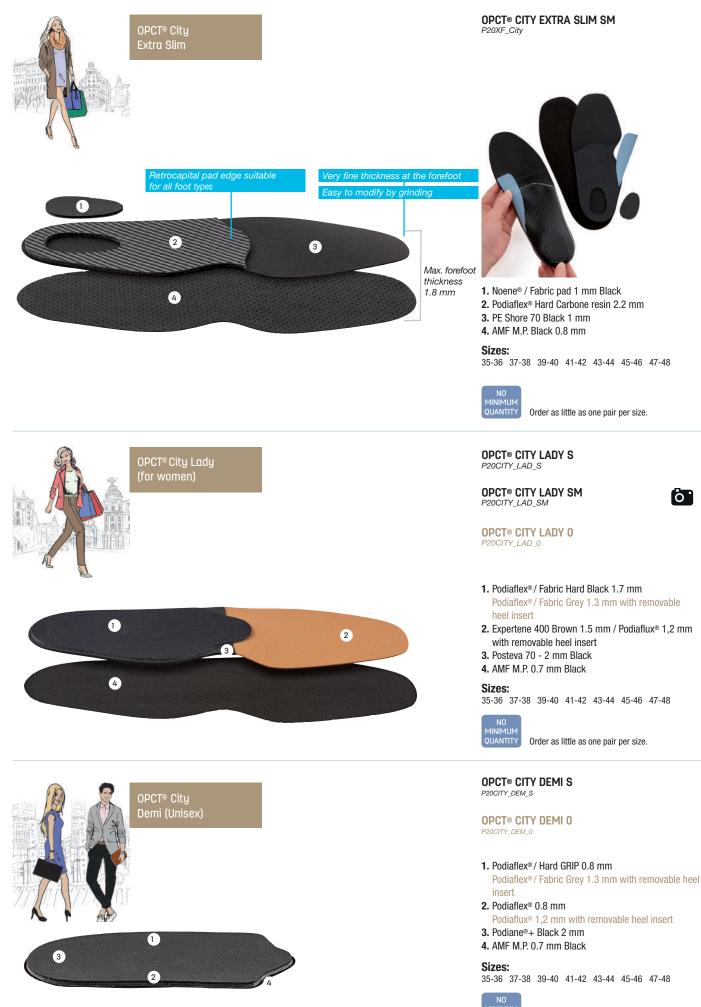

Podiatech reserves the right to change the materials used to make OPCT[®] products **66** www.podiatech.com

OPCT® CITY RANGE











Podiatech reserves the right to change the materials used to make OPCT® products

OPCT[®] CITY MAN S P20CITY_MAN_S

OPCT[®] CITY MAN SM P20CITY_MAN_SM


OPCT® CITY MAN O P20CITY_MAN_O

0. Podiaflex $^{\textcircled{B}}$ / Fabric Grey 1.3 mm with removable heel insert

- 1. Podiaflex R / Fabric Hard Black 1.7 mm
- Podiaflux $^{\textcircled{R}}$ 1.2 mm with removable heel insert 2. Expertene 400 Brown 1.5 mm
- **3.** Posteva 70 Black 2 mm
- 4. Podiane+ Black 2 mm
- 5. AMF M.P. Black 0.7 mm

Sizes:

35-36 37-38 39-40 41-42 43-44 45-46 47-48

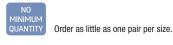

OPCT[®] EVERYDAY SLIM S P20EVRD_FIN_S

OPCT[®] EVERYDAY SLIM SM P20EVRD_FIN_SM

OPCT[®] EVERYDAY SLIM O P20EVRD_FIN_0

OPCT[®] EVERYDAY SLIM OM P20EVRD_FIN_oM

0. Podiaflex[®] / Fabric Grey 1.3 mm with removable heel insert

- 1. Podiaflex[®] / Fabric hard Black 1.7 mm / Podiaflux[®] 1.2 mm with removable heel insert
- 2. Podiane®+ Black 2 mm
- 3. Orthomic Flesh 2 mm
- 4. EMF 0.8 mm Black

Sizes:

35-36 37-38 39-40 41-42 43-44 45-46 47-48



Order as little as one pair per size.

OPCT[®] EVERYDAY SLIM TRANSFLUX® O *P20EVRD_FTX_O*

1. Transflux® Charcoal 1.5 mm


- 2. Podiane®+ Black 2 mm
- 3. Orthomic Flesh 2 mm
- 4. EMF 0.8 mm Black

Sizes: 35-36 37-38 39-40 41-42 43-44 45-46 47-48



Order as little as one pair per size. www.podiatech.com 67

OPCT® CITY RANGE

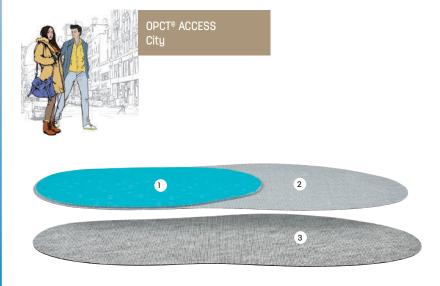

OPCT® MULTI RANGE



OPCT[®] CITY MULTI S P20CITY_MUL_S

Sizes:

Sizes:

MINIMUM OUANTITY

OPCT[®] ACCESS CITY S *P23ACC_VIL_S*

OPCT[®] ACCESS CITY SM *P23ACC_VIL_SM*

1. Podiaflex® Color resin Blue 1.3 mm 2. Podiaflex[®] Soft resin Grey 0.8 mm 3. Sportcover Black 0.8 mm

35-36 37-38 39-40 41-42 43-44 45-46 47-48

Order as little as one pair per size.

OPCT[®] CITY MULTI SM P20CITY_MUL_SM

0.

0.

1. Podiaflex[®] / Fabric Hard Black 1.7 mm

35-36 37-38 39-40 41-42 43-44 45-46 47-48

Order as little as one pair per size.

- 2. Podiaflex[®] / Fabric 0.8 mm
- 3. Podialène M.P. 125 2.5 mm Black



2 1 3 4

OPCT® Everyday

Podiatech reserves the right to change the materials used to make OPCT® products 68 www.podiatech.com

OPCT[®] EVERYDAY STANDARD S P20EVRD_STD_S

OPCT[®] EVERYDAY STANDARD SM P20EVRD_STD_SM



OPCT[®] EVERYDAY STANDARD OM P20EVRD_STD_OM

- 0. Podiaflex[®] / Fabric Grey 1.3 mm with removable heel insert
- 1. Podiaflex[®] / Fabric hard 1.7 mm Podiaflux® 1.2 mm with removable heel insert
- 2. Podiaflex® / Fabric 1.3 mm
- 3. Conforteva M.P. 2.5 mm Grey
- 4. EMF 0.8 mm Black

Sizes:

35-36 37-38 39-40 41-42 43-44 45-46 47-48



Order as little as one pair per size.


OPCT® MULTI SPORT RANGE



OPCT® Extraslim SPORT




OPCT® EXTRASLIM SPORT SM P20XF_SPORT_SM

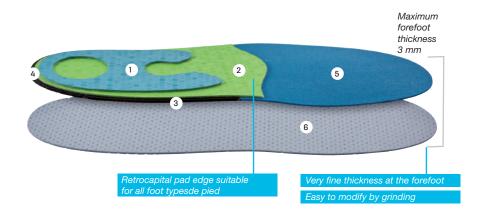

- 1. Podiaflex® Color resin Blue 1.3 mm
- 2. Podiaflex® Color resin Green 1.3 mm
- 3. PE Shore 70 Black 2 mm
- 4. Podiane® pad Black 2 mm
- 5. Viscotène Blue 1.5 mm
- 6. Podiamic 125 M.P. Light Blue 1.5 mm

Sizes:

35-36 37-38 39-40 41-42 43-44 45-46 47-48

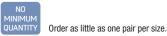

OPCT® ACCESS SPORT S P23ACC_SPOR_S

OPCT® ACCESS SPORT SM P23ACC_SPOR_SM 0.

1. Podiaflex[®] Color resin Blue 1.3 mm

- 2. Podiaflex® Color resin Green 1.3 mm
- 3. Poron® XRD Yellow 3 mm
- 4. Dynatene Orange 3 mm
- 5. Podialene® 200 Black 3 mm
- 6. Podiamic 180 Black 2 mm

Sizes: 35-36 37-38 39-40 41-42 43-44 45-46 47-48



Order as little as one pair per size. www.podiatech.com **69** ORTHOTICS & ELEMENTS

2

6

3

OPCT® Access Sport

OPCT® MULTI SPORT RANGE













Podiatech reserves the right to change the materials used to make ${\sf OPCT}^{\circledast}{\sf products}$ ${\bf 70}$ www.podiatech.com

OPCT® SPORT SLIM S P20MULTI_SLIM_S

OPCT® SPORT SLIM SM P20MULTI_Slim_SM

OPCT® SPORT SLIM O P20MULTI_SLIM_O

OPCT® SPORT SLIM OM P20MULTI_Slim_OM

- $\textbf{0.}\ Podiaflex^{\circledast}/\ Fabric\ Grey\ 1.3\ mm\ with\ removable\ heel\ insert$
- 0. Podiaflux® 1.2 mm with removable heel insert
- 1. Podiaflex® / Fabric Hard Black 1.7 mm
- **2.** Posteva 70 Blue 1 mm
- 3. Podiane®+ 2 mm
- 4. Podiamic 180 Black 2 mm

Sizes:

35-36 37-38 39-40 41-42 43-44 45-46 47-48

NO MINIMUM QUANTITY Order as little as one pair per size.

OPCT® SPORT STANDARD S P20MULTI_STD_S

OPCT® SPORT STANDARD SM

P20MULTI_STD_SM

OPCT® SPORT STANDARD 0 P20MULTI_STD_O

OPCT® SPORT STANDARD OM P20MULTI_STD_OM

- 0. Podiaflex[®] / Fabric Grey 1.3 mm with removable heel insert
- 0. Podiaflux® 1.2 mm with removable heel insert
- 1. Podiaflex® / Fabric Hard Black 1.7 mm (only in S and SM version)
- **2.** Viscotene $^{\ensuremath{\mathbb{R}}}$ Grey 2 mm
- 3. Podieva 70 Black 2 mm
- 4. Podiane®+ Black 2 mm
- 5. Podiamic 180 MP Black 2.5 mm

Sizes:

35-36 37-38 39-40 41-42 43-44 45-46 47-48

N0 MINIMUM

QUANTITY Order as little as one pair per size.

OPCT® SPORT WIDE S P20MULTI_WIDE_S

OPCT® SPORT WIDE SM P20MULTI_Wide_SM

OPCT® SPORT WIDE O P20MULTI_WIDE_O

OPCT® SPORT WIDE OM P20MULTI_Wide_OM


- 0. Podiaflex[®] / Fabric Grey 1.3 mm with removable heel insert
- 0. Podiaflux® 1.2 mm with removable heel insert
- 1. Podiaflex[®]/Fabric Hard Black 1.7 mm (only in S and SM version)
- 2. Jogtene[®] 3 mm Black
- 3. PE 200 Blue 3 mm
- 4. Podiane[®]+ Black 4 mm
- **5.** Podiamic 180 MP Black 2.5 mm
 - FUUIAITIIC TOU IVIF DIACK 2.3 ITII

Sizes: 35-36 37-38 39-40 41-42 43-44 45-46 47-48

NO

Order as little as one pair per size.

OPCT® MULTI SPORT RANGE

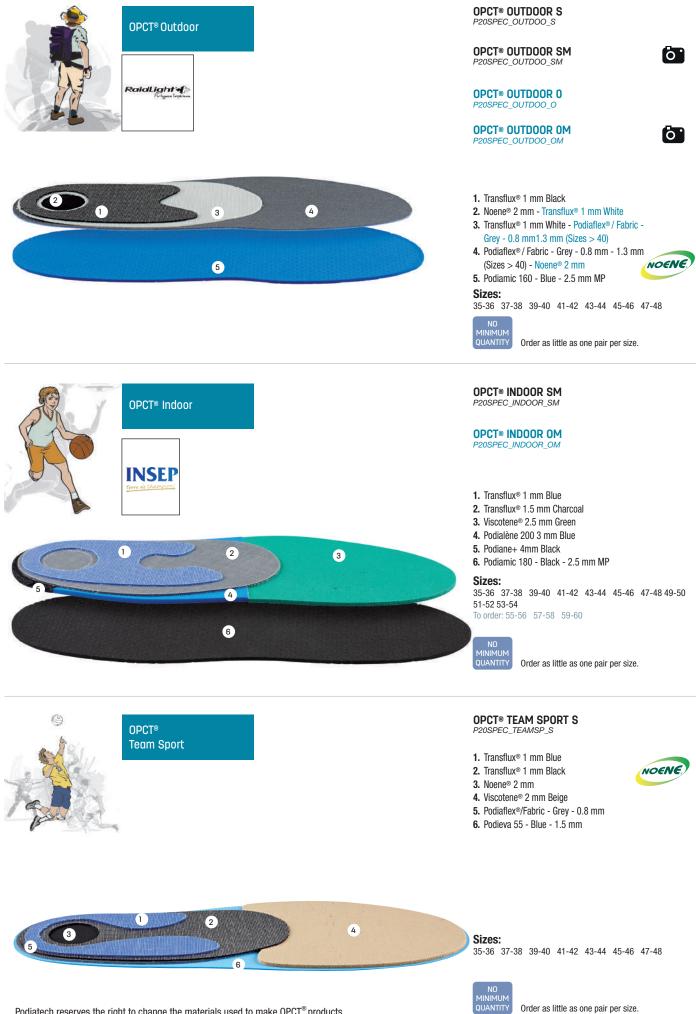


Podiatech reserves the right to change the materials used to make OPCT® products

www.podiatech.com 71

ORTHOTICS & ELEMENTS

OPCT® SPECIFIC SPORTS RANGE

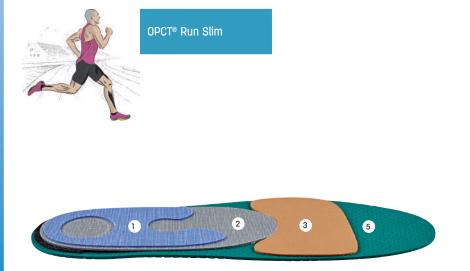

OPCT® RUN SLIM O P20SPEC_RUNSLI_O

- 1. Transflux® Blue 1 mm
- 2. Transflux® Charcoal 1.5 mm
- 3. Expertene Brown 1.5 mm
- 4. Podiane®+ Black 2 mm
- 5. Orthomic® Green 2 mm MP


Sizes:

35-36 37-38 39-40 41-42 43-44 45-46 47-48



Order as little as one pair per size.

ο.

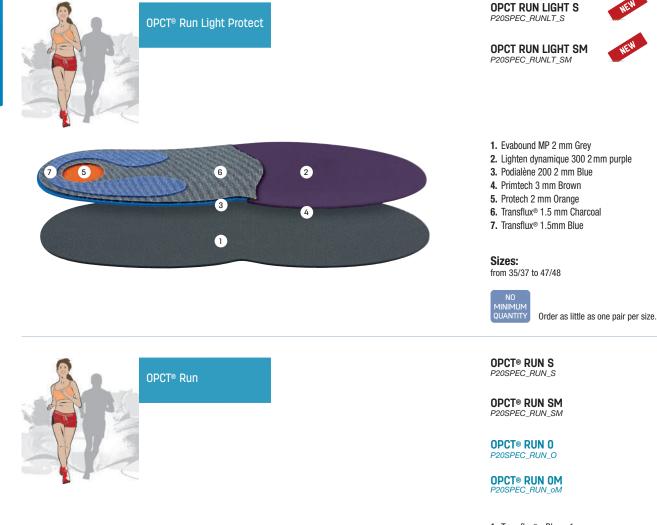

1 2 3 4 6

- 1. Evabound MP 2 mm Grey
- 2. Lighten dynamique 300 2 mm purple
- 3. Podialène 200 2 mm Blue

- 1. Transflux® Blue 1 mm
- 2. Transflux® Charcoal 1.5 mm
- 3. Viscotene® 2.5 mm Green
- 4. Orthomic® 3 mm Green
- 5. Podiane®+ Black 4 mm
- 6. Podiamic 160 Blue 2.5 mm MP

Sizes:

35-36 37-38 39-40 41-42 43-44 45-46 47-48



Order as little as one pair per size.

OPCT® SPECIFIC SPORTS RANGE



www.podiatech.com 73

ORTHOTICS & ELEMENTS

3

5

Podiatech reserves the right to change the materials used to make OPCT® products

4

74 www.podiatech.com




6

7

Sizes: 35-36 37-38 39-40 41-42 43-44 45-46 47-48




Order as little as one pair per size.

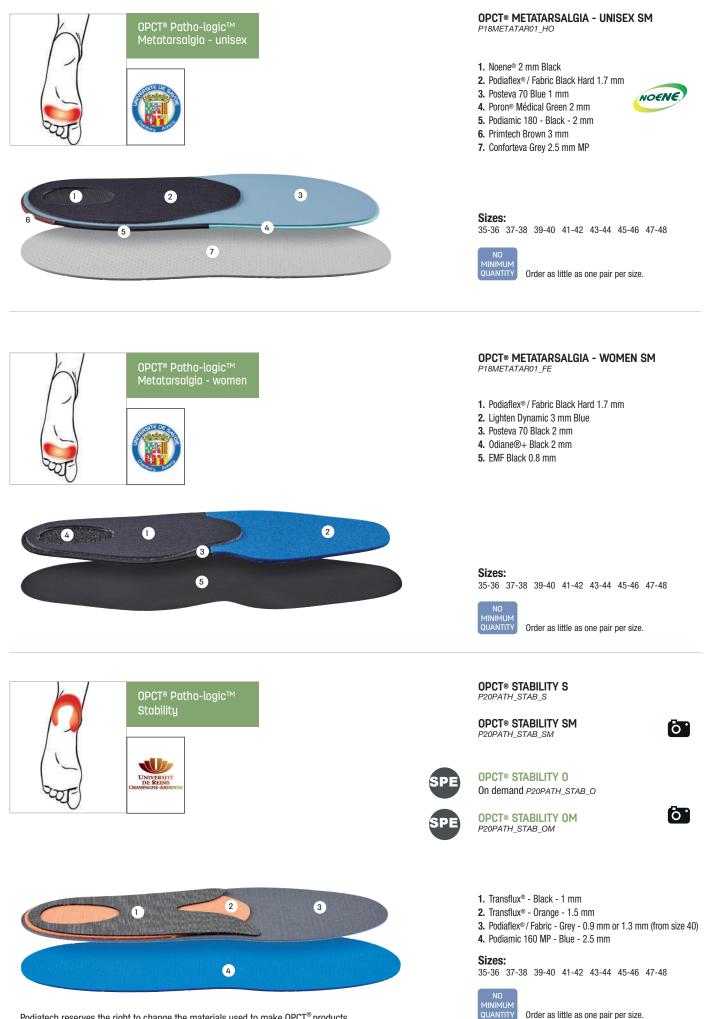

Podiatech reserves the right to change the materials used to make OPCT® products

www.podiatech.com 75



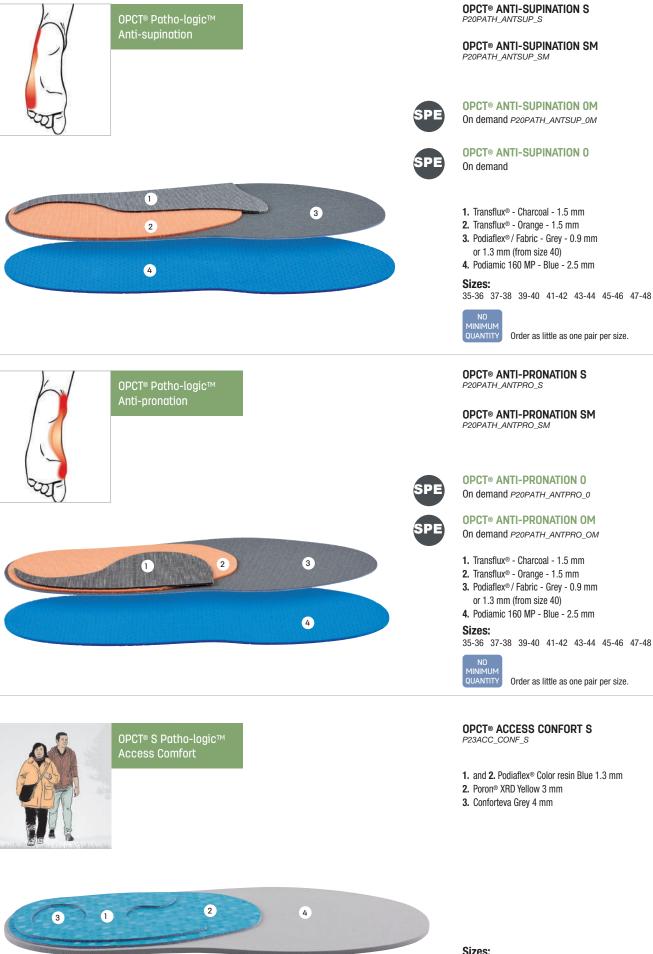


Podiatech reserves the right to change the materials used to make OPCT® products 76 www.podiatech.com

Sizes: 35-36 37-38 39-40 41-42 43-44 45-46 47-48

MUN QUANTITY ο.

0.

ο.

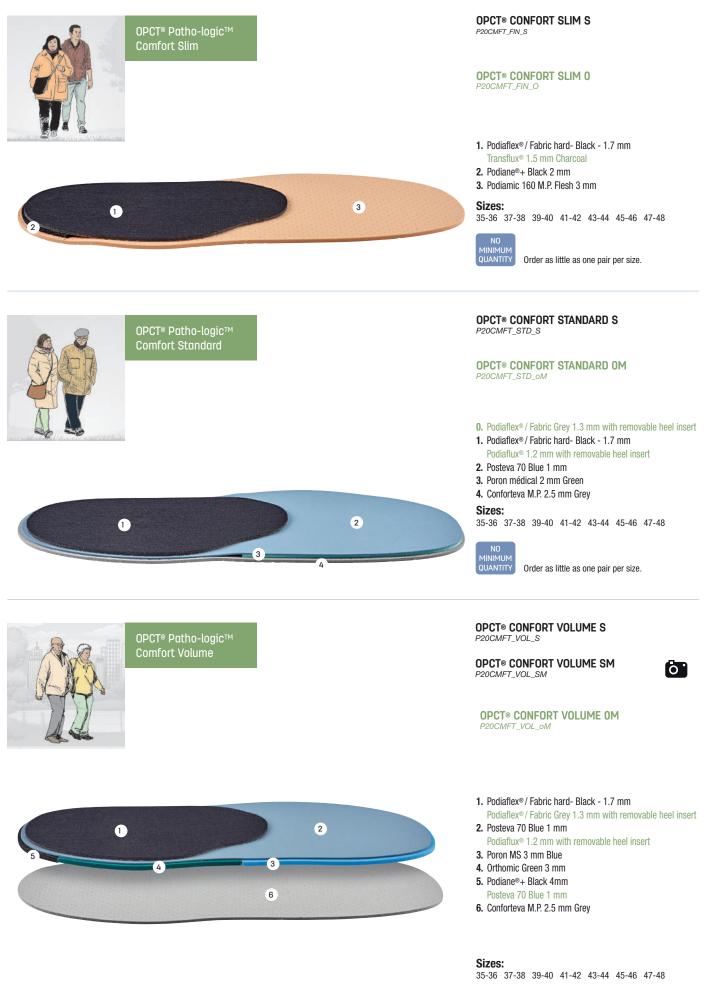

Podiatech reserves the right to change the materials used to make OPCT[®] products

IMUN

Order as little as one pair per size. WWW.podiatech.com **77** ORTHOTICS & ELEMENTS


OPCT® KIDS RANGE







OPCT® Sport

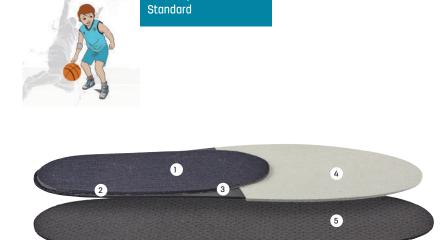

Podiatech reserves the right to change the materials used to make $\mathsf{OPCT}^{\circledast}\mathsf{products}$ 78 www.podiatech.com

OPCT® SEVER S P20PATH_SEVER_S

OPCT® SEVER SM P20PATH_SEVER_SM

ο.

1. Podialène 160 - White - 5mm

- 2. Podiane®+ 2 mm
- 3. Podiaflex® / Fabric Black 1.7 mm
- 4. Podiane®+ Black 2 mm
- 5. Orthomic Flesh 2 mm
- 6. Viscotene® Grey 2 mm
- 7. Podiamic 180 Black 2 mm

Sizes:

33/34 - 35/36 - 37/38 39/40 - 41/42 - 43/44 - 45/46 - 47/48



Order as little as one pair per size.

OPCT® EVERYDAY FIN S P20EVRD_FIN_S

OPCT® EVERYDAY FIN SM P20EVRD_FIN_SM



1. Podiaflex® / Fabric hard - Black - 1.7 mm

- 2. Podiane®+ Black 2 mm
- 3. Orthomic Flesh 2 mm
- 4. EMF 0.8 mm Black

Sizes:

29/30 - 31/32 - 33/34



Order as little as one pair per size.

OPCT® SPORT STANDARD S P20MULTI_STD_S

OPCT® SPORT STANDARD SM P20MULTI_STD_SM

ο.

1. Podiaflex® / Fabric - Black - 1.7 mm

- 2. Podiane®+ Black 2 mm
- 3. Podieva 70 Black 2 mm
- 4. Viscotene® Grey 2 mm
- 5. Podiamic 180 MP Black 2.5 mm

Sizes: 29/30 - 31/32 - 33/34



OPCT® KIDS RANGE



Combination materials in stock. Quotations for custom combination materials available on request



OPCT® SPECIFIC BIOMECHANISM



OPCT® KIDS PINK S

1. Podiaflex® / Fabric Hard Black 1.7 mm

2. Poron® médical 2 mm

3. Marbré Pink 2.5 mm

OPCT® KIDS BLUE S P13GAR01

1. Podiaflex® / Fabric Hard Black 1.7 mm

2. Poron® médical 2 mm

3. Podiamic 125 MP Blue 2.5 mm

Sizes: 29/30 - 31/32 - 33/34



OPCT® SPORT STANDARD O P20MULTI_STD_0

OPCT® SPORT STANDARD OM P20MULTI_STD_OM


1. Podiaflex® / Fabric Grey 1.3 mm with removable heel insert

Order as little as one pair per size.

- 2. Podiaflux® 1.2 mm with removable heel insert
- 3. Podiane®+ Black 2 mm
- 4. Podieva 70 Black 2 mm
- 5. Viscotene® Grey 2 mm
- 6. Podiamic 180 MP Black 2.5 mm

Sizes:

29/30 - 31/32 - 33/34

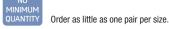

HYPER DRESS P20HYPER_DRESS

- Only available in OM version
- Available on demand
- 1. Podiaflex Blue 1.5 mm
- 2. Podiane 2 mm
- 3. Glass fiber Black 1 mm

HYPER MEDI P20HYPER_MEDI

- Only available in OM version

- Available on demand
- 1. Podialène 125 MP Orange 2.5 mm
- 2. Protech Orange 2 mm
- 3. Podiaflex Blue 1.9 mm
- 4. Glass fiber Black 1 mm
- 5. Glass fiber Black 1 mm

HYPER SPORT P20HYPER_SPORT

- Only available in OM version
- Available on demand
- 1. Lighten dynamique 300 Purple 2 mm
- 2. Protech Orange 2mm
- 3. Glass fiber Black 1 mm
- 4. Podiamic 160 MP Blue 2.5 mm
- 5. Podiaflex Blue 1.9 mm


OPCT® MOUSSE RANGE FOR SENSITIVE FEET: ORTHO CONFORT




OPCT® MOUSSE RANGE FOR SENSITIVE FEET: ORTHO CONFORT

A range of light, pre-designed and pre-moulded comfort orthotics with met bar, customisable to the patient's plantar fat pad condition: atrophied, thin,moderate or normal.



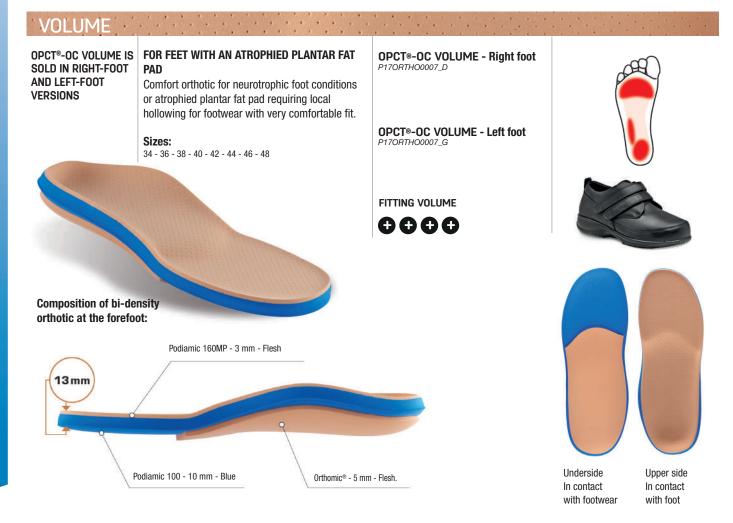

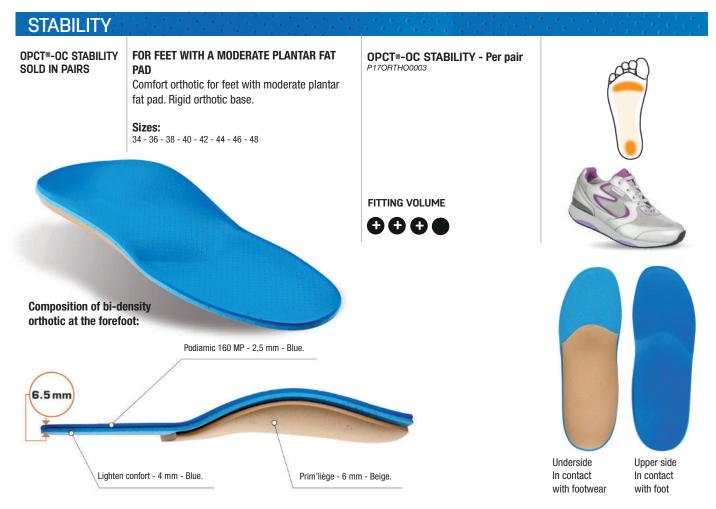



CUSHIONING PREMIUM

OPCT®-OC CUSHIONING PREMIUM SOLD IN PAIRS

Composition of tri-density orthotic at the forefoot:

8.5 mm

FOR FEET WITH A VERY THIN PLANTAR FAT PAD Comfort orthotic for sensitive feet with a very thin plantar fat pad. Variation of the OPCT®-OC COMFORT+ with a very soft forefoot. Rigid orthotic base.

Sizes: 34 - 36 - 38 - 40 - 42 - 44 - 46 - 48

Podiamic 160 MP - 1,5 mm - Brown.

Airmic - 5 mm - Beige.

OPCT®-OC CUSHIONING PREMIUM Per pair P17ORTHO0006

FITTING VOLUME $\Theta \oplus \Theta \bullet$

Poron® médical - 2 mm - Green.

Prim'liège - 6 mm - Beige.




Underside In contact with footwear

with foot

Upper side In contact

FOAM BASE SOLUTIONS

Pre-cut bases - Foam base modules

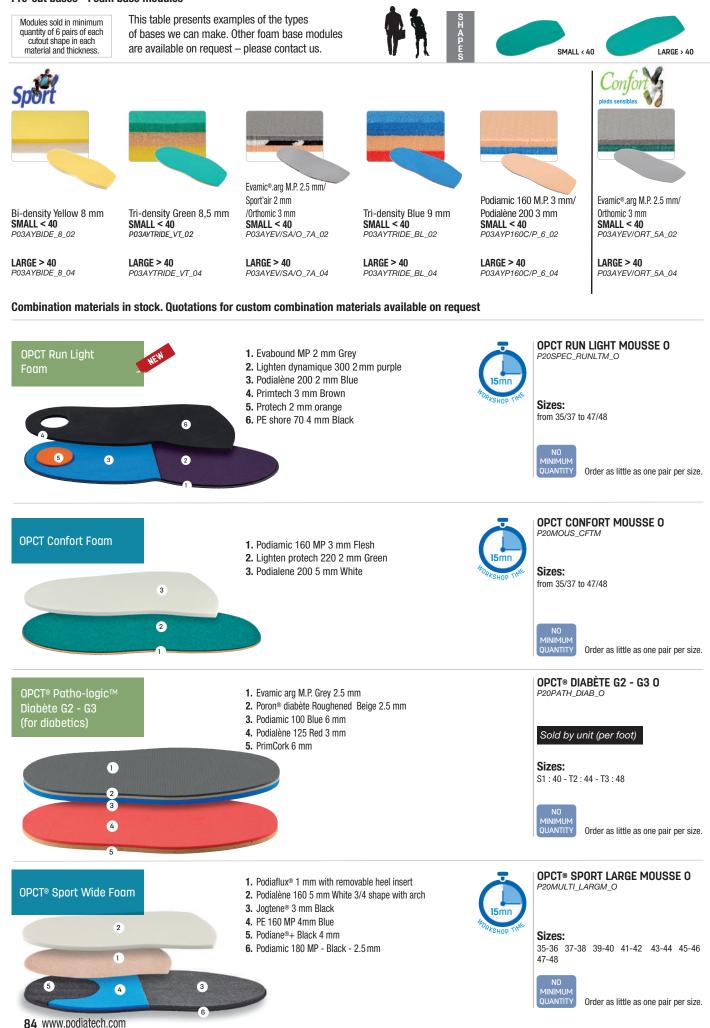

RESIN BASE MODULES

Archless bases Modules sold in min. quantity 6 pairs per cut-out shape and per material and thickness	Full insole		
REF. men's base, archless	BH	BF	BE kids
SIZES	39/40 - 41/42 - 43/44 - 45/46	35/36 - 37/38 - 39/40 - 41/42	25/26 - 27/28 - 29/30 - 31/32 - 33/34
PODIAFLUX® 1 MM	P03AAFLUX_1	P03ADFLUX_1	P03AGFLUX_1
PODIAFLUX [®] 1.2 M	P03AAFLUX_12	P03ADFLUX_12	P03AGFLUX_12
PODIAFLEX® 0.8 MM	P03AAFLEX_08	P03ADFLEX_08	P03AGFLEX_08
PODIAFLEX® 1.2 MM	P03AAFLEX_12	P03ADFLEX_12	P03AGFLEX_12
PODIAFLEX® 1.2 MM (47/48 - 49/50)	P03AAFLEX_12		
PODIAFLEX® 1.9 MM	P03AAFLEX_19	P03ADFLEX_19	P03AGFLEX_19
PODIAFLEX®/TISSU 1.3 MM	P03AATISS_13	P03ADTISS_13	
PODIAFLEX®/TISSU 1.7 MM	P03AATISS_17	PO3ADTISS_17	
Bases with arch	Full insole	3/4 insole	Preformed base with retrocapital p
Modules sold in min. quantity 6 pairs per cut-out shape and per material and thickness	S H A P E S		Primflex reinforcement 1.2 mm

in min. per cut and per thickne	quantity 6 pairs -out shape r material and sss	ŝ		reinforcement 1.2 mm
REF. me	en's / women's with arch	B HF V	B HF V 3/4	
SIZES		35/36 - 37/38 - 39/4	0 - 41/42 - 43/44 - 45/46	
	PODIAFLEX® 0.8 MM	P03AHFLEX_08	P03AIFLEX_08	
	PODIAFLEX [®] 1.2 MM	P03AHFLEX_12	P03AIFLEX_12	
	PODIAFLEX® 1.9 MM	P03AHFLEX_19	P03AIFLEX_19	Primflex base 1.2 mm
	PODIAFLEX®/TISSU 1.3 MM	P03AHTISS_13	P03AITISS_13	PRIMFLEX 01 PREF MOD Men's / women's with arch
	PODIAFLEX®/TISSU 1.7 MM	P03AHTISS_17	P03AITISS_17	35/36 - 37/38 - 39/40 - 41/42 - 43/44 - 45/46 PRIMFLEX 1.2 MM <i>P03PPRIMFL01</i>

PODIAFLEX® RESIN BASES + EMF MICROFIBRE TOP LAYER

DESCRIPTION	THICKNESS	SHOE SIZES (European)	REFERENCE No.	MATERIAL
MEN'S BASE, ARCHLESS	1.2 mm + 0.8 mm	39/40 - 41/42 - 43/44 - 45/46	PO3AAEMF/FLEX-12	PODIAFLEX [®] resin + black EMF*
LADIES' BASE ARCHLESS	0.8 mm + 0.8 mm	35/36 - 37/38 - 39/40 - 41/42	PO3ADEMF/FLEX-08	PODIAFLEX [®] resin + black EMF*
1.				

ORTHOTICS & ELEMENTS

LATEX SHEETS



Green latex foam Density = 270-280 sheet size 25 cm x 35 m



Salmon latex foam Density = 240-250sheet size 25 cm x 35 m 4 MM GREEN
P01LATEX145_4_VT

6 MM GREEN P01LATEX145_6_VT

8 MM GREEN
P01LATEX145_8_VT

10 MM GREEN P01LATEX145_10_VT 4 MM SALMON P01LATEX14S_4_SM

6 MM SALMON P01LATEX14S_6_SM

8 MM SALMON P01LATEX14S_8_SM

10 MM SALMON P01LATEX14S_10_SM

ORTHOTICS & ELEMENTS

PEVA ADHESIVE-COATED BASES

Full insole, adhesive-coated 1 mm, archless

PEVA BE, CHILDREN'S *P03AGPEVA* 25/26 - 27/28 - 29/30 - 31/32 33/34 - 35/36 - 37/38

P03AGPEVA



P03AAPEVA

PEVA BH, MEN'S *P03AAPEVA* 39/40 - 41/42 - 43/44 - 45/46



PEVA BF, LADIES' *P03ADPEVA* 35/36 - 37/38 - 39/40 - 41/42



P03ADPEVA

Option to add top layer to adhesive-coated PEVA bases. Contact us.



Procedure for making an orthotic without using glue.

Pull off the protective film from the area where the element is to be added







Position the element




Elements PRE-SHAPED LATEX ELEMENTS

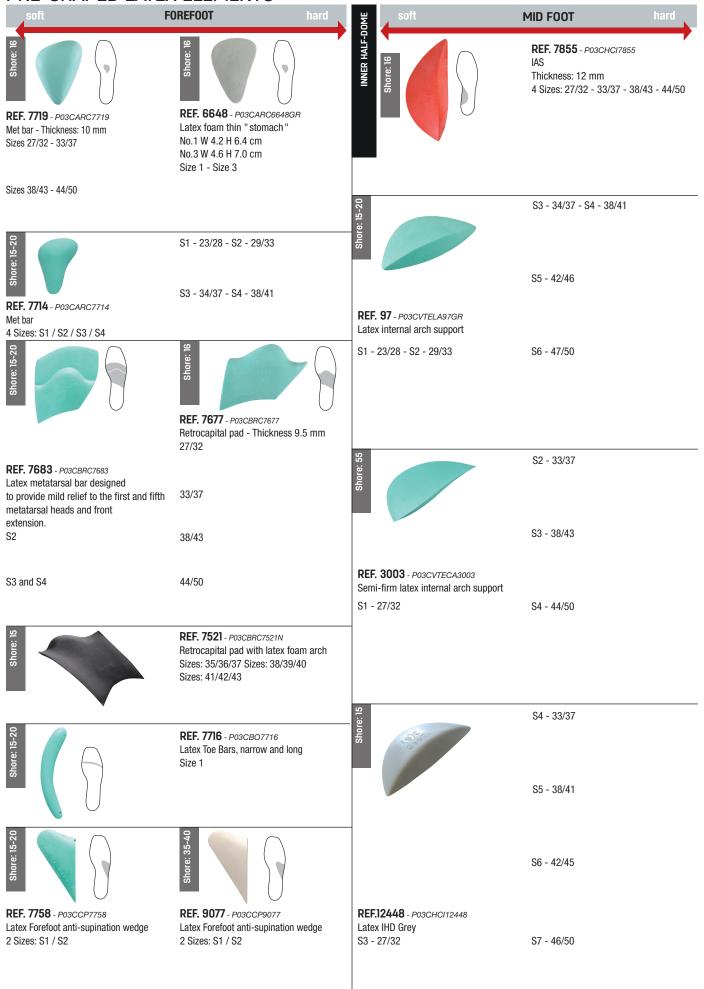

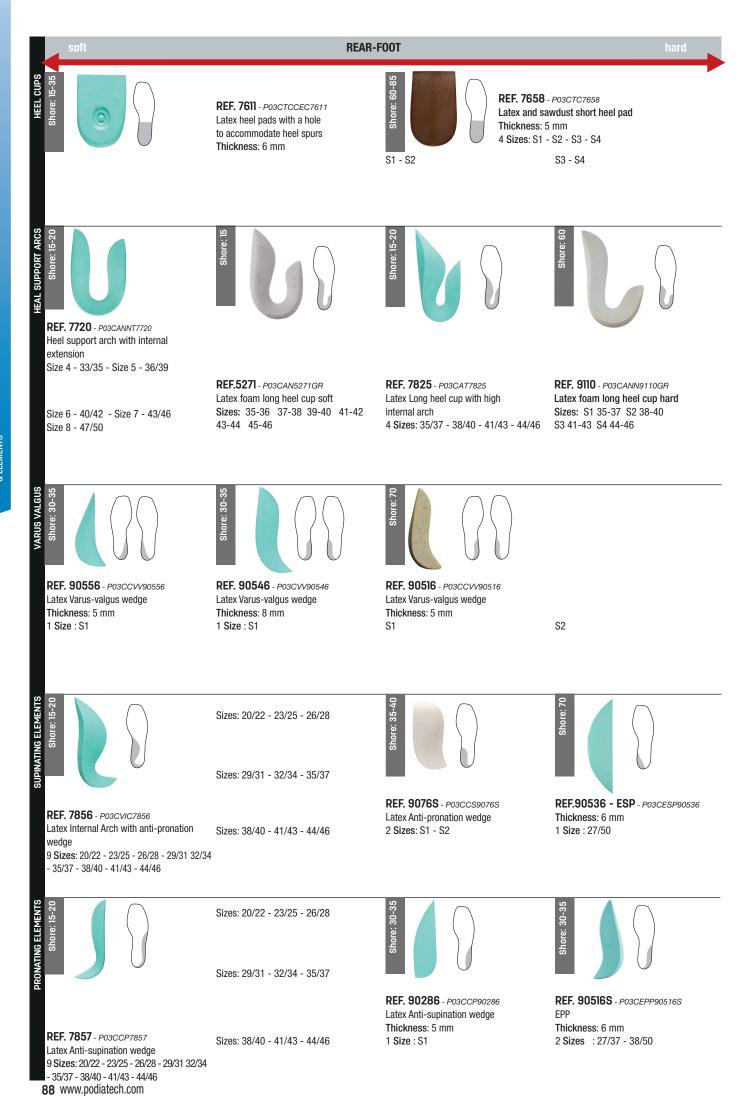

Elements

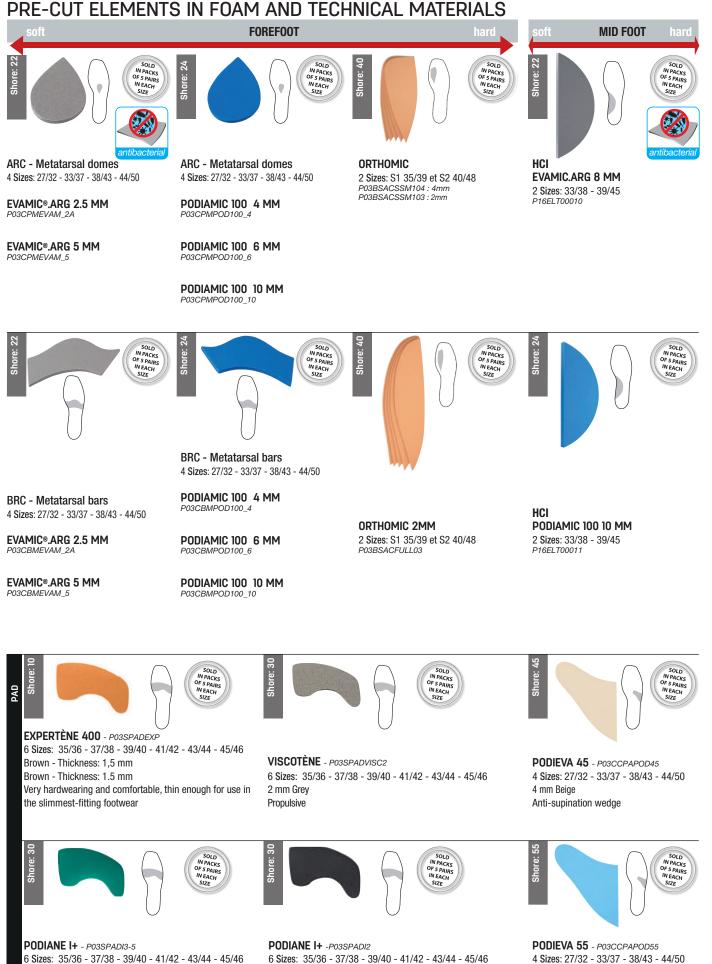

6 Sizes: 35/36 - 37/38 - 39/40 - 41/42 - 43/44 - 45/46 3.5 mm Green


6 Sizes: 35/36 - 37/38 - 39/40 - 41/42 - 43/44 - 45/46 2 mm Black Comfort for forefoot

www.podiatech.com 89

4 mm Light blue

Anterior anti-supination wedge

Elements PRE-CUT REAR-FOOT ELEMENTS

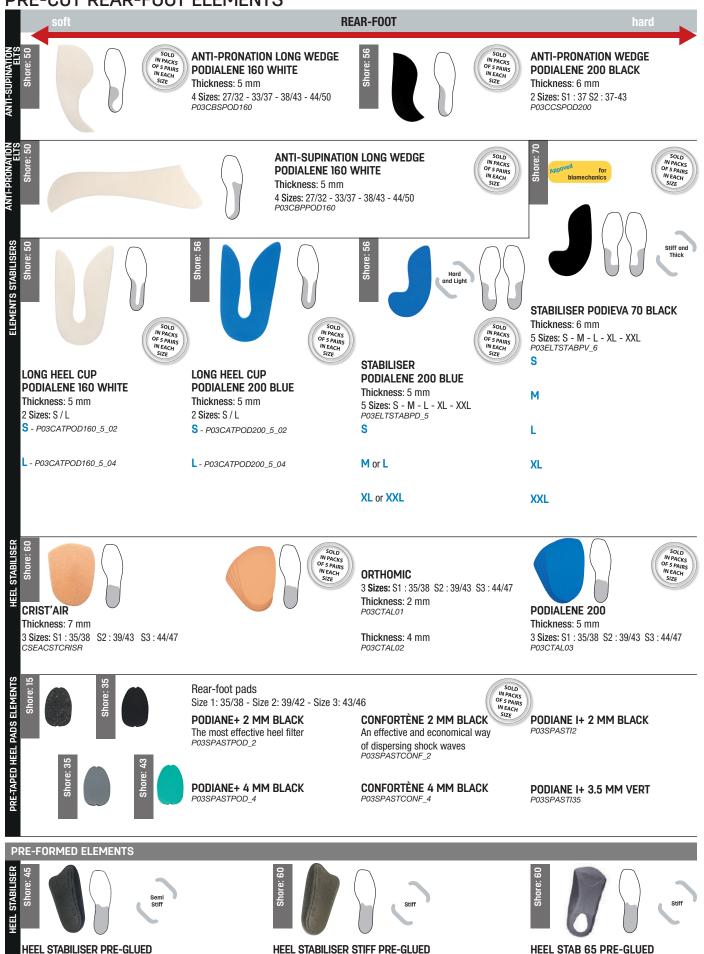

HEEL STABILISER STIFF PRE-GLUED

S1: 34/38 S2: 39/43 S3: 44/45 T4: 46/47

CSEACSTASTIF

HEEL STAB 65 PRE-GLUED 35/36 - 37/38 - 39/40 - 41/42 - 43/44 - 45/46 CSEACSTA65

90 www.podiatech.com

CSEACSTA

S1:35/38 S2:39/43 S3:44/47

Elements PRE-CUT RESIN ELEMENTS

Essential to reinforce Transflux® bases

Rear-foot LONG HEEL CUP 1 MM BLUE Sizes: XS S M L XL Sizes: 33-35 36-38 39-42 43-46 47-50 P03BTXR3/1

1.5 MM CHARCOAL P03BTXR3/1A



EXTERNAL 1 MM BLUE Sizes: XS S M L XL Sizes: 33-35 36-38 39-42 43-46 47-50 P03BTXR1/1

1.5 MM CHARCOAL P03BTXR1/1-5



INTERNAL 1 MM BLUE Sizes: XS S M L XL Sizes: 33-35 36-38 39-42 43-46 47-50 P03BTXR2/1

1.5 MM CHARCOAL P03BTXR2/1A



MULTI-APPLICATION TRANSFLUX® REINFORCEMENT 1 mm White pre-cut

Sizes: XS S M L XL P03BTXMR1/1



Adhesive-coated reinforcements Podiaflex® / Fabric 1.3 mm Sizes: 35-36 37-38 39-40 41-42 43-44 45-46 **EXTERNAL** P03SFTR1_13



INTERNAL P03SFTR2_13

LONG HEEL CUP P03SFTR3_13



SOLD

Anti-pronation reinforcement PODIAFIBER HARD™ 1 MM BLACK Sizes: 35-36 37-38 39-40 41-42 43-44 45-46 47-48

P03BANTPROPF_1_35-37-39-41-43-45-47



Forefoot



MIDFOOT Midfoot specific shape for increased internal and external support. without affecting the rear-foot volume. Sizes: 35-36 37-38 39-40 41-42 43-44 45-46 P03SFTR4_13



SOL

TRANSFLUX 1MM WHITE 3 Sizes: S1 35/39 - S2 40/44 - S3 45/48 P03BPAPILLON01

TRANSFLUX 1.5 MM CHARCOAL 3 Sizes: S1 35/39 - S2 40/44 - S3 45/48 P03BPAPILLON02

SOLD IN PACKS OF 5 PAIRS IN EACH SIZE



TRANSFLUX CARBONE 1.6 mm 2 Sizes: S1 35/39 et S2 40/48 P03BSACFULL01



PODIAFLEX CARBONE ` 2.2 mm Black 2 Sizes: S1 35/39 et S2 40/48 P03BSACFULL02



TRANSFLUX 1mm White 2 Sizes: S1 35/39 et S2 40/48 P03BSACSSM101 TRANSFLUX 1.5 mm Charcoal 2 Sizes: S1 35/39 et S2 40/48 P03BSACSSM102 SOLD IN PACKS OF 5 PAIRS IN EACH

SIZE

Elements PRE-CUT RESIN ELEMENTS



Composition of **"TESTERS" PACK** P81PCKTESTPOSTU

1 pair of each PODIAFLEX® 0.8mm ref. No. in all shoe sizes + 2 pairs of each PODIAFLEX® 1.2mm ref. No. in all shoe sizes

REFERENCE No.	THICKNESS	SHOE SIZES (European)	REFERENCE No.	THICKNESS	SHOE SIZES (European)
1 pair x EAI PODIAFLEX®	0.8 mm	S - M - L	2 pairs x EAI PODIAFLEX®	1.2 mm	S - M - L
1 pair x ECA PODIAFLEX®	0.8 mm	34/39 - 40/45	2 pairs x ECA PODIAFLEX®	1.2 mm	34/39 - 40/45
1 pair x BA PODIAFLEX®	0.8 mm	34/36 - 37/39 40/42 - 43/45	2 pairs x BA PODIAFLEX®	1.2 mm	34/36 - 37/39 - 40/42 - 43/45
1 pair x EPI PODIAFLEX®	0.8 mm	S - M - L	2 pairs x EPI PODIAFLEX®	1.2 mm	S - M - L
1 pair x EMI PODIAFLEX®	0.8 mm	S - M - L	2 pairs x EMI PODIAFLEX®	1.2 mm	S - M - L
1 pair x BPI PODIAFLEX®	0.8 mm	One size	2 pairs x BPI PODIAFLEX®	1.2 mm	One size
1 pair x BI PODIAFLEX®	0.8 mm	34/39 - 40/45	2 pairs x BI PODIAFLEX®	1.2 mm	34/39 - 40/45

Elements PRE-CUT RESIN ELEMENTS

Forefoot	Stimulators	Stimulators developed by	Mr and M	rs Villeneuve (registered designs)
REFERENCE No.	SHOE SIZES (European)	REFERENCE No.	THICKNES	SS SHOE SIZES (European)
P16ELT00130 - 08 BA	Anterior bar	PODIAFLEX [®] resin	0.8 mm	34/36 - 37/39 - 40/42 - 43/45
P16ELT00316NM_1 BA	Anterior bar	WHITEFLEX [®] resin	1 mm	34/36 - 37/39 - 40/42 - 43/45
P16ELT00130 - 12 BA	Anterior bar	PODIAFLEX [®] resin	1.2 mm	34/36 - 37/39 - 40/42 - 43/45
P16ELT00130 - 19 BA	Anterior bar	PODIAFLEX® resin	1.9 mm	34/36 - 37/39 - 40/42 - 43/45
P16ELT00134 - 08 EAI	Interior anterior element	PODIAFLEX [®] resin	0.8 mm	S - M - L
P16ELT00317NM EAI	Interior anterior element	WHITEFLEX [®] resin	1 mm	S - M - L
P16ELT00134 - 12 EAI	Interior anterior element	PODIAFLEX® resin	1.2 mm	S - M - L
P16ELT00134 - 19 EAI	Interior anterior element	PODIAFLEX [®] resin	1.9 mm	S - M - L

Midfoot	Stimulators	Stimulators developed by Mr and Mrs Villeneuve (registered design	IS)
REFERENCE No.	SHOE SIZES (European)	REFERENCE No. THICKNESS SHOE SIZES (European)	
P16ELT00131 - 08	EMI Interior medial element	PODIAFLEX [®] resin 0.8 mm S - M - L	
P16ELT00320NM	EMI Interior medial element	WHITEFLEX® resin 1 mm S - M - L	
P16ELT00131 - 12	EMI Interior medial element	PODIAFLEX® resin 1.2 mm S - M - L	
P16ELT00131 - 19	EMI Interior medial element	PODIAFLEX [®] resin 1.9 mm S - M - L	

Hindfoot		Stimulators	Stimulo	ntors developed by	Mr and M	rs Villeneuve (registered designs)
REFERENCE No.		SHOE SIZES (European)		REFERENCE No.	THICKNES	SS SHOE SIZES (European)
P16ELT00128 - 08	BPI	Interior posterior strip		PODIAFLEX [®] resin	0.8 mm	Unique
P16ELT00315NM_1	BPI	Interior posterior strip		WHITEFLEX® resin	1 mm	Unique
P16ELT00128 - 12	BPI	Interior posterior strip		PODIAFLEX® resin	1.2 mm	Unique
P16ELT00128 - 19	BPI	Interior posterior strip		PODIAFLEX [®] resin	1.9 mm	Unique
P16ELT00129 - 08	BI	Interior strip		PODIAFLEX [®] resin	0.8 mm	34/39 - 40/45
P16ELT00313NM	BI	Bande interne		WHITEFLEX® resin	1 mm	34/39 - 40/45
P16ELT00129 - 12	BI	Interior strip		PODIAFLEX [®] resin	1.2 mm	34/39 - 40/45
P16ELT00129 - 19	BI	Interior strip		PODIAFLEX [®] resin	1.9 mm	34/39 - 40/45
P16ELT00132 - 08	EPI	Interior posterior element		PODIAFLEX [®] resin	0.8 mm	S - M - L
P16ELT00319NM	EPI	Interior posterior element		WHITEFLEX® resin	1 mm	S - M - L
P16ELT00132 - 12	EPI	Interior posterior element		PODIAFLEX® resin	1.2 mm	S - M - L
P16ELT00132 - 19	EPI	Interior posterior element		PODIAFLEX [®] resin	1.9 mm	S - M - L
P16ELT00133 - 08	ECA	Anterior calcaneal element		PODIAFLEX [®] resin	0.8 mm	34/39 - 40/45
P16ELT00318NM	ECA	Anterior calcaneal element		WHITEFLEX® resin	1 mm	34/39 - 40/45
P16ELT00133 - 12	ECA	Anterior calcaneal element		PODIAFLEX [®] resin	1.2 mm	34/39 - 40/45
P16ELT00133 - 19	ECA	Anterior calcaneal element		PODIAFLEX [®] resin	1.9 mm	34/39 - 40/45

Sold in packs of 5 pairs

Sold in packs of 5 pairs

Sold in packs of 5 pairs

3Feet[®] PATHO-LOGIC[®]

The first orthopaedic insole that can be customised to incorporate high-quality SIDAS materials in just 5 minutes. Treat the main pathologies found in your patients by customising your insoles with forefoot, mid-foot and hindfoot elements.

ORTHOTIC-KITS[™] NEUTRE

Orthotic-kits[™] Neutre is your basic insole kit, ready for customising. It comprises a pair of 3Feet PATHO-LOGIC[®] insoles, a technical top layer, plus forefoot and hindfoot elements of the same density as the insoles.

The Orthotic-kits[™] Neutre is available in a High, Mid and Low version, in the following sizes: XS, S, M, L, XL.



ORTHOTIC-KITS[™] NEUTRE HIGH P24PATH NEU HI



ORTHOTIC-KITS[™] NEUTRE MID P24PATH NEU MI




ORTHOTIC-KITS[™] NEUTRE LOW P24PATH NEU LO

ADDITIONAL ELEMENTS

Our range of elements is designed to treat the main pathologies found in your patients and to satisfy customers who are looking for comfort and stability or to improve their sporting performance.

Depending on the needs of your patient, you can customise each insole with forefoot, mid-foot and hindfoot elements.

Each element is pre-glued and has its own easy-to-find location marked out on the insole.

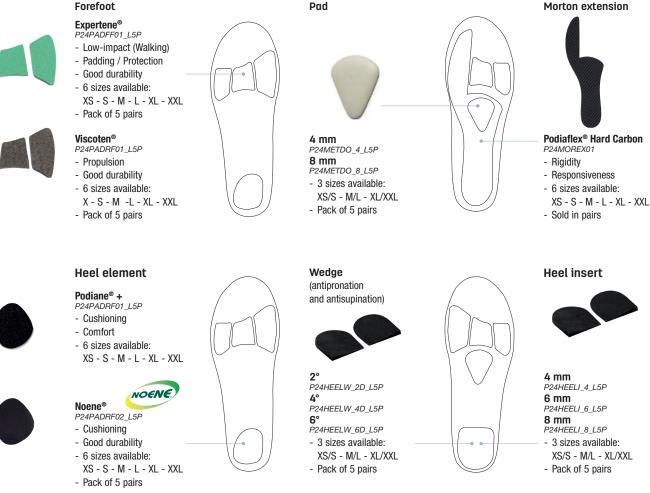





94 www.podiatech.com



Insoles specially designed for high-arched feet. 3Feet[®] Slim insoles are suitable for slim-fitting footwear. They can also be cut for use as half insoles.





Can also be cut for use as half insoles.

3FEET® SLIM HIGH

CSE3FSLIM13_HI

- 1. Microfibre top layer: comfort & antibacterial
- 2. EVA shell with a high arch: foot support & pressure distribution
- 3. Perforations: breathability
- 4. 2mm EVA: thinness & cushioning
- 5. Cutting line to make a half insole

DIMENSIONS (mm)

5



Insoles specially designed for medium-arched feet. 3Feet[®] Slim insoles are suitable for slim-fitting footwear. They can also be cut for use as half insoles.






3FEET® SLIM MID

CSE3FSLIM13_MI

- Sizes: XS-S-M-L-XL-XXL
- 1. Microfibre top layer: comfort & antibacterial
- 2. EVA shell with a medium arch: foot support
- 3. Perforations: breathability
- 4. 2mm EVA: thinness & cushioning
- 5. Cutting line to make a half insole

DIMENSIONS (mm)

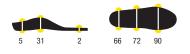

Insoles specially designed for low-arched feet. 3Feet[®] Slim insoles are suitable for slim-fitting footwear. They can also be cut for use as half insoles.






3FEET® SLIM LOW

CSE3FSLIM13_LO

- Sizes: XS-S-M-L-XL-XXL

- 1. Microfibre top layer: comfort & antibacterial
- $\ensuremath{\textbf{2.}}$ EVA shell with a low arch: foot support & stability
- 3. Perforations: breathability
- 4. 2mm EVA: thinness & cushioning

DIMENSIONS (mm)

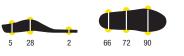





DEMO KIT 3FEET® SLIM

DEMO KIT 3FEET[®] SLIM HIGH CSE3FDKSLIHI13

DEMO KIT 3FEET® SLIM MID CSE3FDKSLIMI13

DEMO KIT 3FEET[®] SLIM LOW CSE3FDKSLILO13



FOOT & SHOES ACCESSORIES

CREAMS

Sidas has begun a new venture, launching a brand-new range of sports creams. These creams containing superior-quality natural ingredients are designed for outdoor sports including running, hiking and trail running.

Anti-rubbing and anti-blister cream. Protects, moisturises and revitalises.

This cream has a protecting action, forming a film that boosts the skin's elasticity and defensive ability.

Apply before sport by rubbing it in at length to areas prone to friction: feet, thighs, arm pits, nipples, etc.



ANTI-FRICTION CREAM 75 ML

CFCCREANTIF17_75ML

- pack of 15 or singly



COMPOSITION

- Shea butter, marine collagen, vitamin E, witch hazel extracts
- Does not contain parabens or PEG, contains natural preservatives. Not tested on animals.
- 99% natural ingredients

3 • •

Provides rapid relief from injury: strained muscles, impact injuries, cramp, dislocations and bruising.

The silicon and targeted essential oils act in synergy to produce a three-fold effect: relaxing, stress-relieving and repairing. Apply locally as soon as possible after injury. Repeat application as necessary.



RECOVERY CRYO GEL 75 ML

CFCCRERECOV17_75ML

- pack of 15 or singly




COMPOSITION

- Organic silicon and essential oils (mint, peppermint, eucalyptus, rosemary, tea tree)
- 77% natural ingredients

FOOT PROTECTIONS

These gel heel pads absorb shocks all day long for optimum comfort. They are easy to fit and will quickly become an essential part of your daily life.

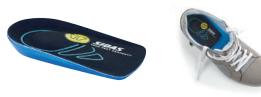

These gel bone spur pads relieve pain affecting the heels specifically. They also absorb shocks all day long for optimum comfort. They are easy to fit and will quickly become an essential part of your daily life.





These wrap-around gel heel cups absorb shocks all day long for optimum comfort. Specially shaped to fit round the heel, they prevent blisters and stabilise your feet, enabling you to enjoy your day or sporting activity to the full.



MARKETING TOOLS

IN-SHOE PADS DISPLAY

CFCFPDISPLAY18

Double-sided display unit designed specifically for the In-Shoe range

- French and English
- Dimensions: 38 cm wide x 125 cm high
- 12 display prongs on each side (36 in all)Up to 5 products per display prong, making
- 120 products in total
- 1 document holder and 1 rivet for attaching a sample product



IN-SHOE PADS PLEXI POS

FR / EN - CFCFPPOSPLEXI12_01 EN / DE - CFCFPPOSPLEXI18_02

Plexiglas® sales aid incorporating the following product: 1 Gel Heel Pad.

- Sales aid enabling consumers to touch a gel product.
- Can be placed on any flat surface.
- French and English
- Dimensions: 20 cm wide x 17 cm high



GEL HEEL PADS

CFCFPHEELPAD18

- S/M (36-41) ; L/XL (42-46)
- Pack: 1 pair
- Supplied in boxes of 10 only

FITTING

Insert the heel cups into the back of your footwear. Clean your heel cups regularly with soap and water to ensure they remain in place.

GEL BONE SPUR PADS

CFCFPBONESPU18

- S/M (36-41) ; L/XL (42-46)
- Pack: 1 pair
- Supplied in boxes of 10 only

FITTING

Insert the heel cups into the back of your footwear. Clean your heel cups regularly with soap and water to ensure they remain in place.

GEL HEEL CUPS

CFCFPHEELCUP18

- S/M (36-41) ; L/XL (42-46)
- Pack: 1 pair
- Supplied in boxes of 10 only

FITTING

Insert the heel cups into the back of your footwear. Clean your heel cups regularly with soap and water to ensure they remain in place.

IN-SHOE PADS WALL POS

EN/ FR - CFCFPWALLPOS18

Hanging display for the In-Shoe Pads range.

Available separately or to complement the main display. - French and English

- Dimensions: 20 cm wide x 87 cm high
- 5 display prongs per column, accommodates a total of 10 products



FOOT PROTECTIONS

Sidas markets these protective devices to ease pain arising from different sporting activities. These products are applied directly onto the skin to protect the wearer from blisters, runner's toenail and chafing.

These gel toe caps enable you to enjoy your sporting activity to the full. Used preventively, they prevent the onset of runner's toenail caused by shocks inside the footwear. If pain is already present, they provide immediate relief by absorbing shocks.

Protective caps made of silicone, washable and reusable.



GEL TOE CAPS

- Sizes: XS (small toes) / S/M (36-40) / L/
- Pack: 1 pair

ASKINGLTOECA17

- Adjustable length
- Available singly or in a box of 20

Sizes Tailles	Big Toe Gros orteil	Middle Toes Orteils moyens	Small Toes Petits orteils
36 - 39 (EU) 4 - 6 (US) 3 - 6 (UK)	S/M	XS	XS
40 - 46 (EU) 7.5 - 13 (US) 6.5 - 12 (UK)	L/XL	S/M	XS

These gel toe wraps enable you to enjoy your sporting activity to the full. They relieve pain caused by the toes rubbing together and prevent corns and calluses.



GEL TOE WRAPS

ASKINGLTOEWR17

- Sizes: S/M (36-40) / L/XL (41-47)
- Pack: Pack of 4
- Adjustable length
- Available singly or in a box of 20

These gel plantar protectors enable you to enjoy your sporting activity to the full. Used preventively, they prevent chafing, friction burns and blisters. If pain is already present, they provide immediate relief by absorbing shocks.




GEL PLANTAR PROTECTORS

ASKINGLPLAPR17

- Sizes: S/M (36-40) / L/XL (41-47)
- Pack: 1 pair
- Anatomically constructed for the right and left feet
- Available singly or in a box of 20

Applying these Silitène[™]-T gel protectors to the toes as soon as rubbing is felt prevents the onset of blisters and relieves pain.

Their biomechanical construction makes them light and almost imperceptible and they will remain in place for a very long time, enabling you to enjoy your activity to the full.

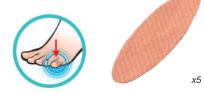

TOE PROTECTORS V2 CFCOSTOEPROT20

- Thickness: 1 mm.

CFCOSFOOTPROT20

- Thickness: 1 mm

- Pack contents: 5 toe protectors

FOOT PROTECTOR V2


- Pack contents: 5 foot protectors

Applying this Silitène[™]-T gel foot protector to the area where chafing occurs as soon as the first rubbing is felt prevents the onset of blisters and relieves pain.

Its biomechanical construction makes this foot protector light and almost imperceptible. It moves with the foot and will remain in place for a very long time, enabling you to enjoy your activity to the full.

Silitene[™]-T protection kit, adapted to the anatomy of the toe and heel. The prevention kit Silitene[™]-T


Their biomechanical construction follows the movement of the foot and allows a very long wear.

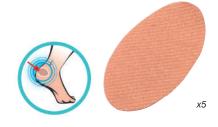

ESSENTIAL PROTECTOR KIT V2

CFCOSESSPROT20

- Thickness: 1 mm
- Pack contents: 5 x toe protectors / 5 x foot protectors

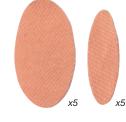

Sheet of Silitène™ to be cut and used as required.

is the essential partner for your sporting activity.

Apply the pads to friction areas for immediate relief and comfort.

Apply to the area where rubbing occurs as soon as the first hotness is felt. You will forget you are wearing this light and almost imperceptible Silitène[™] sheet, enabling you to enjoy your activity to the full.



FOOT PROTECTOR SHEET

FOOT PROTECTOR SHEET 1 MM ASKINFOOTPROT14_01SH

FOOT PROTECTOR SHEET 2 MM ASKINFOOTPROT14_02SH

- Dimensions: 90x100mm (1mm) / 100x115mm (2mm)

MARKETING TOOLS

ON SKIN BOX

Product box can also be used as a display unit. All On Skin products are supplied in 20s in this box.





CEDAR WOOD DRYER BAG

CSHACDRYBAWOOD19

This cedar wood dryer bag is the natural way to dry out all types of damp footwear.

- Absorbs the humidity from footwear four times more quickly than air drying.
- Antifungal action removes unpleasant odours.
- Made of 100% natural cedar wood and a polyester fabric.
- Re-usable, hardwearing and long-lasting product.
- Supplied in boxes of 10. Box can be used as a display unit.

www.podiatech.com 99

FOOT AND SHOES SPECIALIST

PU RESINS



MOULD RELEASE AGENTS AND PU ACCESSORIES



Accessories: PRESTOFORM® ULTIMATE



Accessories

Magnetic orthotic stand PORSOCLEP_S

Top piece for orthotic stand PORCOQUILLE

Articulated seat + M12 pin for orthotic stand PORPATINMOBI

Set of 6 pins for orthotic stand (Ø 7, 8, 9, 10, 11 and 12 mm) $\ensuremath{\textit{PORTIGEP}_S}$

Pack of 20 orthotic stand rectangles PORLOT20R

Pack of 20 shaft circles PORLOT20C

Cooling pouch size S PORPOCHREF10

Cooling pouch size S PORPOCHREF30

Strap for cooling pouch PORSANGLPOCH

Consumable accessories

Latex membrane 10/10° 600 x 700 mm P99MLATORT1

Latex membrane 20/10° 600 x 700 mm P99MLATORT2

Large plate membrane 10/10° 800 x 1000 mm P99MLATGPL1

Large plate membrane 20/10° 800 x 1000 mm P99MLATGPL2

Parasuper membrane 1 mm 600 x 700 mm P99MPARA1MM

Parasuper membrane 2 mm 600 x 700 mm P99MPREST3

Ultimate membrane 3 mm 600 x 700 mm P99MULTIMATE_30

Consumables: PRESTOFORM® ULTIMATE



Podialène 160 white 5 mm or 10 mm



Podialène 160 sand 5 mm or 10 mm



Podialène 200 blue 5 mm or 10 mm

Other custom combination materials can be manufactured on request, contact us.

Second layer module (1 piece) Sizes S and L

First layer module

(1 piece) Sizes S and L

Surrounds (AEPs)

(1 piece) Sizes S and L



Podialène 160 white 5 mm, 10 mm or 15 mm

Podialène 160 sand

5 mm, 10 mm, 15 mm or 20 mm





AEP Podialène AEP 160 White 5 mm S P15APOD1605_BC_02

AEP Podialène AEP 160 White 5 mm L P15APOD1605_BC_04

AEP Podialène AEP 160 White 10 mm S P15APOD16010_BC_02

AEP Podialène AEP 160 White 10 mm L P15APOD16010_BC_04

AEP Podialène AEP 160 Sand 5 mm S P15APOD1605_SB_02

AEP Podialène AEP 160 Sand 5 mm L P15APOD1605_SB_04

AEP Podialène AEP 160 Sand 10 mm S P15APOD16010_SB_02

AEP Podialène AEP 160 Sand 10 mm L P15APOD16010_SB_04

AEP Podialène AEP 200 Blue 5 mm S P15APOD2005_BL_02

AEP Podialène AEP 200 Blue 5 mm L P15APOD2005_BL_04

AEP Podialène AEP 200 Blue 10 mm S P15APOD20010_BL_02

AEP Podialène AEP 200 Blue 10 mm L P15APOD20010_BL_04

Bi-densité 8 mm S - left or right foot P15M1CBIDE8

Bi-densité 8 mm L - left or right foot P15M1CBIDE8

Bi-densité 13 mm S - left or right foot *P15M1CBIDE13*

Bi-densité 13 mm L - left or right foot P15M1CBIDE13

Podialène 160 White 5 mm S P15M1CP1605_BC_02

Podialène 160 White 5 mm L P15M1CP1605_BC_04

Podialène 160 White 10 mm S P15M1CP16010_BC_02

Podialène 160 White 10 mm L P15M1CP16010_BC_04

Podialène 160 Sand 10 mm S P15M1CP16010_SB_02

Podialène 160 Sand 10 mm L P15M1CP16010_SB_04

Podialène AEP 160 White 5 mm S P15M2CP1605_BC_02

Podialène AEP 160 White 5 mm L P15M2CP1605 BC 04

Podialène AEP 160 White 10 mm S P15M2CP16010_BC_02

Podialène AEP 160 White 10 mm L P15M2CP16010_BC_04

Podialène AEP 160 White 15 mm S P15M2CP16015_BC_02

Podialène AEP 160 White 15 mm L P15M2CP16015_BC_04

Podialène AEP 160 Sand 5 mm S P15M2CP1605 SB 02

Podialène AEP 160 Sand 5 mm L P15M2CP1605 SB 04

Podialène AEP 160 Sand 10 mm S P15M2CP16010_SB_02

Podialène AEP 160 Sand 10 mm L P15M2CP16010_SB_04

Podialène AEP 160 Sand 15 mm S P15M2CP16015_SB_02

Podialène AEP 160 Sand 15 mm L P15M2CP16015_SB_04

Podialène AEP 160 Sand 20 mm S P15M2CP16020_SB_02

Podialène AEP 160 Sand 20 mm L P15M1CP16020_SB_04

MATERIALS FOR FITTING AND REPRODUCTION - EXTRUDED FILMS







ORTHOCONTROLE PEVA 0.9 MM - *POGCTRLVA* **Fitting, soft enough for walking.**

- Sheet dimensions: 500 mm x 600 mm

- Orthocontrôle PEVA, 0.9 mm
- Kit of 50 sheets

EXTRUDED FILMS - Thickness: 0.9 mm

EXTRUDED FILMS - Thickness: 0.4 mm

- Sheet dimensions: 500 mm x 600 mm

ORTHOCONTROLE PVC - POGCTRLVC Fitting, semi-rigid, good transparency. - PVC, 0.4 mm

- Orthocontrôle PVC, 0.4 mm
- Kit of 50 sheets
- KIL UI JU SHEELS

ORTHOLAST SHEETS Reproduction for heavy use. - Thickness: 2 mm and 4 mm

ORTHOLAST KIT, 2 MM - PO6KITLAST2 - Kit of 20 sheets

- Sheet dimensions: 500 mm x 600 mm

ORTHOLAST KIT, 4 MM - PO6KITLAST - Kit of 10 sheets

		Advantageo		
	TEPS LYCRA Grey	Advantages: - Polyester fabric stretchable in every direction.		
	Pre-glued	- Very hardwearing.	Diabetes	
S	- Sheet dimensions: 1500 mm x 1000 mm	 Will fit on even the most difficult shapes. Easier to slide foot in footweart. 	DIABETES: Mate	erial recommended
	 Thickness: 0.6 mm Colour: Grey 	- Pleasant to the touch, does not induce	for the feet of d	iabetics.
		overheating or have a characteristic smell. - Washable.		
		Precautions for use:		
		- Apply a film of neoprene cement to the other material		
		when adding this material (avoid thick glues). Remove the protective film on the Lycra and re-activate it		
		either in an oven or with a heat gun. Then place on	P01BTEPSLYCR	
		top of the other material.		
	TEPS CHAMOISE	Advantages: - Extensible and very comfortable.		
	- Sheet dimensions:	- Can be used to cover complex surfaces when making		
V	1500 mm x 1000 mm - Thickness: 0.9 mm	complex devices. - Well tolerated by the skin.		
S	- Colour: Desert	- Wein tolerated by the skin. - Washable.		
		Technical applications and observations:		
		 Covering for plantar orthotics for elderly patients, arthritis, rheumatoid foot, diabetic foot (first and 		
		second degree).		
		- Covering for complex sockets and orthotics; shoe and		
		slipper lining; - Interior covering for special appliances.		
		- This material should be turned over at the edges to	P01BTEPSCHAN	1
		prevent delamination.		
	LIGHTEN CONFORT 150	Light, durable and comfortable. For use in many		
	- Sheet dimensions:	applications.		
V	490 mm x 550 mm - Thickness: 7 mm			
S	- Colour: Blue			
		25 150 👗		
		Passive padding	P01HLIGHTENC	_7_BL
	Ring uppers for Se	For use with a slow forward flex movement		
MANUFACIU				
	DIATEPS	Three-layer material comprising: - 1 layer TEPS sueded		A AT
あんたたまましい	 Sheet dimensions: 1400 mm x 950 mm 	- 1 layer Podiamic 180 MP black, 2.5 mm		The second second
	- Thickness: $0.5 \text{ mm} + 2.5 \text{ mm}$	- 1 layer Extensouple carbone		State State
	+ 0.9 mm = 3.9 mm	Technical applications and observations: combination material that becomes thermo-stretchable		
TETRA MARTINETERS AND	- Colour: Black / Black / Desert	when subjected to 100°C heat in an oven.		
		Manufacturing uppers for sensitive feet.	P02DIATEPS4	
EXTERIOR CO	VERING FOR UPPERS	S AND ORTHOTICS		
Contraction States Section 5	EXTENSOUPLE CARBONE	Advantages:		
	- Sheet dimensions:	Fabric extensible in all directions.Carbon-appearance finish.		
	1400 mm x 1000 mm	- Will fit on even the most complex shapes.		
S CONTRACTOR CONTRACTOR	 Thickness: 0.5 mm Colour: Black 	Technical applications and observations:		
		- Exterior covering for uppers and orthotics.		-
and the second start of the			P06EXTENSCAR	R
Please contact us it	f you wish to order CSL FIN, Al	NF MP or SPC MP in rolls (large sizes)		
TOE-ENDS				
BOUTFLEX V25P PR	E-CUT TEMPLATES - P06BDURF	BOUTFLEX V45P PRE-CUT TEMPLATES - POG	BDURH	Advantages :
B				- Very malleable when hot
J.		V		- easy to heat shape
	- For a thinner upper	- For a firmer upper		Heating temperature=
	- 0.8 mm Beige - 1 pair	- 1 mm Beige - 1 pair		$\bigcap_{i=1}^{n}$
•	н.».			C° T°. 100° C

- For a thinner upper
 0.8 mm Beige
 1 pair

www.podiatech.com 105



PATIENT CHAIRS





PODIATECH, exclusive partner of BTC




1


ъ

11

1

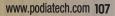

CONFIGURE YOUR CHAIR With PODIATECH you can order a custom

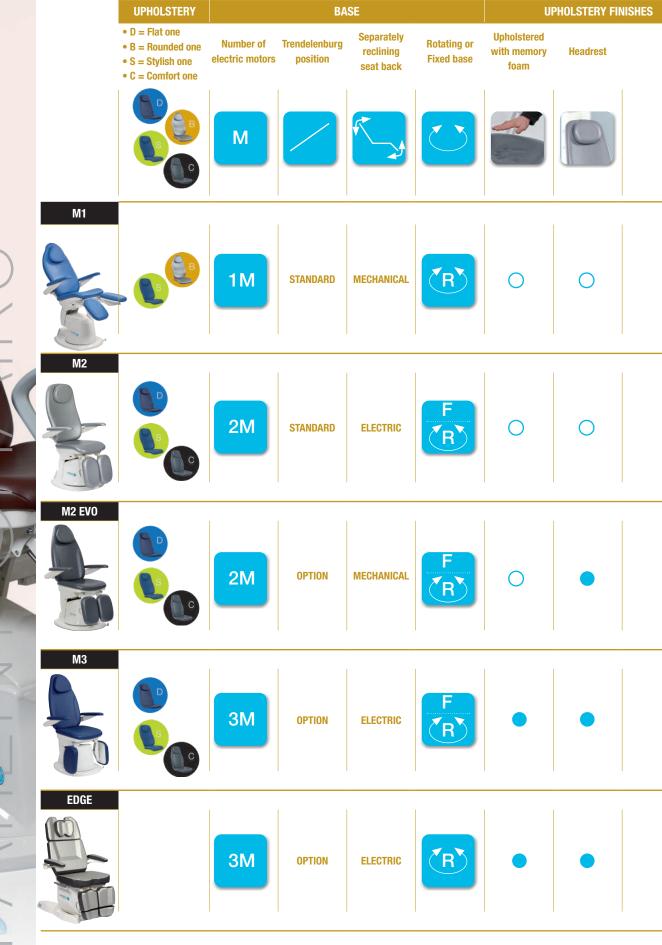


CHAIRS AND SEATS

PODIATECH



CHAIRS AND SEATS



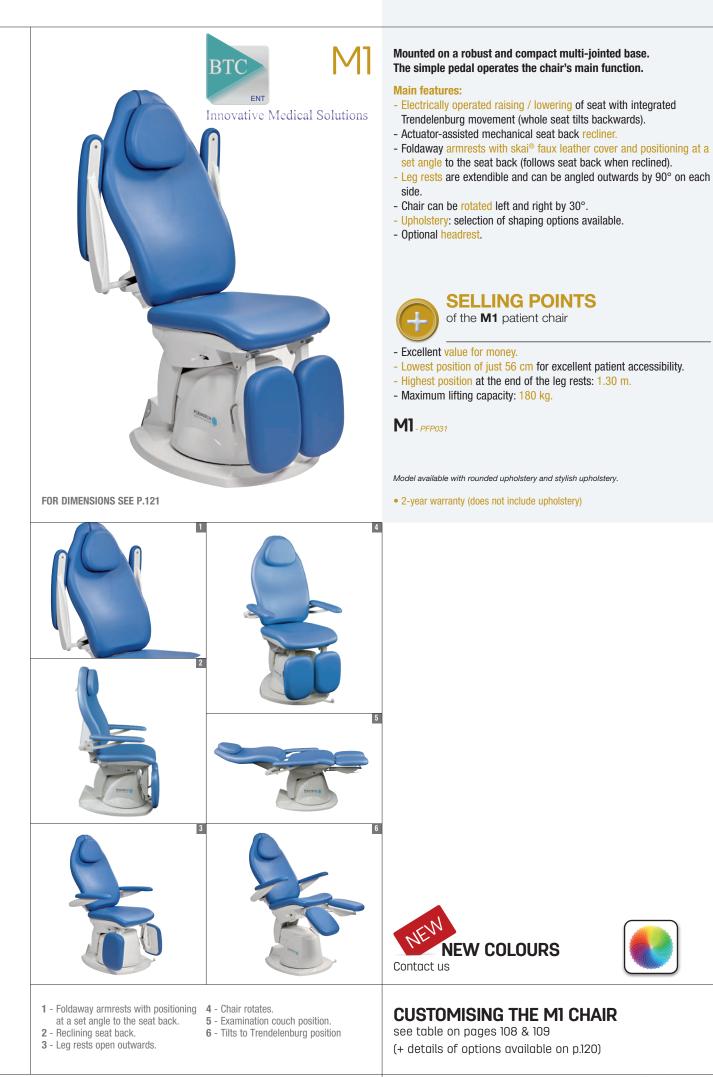

CHAIRS AND SEATS







Innovative Medical Solutions

Mounted on a robust and compact multi-jointed base. The main movement functions are controlled by a simple pressbutton pedal.

Main features:

- Maximum height: 130 cm.
- Minimum height of seat: 56 cm.
- Electrically operated raising / lowering of seat with integrated Trendelenburg movement (whole seat tilts backwards).
- Actuator-assisted mechanical seat back recliner.
- Foldaway armrests with skai® faux leather cover and positioning at a set angle to the seat back (follows seat back when reclined).
- Leg rests are extendible and can be angled outwards by 90° on each side.
- Chair can be rotated left and right by 30°.
- Fully reclining.
- Automatic chair reset to the starting position.
- Optional headrest.
- Upholstery: selection of shaping options available.

M2 - PFP031

Model available with rounded upholstery and stylish upholstery.

• 2-year warranty (does not include upholstery)

M2 - Trendelenburg position

FOR DIMENSIONS SEE P.121




- 2 Leg rest protector kit (one pair).
- 3 Memory foam.
- 5 Kit: fastenings + collection trays (one pair
- 6 Leg rests open outwards.

CUSTOMISING THE M2 CHAIR

see table on pages 108 & 109 (+ details of options available on p.120)



CHAIRS AND SI

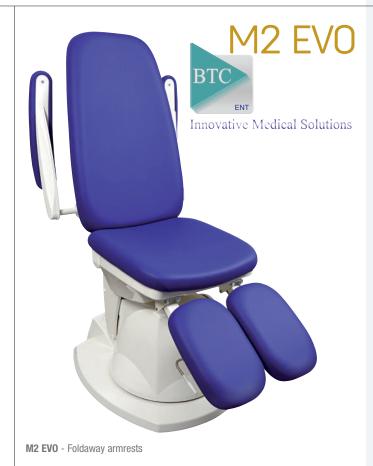

Its new shape and recent kinematic upgrade make this a highly sought-after chair on the market in terms of functionality and elegance.

Main features:

- Maximum height: 145 cm.
- Minimum height of seat: 56 cm.
- Electrically operated seat raising / lowering.
- Independently electrically operated backwards tilting of the whole chair (Trendelenburg position).
- Actuator-assisted mechanical seat back recliner.
- Foldaway armrests with skai® faux leather cover and positioning at a
- set angle to the seat back (follows seat back when reclined). Headrest.
- Leg rests are extendible and can be angled outwards by 90° on each side.
- Chair can be rotated left and right by 115°, or fixed-base option.
- Adjustable couch position.
- Automatic chair reset to the starting position.
- Upholstery: selection of shaping options available.

M2 EVO - PFP032 + PFP026

Price includes fixed base, skai® faux leather armrests and back handles.

FOR DIMENSIONS SEE P.121

• 2-year warranty (does not include upholstery)

• 3-year warranty for raising/lowering actuator



1 - Starting position.

2 - Maximum amplitude.

- 3 Memory foam.
- 4 High couch position.
 5 Kit: fastenings + collection trays (one pair).
 6 Leg rests open outwards.





CUSTOMISING THE M2 EVO CHAIR see table on pages 108 & 109

(+ details of options available on p.120)









CHAIRS AND SEATS

PATIENT CHAIR OPTIONS AND ACCESSORIES





LEG REST BAG HOLDER - PFP045



LEG REST PROTECTOR KIT (ONE PAIR) - PFP007 Washable.



FLASH GERM' MULTIPURPOSE **SPRAY 750 ML** - *P71D035*



KIT: FASTENINGS + COLLECTION TRAYS (ONE PAIR) - PFP005

Includes fastening mechanisms and pair of removable collection trays.



FLEXIBLE COLLECTION TRAY

PFP041

MEMORY FOAM - PFP018

In addition to the HR (high-resilience) foam on the chair, we offer a combination material featuring MEMOREX memory foam for enhanced patient well-being.




M1 & M2 HEADREST KIT - PFP029

ROLL HOLDER - PFP027



KIT OF 2 FOAM INSERTS (EXAMINATION COUCH POSITION) PFP028

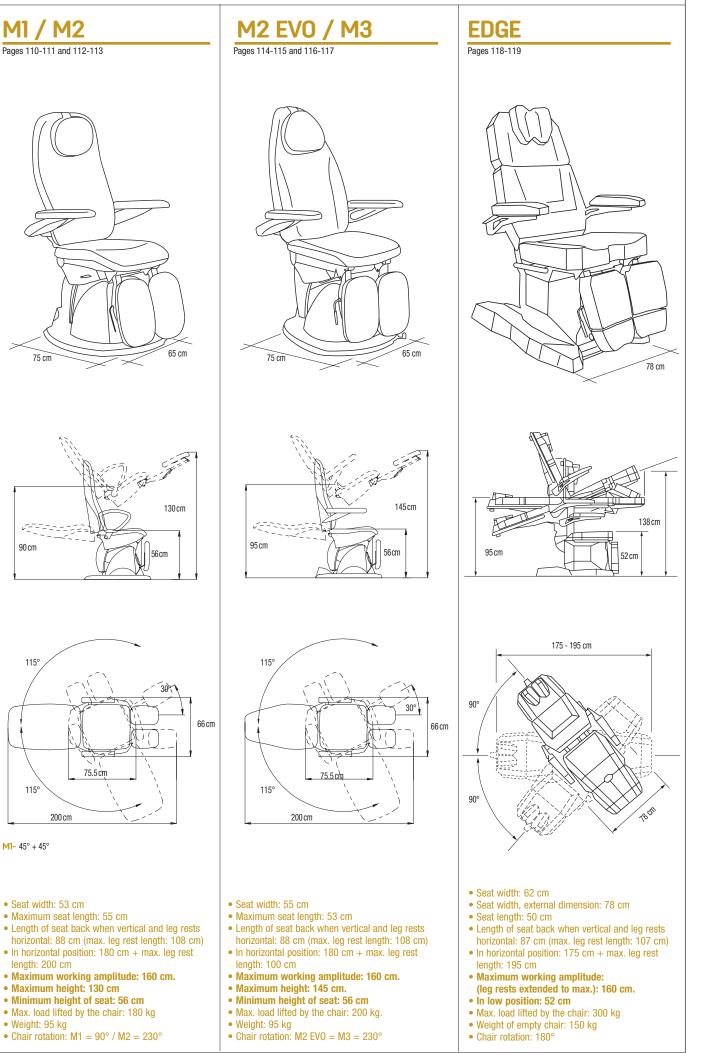



• 1-year warranty (does not include upholstery)



(does not include upholstery)



• 1-year warranty (does not include upholstery)



(does not include upholstery)

S2 OPERATOR SEAT

Practitioner chair with saddle-type seat.

Main features:

- Seat tilts forwards and backwards, with locking mechanism.
- Rubber casters suitable for all floor types.
- Footprint: 65 cm x 65 m.

Seat height adjustment:

- Low position: 53 cm.
- High position: 66 cm.
- Weight : 10 kg

S2 - PSIE002

S2 EVO OPERATOR SEAT

Practitioner chair with saddle-type seat and ergonomic back.

Main features:

- Seat tilts forwards and backwards, with locking mechanism.
- Rubber casters suitable for all floor types.
- Footprint: 65 cm x 65 cm.

Seat height adjustment:

- Low position: 56 cm.
- High position: 77 cm.
- Weight : 10 kg

S2 EVO - PSIE009

S3 OPERATOR SEAT

Practitioner chair with circular memory foam seat.

Main features:

- Seat back is height-adjustable and can be positioned to the left or right.
- Rubber casters suitable for all floor types.
- Footprint: 65 cm x 65 cm.

Back rest height adjustment:

- Low position: 83 cm.
- High position: 90 cm..

Seat height adjustment:

- Low position: 53 cm.
- High position: 66 cm.
- Weight : 15 kg

S3 - PSIE001

S4 EVO OPERATOR SEAT

Practitioner chair with memory foam seat and ergonomic seat back

Main features:

- Seat tilts forwards and backwards, with locking mechanism.
- Rubber casters suitable for all floor types.
- Footprint: 65 cm x 65 cm.

Seat height adjustment:

- Low position: 50 cm.
- High position: 64 cm.
- Weight : 20 kg

S4 EVO - PSIE010

NEW

COLOURS

Contact us






TREATMENT UNITS



CHOOSE YOUR TREATMENT UNIT

		ERGONOMIC OPERATION		TREATMENT Technique						1	FECHNICAL	L FEATURES
		Retractable cords	Hanging cords	Dry (motor + suction)	Spray (motor + spray)	Control pedal	ON/OFF Pedal	Instrument detection	Induction micromotor	Carbon brush micromotor	Integrated suction with adapted filter and tube	Suction pump in cabinet with adapted tube
	ACCESS											
			•									
)	UPLAY											
	CLASSIC						•		-		•	
7	CLASSIC		•			0	•			Dry	Dry	
-	COLIBRI		-			-	-		Spray	Spray	Spray	Spray
			•			0	•		Dry O Spray	Dry Spray	Dry Spray	Dry O Spray
- NEL	SIROCCO									0	0	
		•	•	•		0	•	•	Dry O Spray	Dry Spray	Dry Spray	Dry O Spray
	ÉLÉGANCE								Dry	Dry	Dry	Dry
			•			0			Spray	Spray	Spray	Spray
			•	•		•						
		ON KIT							Dry	Dry	Dry	Dry
			•			0			Spray	Spray	Spray	Spray
	KART UK					\cap			Dry	Dry	Dry	Dry
			Optional	mini a		Choice of			Spray	Spray	Spray	Spray
intach agr	(1) Classic or Pi	remium. (2)	Optional	mm-comp	pressor. (3) (unoice of	wo colou	15.				

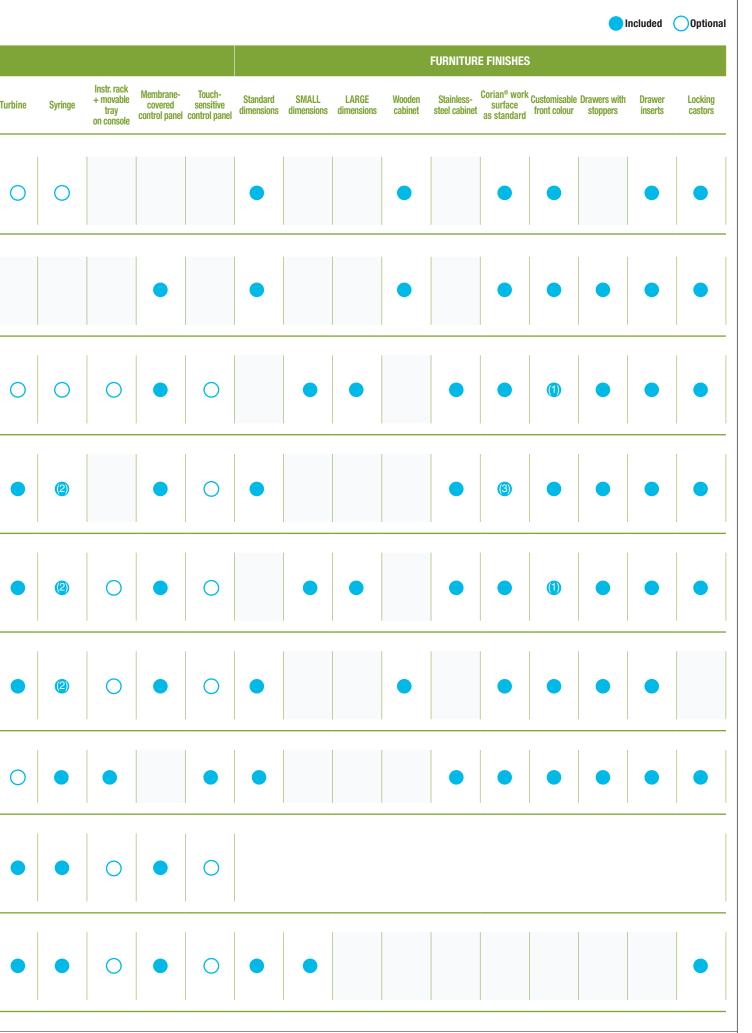







TREATMENT UNITS

U Play

C spinned

SIDAS

E

ő





SIDAS (0)

U Play 1. Cabinet

Cabinet with wood finish.


- 1 storage module
- 1 Corian® work surface.
- 4 removable compartmentalised drawer inserts as standard.
- Stainless-steel flush handles.
- Locking casters.
- 1 technical module.

2. Ergonomics

- Hanging cords.
- ON/OFF control pedal.
- Front-mounted control keypad.
- 1 instrument rack for 2 instruments.
- The motor unit is accessible via an opening at the front.

Left-handed version: contact us.

FOR DIMENSIONS SEE P.146

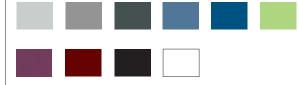

COMPARTMENTALISED DRAWERS

pas @

1 - U Play instrument rack.

- 4 Integrated suction pump.
- 2 Detachable control panel.
- 3 Water reservoir.
- 5 ON/OFF pedal.

Customising the U Play treatment unit (options) LIGHTING (See peripherals on pages 156-157)





Visano lamp PECI010

Opticlux LED magnifying lamp (Tévisio) PECI015

Standard LED magnifying lamp PECI017

3. Working technique

Motor unit with removable internal spray. For treatments within and outside the practice.

- 1 one-piece induction micromotor with internal spray, E-type attachment, speed 1,000 to 40,000 rpm, choice of 3 working techniques.
- 1 integrated Cattani suction pump with motor, collection tray, bacteriological filter and special tube.
- Choice of 3 working techniques.

4. **Options and accessories**

- Adapter pack + tray
- Selection of lighting options p.156-157

U Play - PUSU1PLAY

- Unit supplied without options
- 2-year warranty (not incl. rotary instruments)

COLOURED FRONT OPTION (available in all colours)



CLASSIC 1. Cabinet

Cabinet with stainless-steel finish.

1 storage module

- 1 Corian[®] work surface with integrated front handle, available in black or white. Other colours: contact us.
- 4 drawers + 1 service cupboard.
- Drawers with compartmentalised inserts and push-pull opening system.
- Customisable drawer fronts: classic or premium.
- Locking casters.
- 1 technical module.

2. Ergonomics

- Hanging or retractable cords.
- Swivel console with auto-detect instrument rack for 2 or 4 instruments (depending on customisation option).
- ON/OFF or variable control pedal.
- Hygienic membrane-covered control keypad.

Left-handed version: contact us.

3. Working technique

DRY

(î)

COMPARTMENT

AWFRS

3 - Removable drawer inserts.

4 - Integrated suction pump.5 - Adapter kit + EVO tray.

- 1 carbon brush micromotor without spray.
- 1 handpiece without spray, EX6B NSK.
- 1 integrated Cattani surgical suction pump with adapted tube.
- 1 ON/OFF pedal.

SPRAY

3

- 1 BIEN AIR MCX induction micromotor with internal spray.
- 1 BIEN AIR micro-series handpiece with internal spray.
- 1 BIEN AIR Black Pearl turbine with spray reservoir. *(Supplied without compressor).*
- 1 three-way syringe: air, water, mist. (Supplied without compressor).
- 1 ON/OFF pedal. Can be used with variable control pedal.

This kit requires a compressor.

4. Options and accessories

- Mini-compressor syringe kit.
- Integrated surgical suction pump.
- Touch screen
- Adapter pack + tray
- Selection of lighting options p.156-157

CLASSIC LARGE with hanging cords

- Units supplied without options
- 2-year warranty (not incl. rotary instruments)

RETRACTABLE-CORDS VERSION Contact us

Visano lamp

1 - Instrument rack.

2 - Membrane-covered control panel.

BTC

ENT

FOR DIMENSIONS SEE P.146

Innovative Medical Solutions

0

Halux N30 LED lamp - Flexible arm PECI019

Customising the CLASSIC treatment unit (options)

LIGHTING (See peripherals on pages 156-157)

Column-mounted scialytic lamp PECL004

Bien Air

Dental

COLIBRI

Handleless drawer fronts





TREATMENT L

SIROCCO





SOCCO 1. Cabinet

Cabinet with stainless-steel finish.

1 storage module

- 1 Corian[®] work surface with flush front handle. Other colours: contact us.
- Customisable drawer fronts: classic or premium.
- Locking casters.

1 technical module.

Choice of 'straight' or 'wave' front: contact us.

2 versions: small and large

Ergonomics 2.

- Retractable or hanging cords.
- Swivel console on articulated arm.
- ON/OFF or variable control pedal.
- Membrane-covered control keypad.
- Choice of tray heights on installation.

Left-handed version: contact us.

3. Working technique

DRY

- 1 micromotor without spray.
- 1 handpiece without spray, EX6B NSK.
- 1 integrated surgical Cattani suction pump with a bacteriological filter and adapted tube.
- 1 ON/OFF pedal.



Denta

NSK

- 1 BIEN AIR induction micromotor with internal spray.
- 1 BIEN AIR handpiece with internal spray.
- 1 BIEN AIR Black Pearl turbine with spray reservoir. (Supplied without compressor).
- 1 three-way syringe: air, water, mist. (Supplied without compressor).
- 1 ON/OFF pedal. Can be used with variable control pedal.
- This module requires a compressor.

4. **Options and accessories**

- Mini-compressor syringe kit.
- Integrated surgical suction pump For Cattani suction pump contact us.
- Touch screen.
- Electrically assisted raising/lowering.
- Touch-sensitive keypad: contact us.
- Adapter pack + tray.
- Selection of lighting options p.156-157
- Sirocco motorised arm. For price, contact us.

SIROCCO SMALL with hanging cords

SIROCCO LARGE with hanging cords

 Unit supplied without options • 2-year warranty (not incl. rotary instruments)

RETRACTABLE-CORDS VERSION Contact us



S arm PECL015 scialytic lamp PECL004

scialytic LED lamp PECLOOS

SPRAY

(()

BT(

ENT

FOR DIMENSIONS SEE P.147

Innovative Medical Solutions





INGENii





INGENII 1. Cabinet



FOR DIMENSIONS SEE P.147



Epoxy powder painted steel cabinet. Corian® work surface. Fitted with 8 medical-grade casters with the 4 front casters lockable. Comprises 2 modules that can be connected, rotated or separated.

1 storage module

- 4 drawers.
- 1 Corian® shelf
- Courtesy light.

1 technical module (console and tray / suction pump / compressor / water reservoir).

2. Ergonomics

- Hanging cords.
- Comprising 2 cabinets: one technical compartment and one drawer section, which can be separated and rotated (0 to 30°).
- The 2 cabinets are joined by electromagnets.
- Push-pull opening.
- Rotating and tilting instrument rack.
- UV drawer with LED light.
- Console with 10-inch touch screen control panel on a triple-jointed arm.
- Detachable tilting touch screen.
- All functions of the unit, chair and multimedia support functions can be controlled from the touch screen.
- Management system for up to 5 users.
- Integrated diary function that synchronises with Google Calendar.

Left-handed version: contact us.



INTEGRATED UNIT

Custom treatment-unit arrangement solution tailored to your activity and your practice.



ADVICE AND ASSISTANCE

Contact your technical sales advisor.

S/he will assist you through each step, from choosing to configuring your treatment unit.

VIEW YOUR LAYOUT

With the 3D integration software developed by our Research and Development department, you can see your treatment unit 'installed' in your practice.

TREATMENT UNIT INTEGRATION KIT

Working technique

DRY

- 1 carbon brush micromotor without spray.
- 1 handpiece, without spray, EX6B NSK.
- 1 integrated surgical Cattani suction pump with a bacteriological filter and adapted tube.
- 1 ON/OFF pedal.
- Offset arm.

SPRAY




NSK

Dynamic instruments used for spray-technique treatments, comprising:

- 1 Bien Air MCX brushless micromotor with internal spray, E-type attachment, speed 1,000 to 40,000 rpm.
- 1 Bien Air micro-series handpiece with internal spray, E-type attachment.
- 1 Bien Air Black Pearl turbine with push-button system. Speed 310,000 rpm.
- 1 three-way straight syringe: air, water, mist.
- 1 ON/OFF pedal.
- Offset arm.

This module requires a compressor.

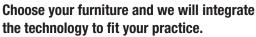

A custom solution for a unique treatment unit tailored to your preferred layout.



TREATMENT UNIT INTEGRATION KIT

Hanging cords with offset arm

• 2-year warranty (not incl. rotary instruments)



CONSULT US FOR AVAILABILITY IN YOUR COUNTRY.


To ensure satisfactory kit integration, furniture must meet SIDAS

142 www.podiatech.com


PLEASE NOTE

specifications.

CONTACT US FOR ANY SPECIAL REQUESTS

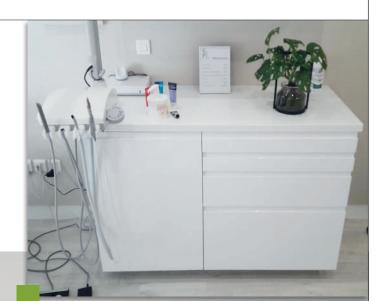

TREATMENT UNITS

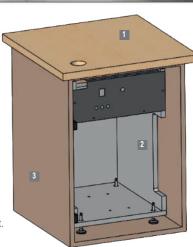

1 - Customer work surface.2 - Podiatech integration kit.

3 - Customer carcass.



KART UK NANO 1. Nano version



Sold without lamp FOR DIMENSIONS SEE P.147



- 1 Membrane-covered control panel.
- 2 Touch-screen control panel.
- 4 Suction pump inside unit. 5 - Adapter + EVO tray kit. 3 - Multifunction variable control pedal





Ergonomics

- Retractable or hanging cords.
- Swivel console with auto-detect instrument rack for 3 or 4 instruments (depending on customisation option).
- ON/OFF or variable control pedal.
- Membrane-covered control keypad.

3. Working technique



- 1 micromotor without spray.
- 1 handpiece, without spray, EX6B NSK.
- 1 Cattani surgical suction pump inside unit.
- 1 ON/OFF pedal.

SPRAY

DRY



Dynamic instruments used for spray-technique treatments, comprising:

- 1 Bien Air MCX brushless micromotor with internal spray, E-type attachment, speed 1,000 to 40,000 rpm.
- 1 Bien Air micro-series handpiece with internal spray, E-type attachment.
- 1 Bien Air Black Pearl turbine with push-button system. Speed 310,000 rpm.
- 1 three-way straight syringe: air, water, mist.
- 1 multifunction variable control pedal: operates the micromotor and turbine and activates the spray and chip-blower.

Options and accessories

- Selection of lighting options p.156-157.

KART UK SUCTION PUMP IN HOUSING - PASP006

TOUCH-SCREEN CONTROL PANEL (UNIT + CHAIR) PUS006

Three-setting rotation speed memory for each instrument plus chair controls when chair pedal is not used (depending on model).

ADAPTER AND EVO TRAY PACK - PUS037


KART UK spray technique

- Kart supplied without options
- 2-year warranty (not incl. rotary instruments)
- · Left-handed version: contact us

DRY TECHNIQUE VERSION WITH SUCTION PUMP - Contact us

RETRACTABLE-CORDS VERSION Contact us



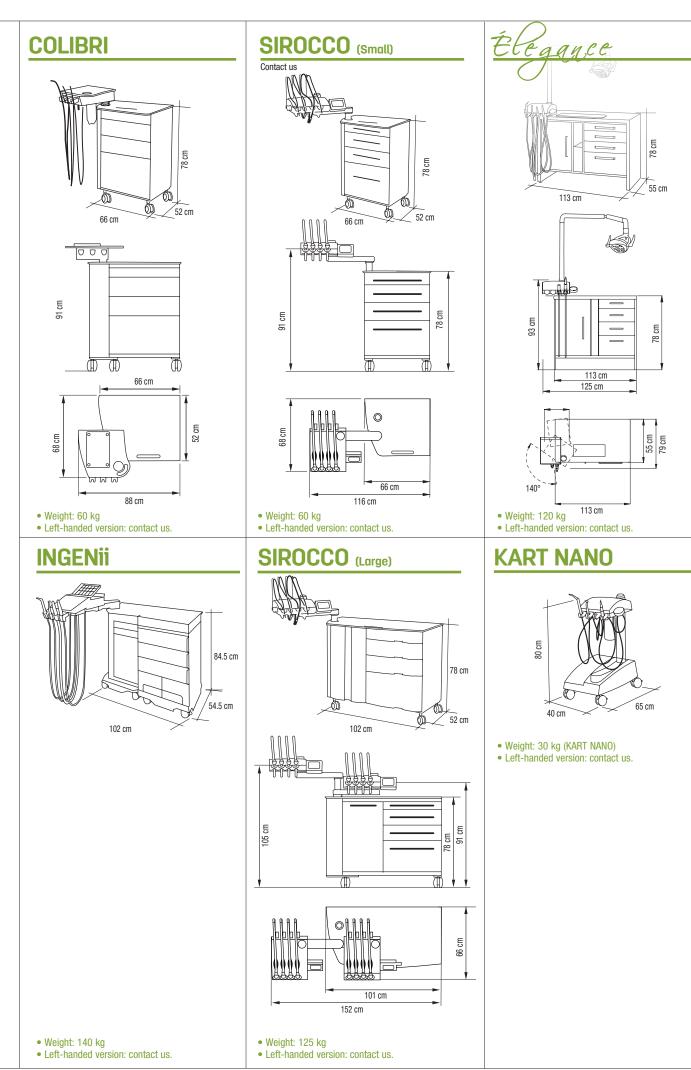

TREATMENT UN







AIR INSTRUMENTS AND ELECTRICAL INSTRUMENTS




This Swiss brand has been the leading make on the podiatry market since 1959.



Hardwearing, superior-quality, high-tech instruments that make the practitioner's work easier.

INSTRUMENT TYPE			MICROMOTORS				
		Post Seller Seller Seller	Dal zman Iteration			Con Arter Bandar	COLUMN TRANSPORT
	Bien Air	NSK			Bien Air		
Model	BLACK Pearl	PANA-MAX	BORALINA	BORA	BORA LED	AQUILON	MC2 IR
PODIATECH ref.	PTUR002	PTUR011	PTUR009	PTUR001	PTUR010	PMIC014	PMIC003
Technology							Carbon brush
UNIFIX rotating coupling		NO	YES (not supplied)	YES (not supplied)	YES (not supplied)		
SPRAY	SINGLE	SINGLE	SINGLE	TRIPLE	TRIPLE	YES INTERNAL	YES INTERNAL
Light					YES		
Bearings	s. steel	s. steel	s. steel	ceramic	ceramic		
Cooling						air	air
Maximum speed (rpm)	310,000	450,000	310,000	320,000	320,000	20,000	40,000
Max. power output						16W	55 W
Weight	55 g	46 g	61 g	75 g	75 g	90 g	74 g
Parts and labour warranty (not incl. wearing parts)	12 months	12 months	12 months	24 months	24 months	24 months	18 months

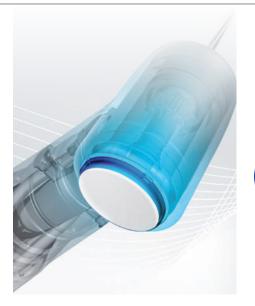



After-sales: Ask us about repairing your BIEN AIR & NSK instruments.

MICROMOTORS					HANDPIECES								
Bost Seller	Bost Soller	(Complete Complete	(Second		R Bernde	E MOX E Ber Ar	MCX Ben Ar restore	Rost Soller	Bost Soller	There are a second and a second		Dende	C C R TLANK X/S
NS	K	Bier	thrologies	NSK	-	Bien Air		NS	SK	Bien Air Medical Technologies	PODIATECH	Bien Air Medical Technologies	NSK
TI-MAX M40 XS N	TI-MAX M40 XS LUM	MC2 GTIR	MC2 GTAV	VOLVER i7	МХІ	MCX LED	MCX BIEN AIR WITH SPRAY	PAM EX6B	PAM M65	PM 1:1 Micro Séries	PM WITHOUT SPRAY	PM 1:1 Spray	TI-MAX X65L NSK
PMIC028	PMIC031	PMIC013	PMIC002	PMIC036	PMIC015	PMIC024	PMIC001	PPMA008	PPMA009	PPMA002	PPMA005	PPMA001	PPMA007
Carbon brush	Carbon brush	Carbon brush	Carbon brush	Carbon brush	Induction	Induction	Induction						
YES INTERNAL	YES	YES	NO	NO	NO	YES INTERNAL	YES INTERNAL	EXTERNAL	YES INTERNAL	YES INTERNAL	NO	INTERNAL	YES INTERNAL
NO	YES					YES							YES
			self-	self-	self-								
air		air	ventilated	ventilated	ventilated	air	air	air	air				
40,000	40,000	40,000	40,000	30,000	40,000	40,000	40,000	40,000	40,000				
		75W	74W		140 W	110W	110W						
97 g	104 g	75 g	110g	192 g	119g	79 g	78 g	54 g	80 g	70 g	87 g	80 g	56 g
12 months	12 months	18 months	18 months	12 months	36 months	36 months	36 months	12 months	12 months	24 months	12 months	24 months	12 months

MICROMOTORS WITHOUT SPRAY



1

4

5

6- LUBRIMED GREASE GUN - PNIR005 - Sold singly.

MICROMOTORS WITH SPRAY



HANDPIECES



154 www.e-podiatech.com



- Triple spray turbine.

- Separate air/water cooling.
- Push-button bur chuck with heat-guard technology.
- Fast connection to UNIFIX rotating coupling (not supplied)
- Ceramic ball bearings for longer service life.
- LED light comprising single-piece optical glass bars.
- Autoclave-compatible at 134°C.
- 24-month parts and labour warranty (does not include wearing parts).






www.e-podiatech.com 155

LIGHTING



LIGHTING

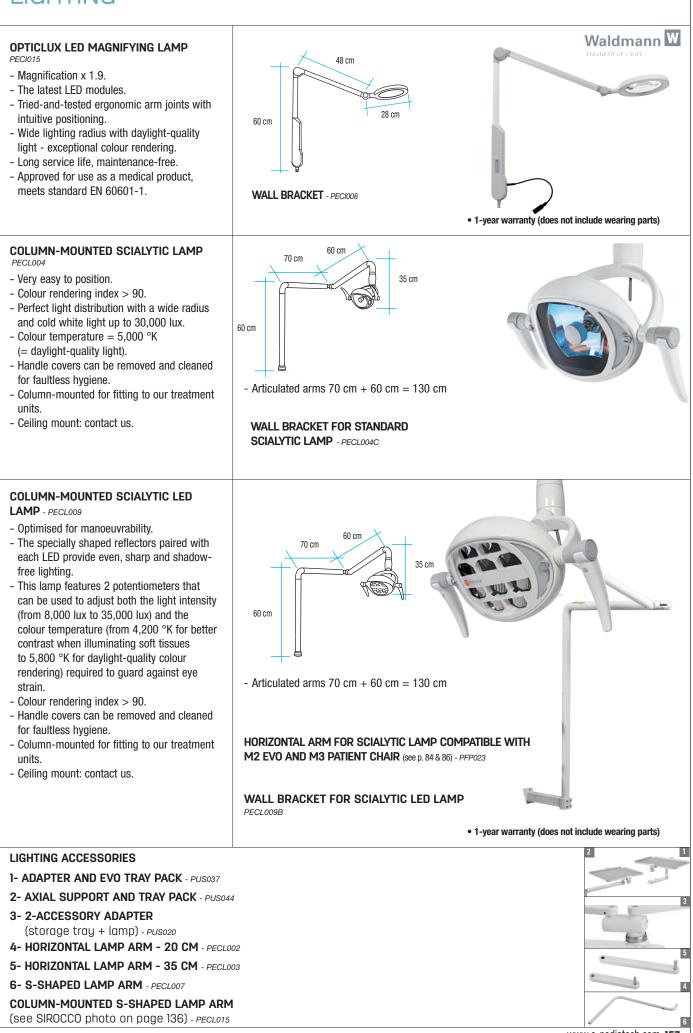

www.e-podiatech.com 157

COMPRESSORS



DRY AIR COMPRESSORS

NON-ENCLOSED SINGLE-CYLINDER DRY-AIR COMPRESSOR - PCOMKIT001

SINGLE-CYLINDER DRY-AIR MOTOR with dryer and automatic purging system with condensate tank.

- Hour meter.
- No purging required.
- 30-litre tank with internal food-grade resin lining.
- Maximum pressure: 8 bar.
- Effective air yield at 5 bar: 67.5 l/min.
- Operating pressure: 5.5 to 7.5 bar.
- Noise level 70 dB.
- Soundproofing cabinet: reduces compressor noise by approximately 20 dB.
- Weight: 36 kg
- 230 V 550 W 3.8 A




• 3-year warranty (does not include wearing parts)

SURGICAL ASPIRATORS

UNI-JET 75 MOTOR KIT- PASP001

ASPIRATION KIT: UNI-JET 75 motor - 75 m3/h + DUST COLLECTOR + BACTERIOLOGICAL FILTER Operating head: 1500 mm water. Comprises: - 1 Uni-Jet 75 motor: length 232 mm, depth 246 mm, ht. 251 mm. 62 dB. - 230 V - 400 W - 3.1 A - 1 nail dust collector: length 185 mm, depth 125 mm, height 371 mm. - 1 bag for the nail dust collector. - 1 active-charcoal bacteriological filter: diam. 137 mm, height 225 mm. - 1 set of tubes to connect the motor to the collector and filter.	CATTANI Dir technology
Compatible with our treatment units.	- Communition
Can also be installed in a separate cabinet (maximum distance 10 m).	
UNI-JET 75 MOTOR IN MONO-LABOR HOUSING - PASP006	
UNI-JET 75 MOTOR in metal housing, with switch and outlet for equipment connection tube.	
Comprises:	CATTANI
- 1 Uni-Jet 75 motor. - 1 nail dust collector.	air technology
 - 1 bag for the nail dust collector. - 1 set of tubes to connect the motor to the collector and filter. 	
- 58 dB	•
- 230 V - 400 W - 3.1 A	• 1-year warranty (does not include wearing parts)
MAINTENANCE AND ACCESSORIES FOR SURGICAL ASPIRATORS	
BACTERIOLOGICAL FILTER - PASP003	
PACK OF 5 ASPIRATION BAGS - PASP004	
MAINTENANCE ADVICE FOR COMPRESSORS AND ASPIRATORS	
Our experienced professionals are on hand to provide clearly explained answers to your technical requi	ests and questions:
diagnoses over the phone, solutions. Contact the PODIATECH technical department from Monday to Friday, 8.30 to 12.30 and 13.30 to 17.00), on + 33 (0)4 76 93 46 06 .
	www.e-podiatech.com 159

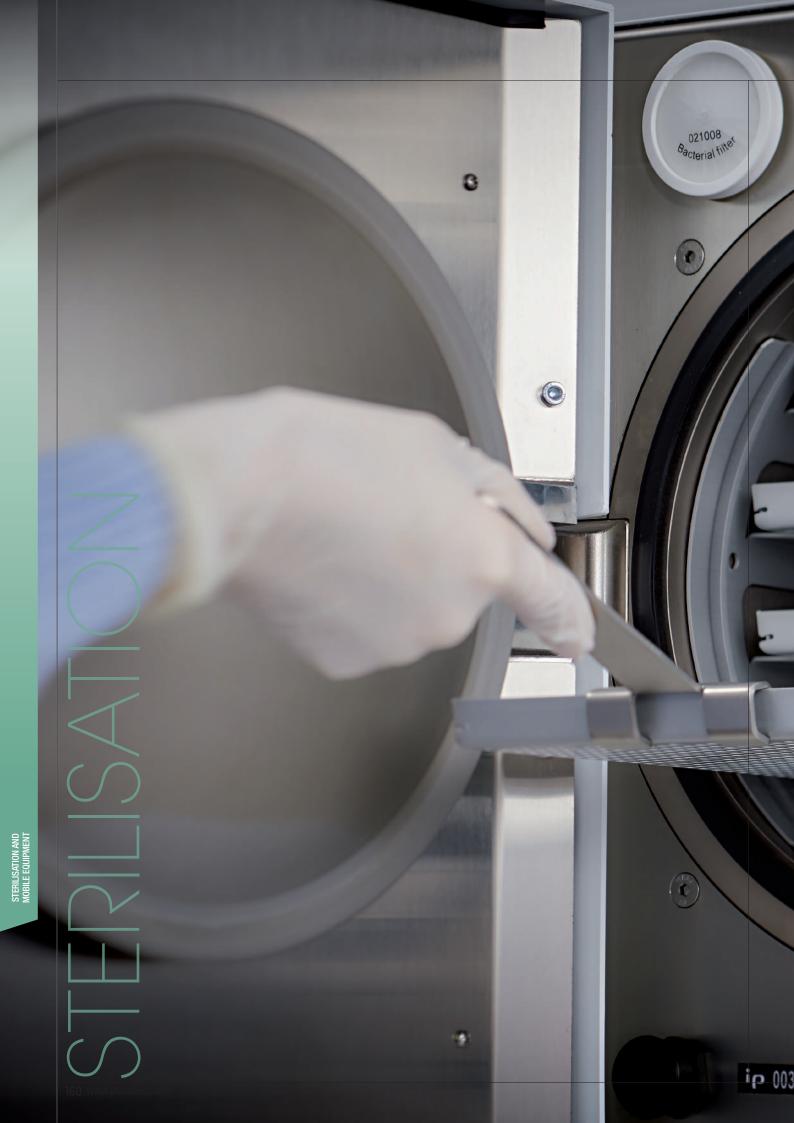







STERILISATION - 3 ESSENTIAL STEPS

TRACEABILITY

DID YOU KNOW? French ministerial circular No. 138 of 14 March 2001 states that "the traceability of actions, equipment and procedures for handling medical devices is compulsory [...]".

PACKAGING

DID YOU KNOW? The packaging stage prior to instrument sterilisation is essential for ensuring that the sterile state is maintained until next use with a patient.

MAINTENANCE

DID YOU KNOW? In article D665-5-1 of the French Public Health Code issued by decree No. 2001-1154 of 5 December 2001 - art.1 JORF of 7 December 2001, repealed by decree 2004-802 2004-07-29 art.4 JORF of 8 August 2004: "The operator is required to ensure the continuation of this performance and the maintenance of all medical devices. This obligation shall be verified by a quality control process defined by decree, the cost of which must be borne by the device operator."

SOAKING TANKS

SOAKING TANK, 3 L-*P71D009* Large tank with perforated basket and integrated lid.

- Dimensions: length 35 cm, width 23 cm, height 12 cm.
- Basket dimensions: length 34 cm, width 22 cm, height 3 cm.

SOAKING TANK, 2 L - P71D026

SOAKING TANK, 1,25 L - P71D010

- Practical to use and compact.
- Dimensions: length 26 cm, width 11 cm, height 8 cm.
- Basket dimensions: length 24 cm, width 9,5 cm, height 3 cm.

AUTOCLAVE-COMPATIBLE BRUSH - *P65BR001* For cleaning instruments.

SOAKING POT FOR BURS - *P71D011* For soaking and draining burs.

ULTRASONIC CLEANER, - PBAC006

- Frequency of transducers: 42 kHz.
- Power output: 170 W.
- External dimensions: 340 x 255 x 225 mm.
- Supplied with basket and lid.

NEWSEAL PLUS HEAT SEALER

The Newseal Plus heat sealer is the easy way to seal bags before sterilisation.

- Automatically adjusted sealing from 12 mm.
- Retractable-blade cutter.
- Visual and audible sealing indicator.
- Roll holder.
- Pre-adjusted for wall mounting.

FLAT STERILISATION SLEEVES (rolls sold singly)

200 m, width 50 mm - P72S004
200 m, width 75 mm - P72S001
200 m, width 100 mm - P72S002
200 m, width 150 mm - P72S003

SELF-SEAL BAGS (Box of 200 bags)

60 mm x 100 mm - P72S005 90 mm x 165 mm - P72S017 90 mm x 250 mm - P72S006 140 mm x 260 mm - P72S007









RECTANGULAR ALUMINIUM BOXES

- For storing instruments for home visits.
- Sold singly.

Rectangular aluminium box 18 x 9 x 3 cm - P40B001 Perforated aluminium box, grey, 17 x 7 x 3 cm - P40B011 Aluminium blue box 18 x 9 x 3 cm - P40B012

ALUMINIUM TRAY

- 28.4 cm x 18.3 m. Height 1.7 cm.
- Designed to hold an instrument rack or podiatry tray, tray base of the same dimensions.
- Sold singly.

Aluminium tray - P40P001

ALUMINIUM INSTRUMENT RACK

Designed to be placed inside the tray. Holds 17 instruments.Sold singly.

Aluminium instrument rack - P40P002

PODIATRY TRAY

- Single-use cover for aluminium tray, made of disposable plastic for faultless hygiene. Use a new cover for each patient.
- Sold in boxes of 400.

Podiatry tray - P40P003

BOTTLES

- 3 sizes of bottle: 60, 100, 250 cc, with red perforated cap.

- Sold singly.

60cc - 100cc - 250cc bottle - P40F003-04-05

WASH BOTTLES

- Made of flexible plastic for medical use, with curved nozzle.
- 2 sizes of wash bottle: 100 cc, 250 cc.
- Sold singly.

100cc wash bottle - *P*40F001 **250cc wash bottle** - *P*40F002





MICROMOTORS



MICROMOTORS - CASES - HYGIENE - BOXES



 \mathcal{O} _

168 www.podiatech.com



CONSUMABLES

BLADES & HANDLES

CLIPPERS SCISSORS & BURS p.172







TREATMENTS & CREAMS p.188

p.196


BLADES AND HANDLES

STERILE GOUGE BLADES

Box of 50 sterile stainless-steel gouge blades, packaged individually in foil.

No.1 - No.2 - No.3 - No.5 - No.8 - P41G_FIRST_1 to 8 No.10 - P41G_FIRST_10



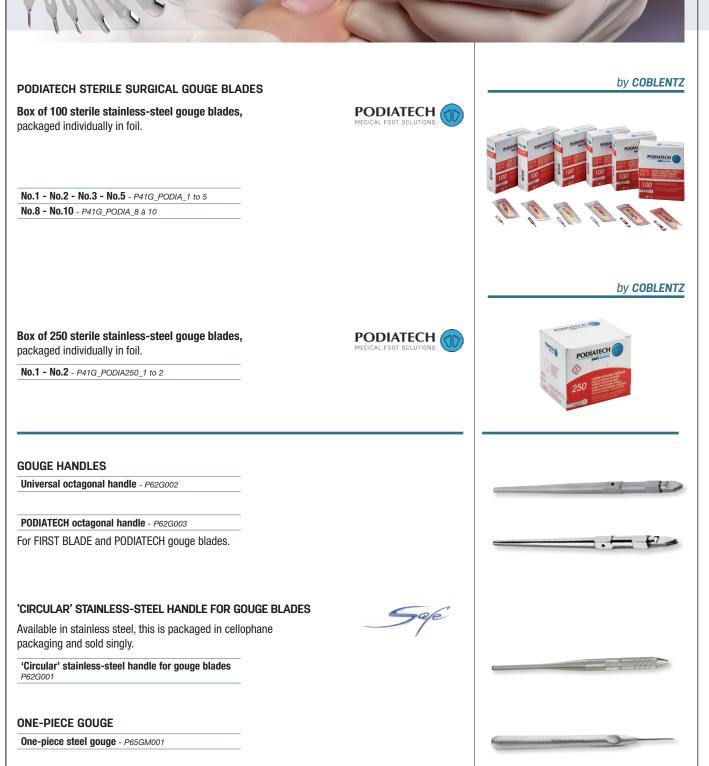

BLADES AND HANDLES

ESSENTIALS

BLADES AND HANDLES		
SWANN-MORTON STERILE SCALPEL BLADES Box of 100 sterile SF 100 stainless-steel blades,	Swann-Morton	
packaged individually in foil.		CARDON STEEL SANDAR STEEL SANDA
No.10 - No.11 - No.15 - No.20** - P41B_10 to 23		
SWANN-MORTON SCALPEL HANDLES Sold singly.	Swann-Morton	
Handle No. 3 for blades No. 10, 11, 15 - $P62B_3$ Handle No. 4 for blades No. 20 - $P62B_4$ Handle No. 9 for blades No. 23 - $P62B_7$		
B BRAUN STERILE SCALPEL BLADES Box of 100 sterile carbon-steel blades, packaged individually in foil.	BBRAUN	
No.10 - No.11 - No.15 - No.20 - P41B_BB510 to 20		
PODIATECH SCALPEL HANDLES	PODIATECH	
Handle No. 3 for blades No. 10, 11, 15 - <i>P62B_5</i> Handle No. 4 for blades No. 20 - <i>P62B_6</i>		
SWANN-MORTON STERILE CHISEL BLADES Box of 25 sterile SF 100 stainless-steel blades, packaged individually	Swann-Moston	
1 narrow tip - SM61S - P41F_SM61S 2 narrow tips - SM61 - P41F_SM61 1 wide tip - SM62S - P41F_SM62S		
2 wide tips - SM62 - P41F_SM62		And the second s
SWANN-MORTON CHISEL HANDLE	Swann-Morton	
SHARPS BOXES		6
An essential safety product, these boxes are for the safe collection of used and soiled blades, preventing risk of injury or secondary contamination.		
250 ml box - <i>P40B003-250</i> 400 ml box - <i>P40B003</i>		
CORN-CUTTER BLADES AND CREDO STERILISABLE STAINLE Box of 100 stainless-steel blades, packaged individually in packs of 10.	SS-STEEL HANDLE	
Corn-cutter blades - P41C001		110
Available in stainless steel, the handle is individually packaged. Sold singly.		Cose C
CREDO sterilisable stainless-steel handle - P62C001		

NAIL CLIPPERS

AESCULAP[®]

- The world's leading manufacturer of patient-care instruments. - The go-to manufacturer of medical and surgical instruments. - Packaged in a box with a spare spring and sold singly.



- Long-lasting. precision-cutting instruments.
- Ergonomically designed for user comfort.
- Packaged in a box with transparent lid and sold singly.

Ingrown nail clippers

INGROWN NAIL CLIPPERS. 11.5 CM AESCULAP

Without clasp, straight flat and slim jaws curved at the tips. Cutting length 12 mm, straight cut. Supplier code: HF475R.

Ingrown nail clippers. 11.5 cm - P61CD007

INGROWN NAIL CLIPPERS. 11.5 CM SAM

Without clasp, flat and slim jaws. Cutting length 14 mm, straight cut. Supplier code: 613Al115SA.

Ingrown nail clippers. 11.5 cm - P61CD002

INGROWN NAIL CLIPPERS. 12 CM AESCULAP

Double spring without clasp, flat narrow jaws, cutting length 16 mm, straight cut. Supplier code: HF482R.

Ingrown nail clippers. 12 cm - P61CD012

INGROWN NAIL CLIPPERS. 13 CM AESCULAP

Without clasp, straight flat and slim jaws curved at the tips. Cutting length 15 mm, straight cut. Supplier code: HF483R.

Ingrown nail clippers. 13 cm - P61CD006

INGROWN NAIL CLIPPERS. 13 CM AESCULAP

Without clasp, straight flat and very slim jaws curved at the tips. Cutting length 16 mm, straight cut. Supplier code: HF478R.

Ingrown nail clippers. 13 cm - P61CD008

INGROWN NAIL CLIPPERS. 13 CM SAM

Without clasp, ovoid flat and slim jaws. Cutting length 15 mm, straight cut. Supplier code: 613Al130SA

Ingrown nail clippers. 13 cm - P61CD003

INGROWN NAIL CLIPPERS. 13 CM SAM

Double spring, without clasp, ovoid flat and slim jaws. Cutting length 17 mm, straight cut. Supplier code: 612AI130SARD

Ingrown nail clippers. 13 cm - P61CD010

INGROWN NAIL CLIPPERS. 13.5 CM SAM

With clasp, straight flat and slim jaws. Cutting length 18 mm, straight cut. Supplier code: 614Al130SA


Ingrown nail clippers. 13.5 cm - P61CD004

INGROWN NAIL CLIPPERS. 13.5 CM SAM

With clasp, straight flat and very slim jaws. Cutting length 18 mm, straight cut. Supplier code: 615Al135SA

Ingrown nail clippers. 13.5 cm - P61CD001

172 www.podiatech.com

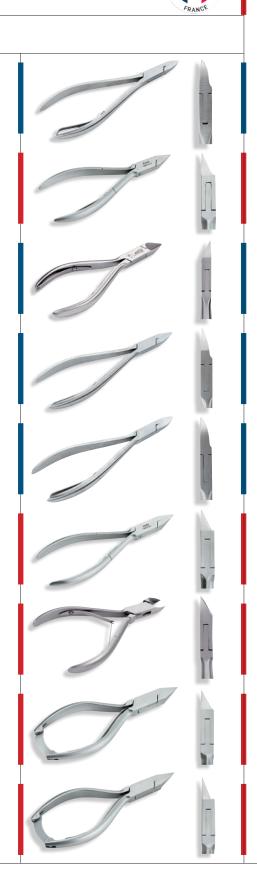



NAIL CLIPPERS

SHARPENING SERVICE, ALL BRANDS

The conditions of circulation of medical items are subject to certain safety rules. We accept clippers and gouges for sharpening only when sent in sterile packs.

P66A001

P66A002

CLIPPERS SHARPENING, ALL BRANDS

ONE-PIECE GOUGE SHARPENING

CLEANING, STERILISATION AND MAINTENANCE - PAGE 174

Standard nail clippers

NAIL CLIPPERS, 12.5 CM AESCULAP

Without clasp, concave curved jaws, cutting length 18 mm, curved cut. *Supplier code: HF201R*

Nail clippers, 12.5 cm - P61CA005

NAIL CLIPPERS, 13 CM AESCULAP

Without clasp, concave jaws curved at the tips. Cutting length 19 mm, straight cut. *Supplier code: HF479R*

Nail clippers, 13 cm - P61CD015

NAIL CLIPPERS, 13 CM AESCULAP

Without clasp, concave jaws curved at the tips. Cutting length 20 mm, curved cut. *Supplier code: HF203R*

Nail clippers, 13 cm - P61CA003

NAIL CLIPPERS, STANDARD 13.5 CM SAM

Double spring, with clasp, standard straight curved flat jaws. Cutting length 15 mm. *Supplier code: 609Al135SARD*

Standard nail clippers, 13.5 cm - P61CD011

NAIL CLIPPERS, STANDARD 13.5 CM SAM

With clasp, standard straight curved flat jaws. Cutting length 15 mm. *Supplier code: 609Al135SA* Standard nail clippers, 13.5 cm - *P61CD005*

NAIL CLIPPERS, 13.5 CM SAM

With clasp, concave curved jaws, curved cut. Cutting length 16 mm. *Supplier code: 610AI135SA*

Standard nail clippers, 13.5 cm - P61CA001

NAIL CLIPPERS, 14 CM SAM

With clasp, oblique jaws, curved cut. Cutting length 18 mm. *Supplier code: 618Al140SA*

Nail clippers, 14 cm - P61CA002



We offer

MORE

SHARPENING

OTHER NAIL CLIPPERS

CLEANING, STERILISATION AND MAINTENANCE OF TWEEZERS AND SCISSORS

1. Pre-disinfection: decontamination of tweezers by immersion in a bath.

2. Manual cleaning with brushes or in an ultrasonic cleaner.

3. Rinsing and drying, manually or in a drier.

4. Lubrication and bagging.

IMPORTANT: Stainless-steel and nickel-plated instruments must not be mixed because they may become contaminated by brown marks.

5. Sterilisation and storage. Prion Cycle at 134 °C.

Special nail clippers

NAIL NIPPER, 11 CM OR 14 CM SAM

With clasp, standard straight jaws. 11 cm nipper: cutting length 15 mm; 14 cm nipper: cutting length 17 mm.

Nail nipper, 11 cm - P61SE001

Nail nipper, 14 cm - P61SE002

HANGNAIL CLIPPERS, 12 CM SAM

With clasp, 3mm oblique straight jaws.

Hangnail clippers, 12 cm - P61P001

PÉAN TWEEZERS

Available in stainless steel, 14 cm, can be used to grasp cotton wool. **Sold singly.**

Péan tweezers - P61PR002

KOCHER TWEEZERS

Available in stainless steel, 14 cm, can be used to grasp gauze. **Sold singly.**

Kocher tweezers - P61PR003

SPLINTER TWEEZERS

Available in stainless steel, 10 cm, can be used to grasp a cut corner of nail or a splinter. **Sold singly.**

Splinter tweezers - P61PR001

ADSON TWEEZERS

Available in stainless steel, 12 cm, grooved jaws enabling it t o be used to grasp a cut corner of nail. **Sold singly.**

Adson tweezers - P61PR004



SCISSORS & CURETTES

Τ

AESCULAP LUBRICATING OIL	
Recommended for maintaining the proper functioning of joints, spring tab operation on the handle and clasp lock, before sterilisation. 50ml bottle, sold singly.	
AESCULAP lubricating oil - P61N001	
	C
Scissors & Curettes	
BLUNT-TIP SCISSORS	0
Available in stainless steel, 14 cm, for cutting bandages, plasters and gauze. Sold singly.	0
Blunt-tip scissors - P65C001	0
POINTED-TIP SCISSORS	\bigcirc
Available in stainless steel, 14 cm, for cutting bandages with a central depression, plasters and gauze. Sold singly.	
Pointed-tip scissors - P65C002	
SLIM POINTED-TIP SCISSORS, STRAIGHT	\bigcirc
Available in stainless steel, 10 cm, for cutting orthoplastic devices. Sold singly.	
Slim pointed-tip scissors, straight - P65C003	\bigcirc
SLIM POINTED-TIP SCISSORS, CURVED	
Available in stainless steel, 10 cm, for cutting and sculpting orthoplastic devices. Sold singly.	0
Slim pointed-tip scissors, curved - P65C004	
EXCAVATOR CURETTE, 15 CM	
Excavator curette, 15 cm - P65EX002	
STAINLESS-STEEL RIFFLER FILE, 20 CM	
Stainless-steel riffler file, 20 cm - P42R010	

HANDPIECE BURS

edenta 🖯

ZIRCONIUM BURS (ZrN – zirconium nitride)

Hard and long-lasting - quiet grinding - reduced heat generation

	g quiet gimanig rouad	ou nout gonorauon			
		AC flutes AC	46AC Coarse flutes AC	75AC Plain cross-cutting flutes AC	53AC Fine and coarse flutes AC
Reference	45AC	45AC	46AC	75AC	53AC
ø in mm	031	040	040	040	060
Reference	P63CZ001	P63CZ002	P63CZ003	P63CZ005	P63CZ004
Maximum speed	30,000 rpm	30,000 rpm	30,000 rpm	30,000 rpm	20,000 rpm
Optimum speed	15,000 rpm	15,000 rpm	15,000 rpm	15,000 rpm	15,000 rpm

'STANDARD' TUNGSTEN CARBIDE BURS

Optimum comfort and effectiveness for all types of grinding. Stainless-steel stem. Vibration-free rotation.

			ŧ				
Reference	0720	6220	5620		5220	5420	5520
ø in mm	023	040	045		060	060	060
Reference	P63CT001	P63CT002	P61CT0		SCT003	P61CT016	P61CT017
Maximum speed	40,000 rpm	30,000 rpm	30,000 r	-	00 rpm	20,000 rpm	20,000 rpm
Optimum speed	15,000 rpm	15,000 rpm	15,000 r	rpm 12,0	00 rpm	12,000 rpm	12,000 rpm
Reference	0110	6410 001	5710	5610 045	5210 060	5410 060	5510 060
ø in mm	023	031	040	045	060	060	060


			Â				4
Reference	5650	5350	5550	Reference	5455	5855	6055
ø in mm	045	060	060	ø in mm	060	060	070
Reference	P63CT008	P61CT022	P61CT023	Reference	P63CT009	P63CT010	P63CT011
Maximum speed	30,000 rpm	20,000 rpm	20,000 rpm	Maximum speed	20,000 rpm	20,000 rpm	20,000 rpm
Optimum speed	15,000 rpm	12,000 rpm	12,000 rpm	Optimum speed	12,000 rpm	12,000 rpm	10,000 rpm

BLUE RING, standard crisscross flutes for filing normal nails. GREEN RING, coarse crisscross flutes for filing thick nails and calluses. RED RING, fine crisscross flutes for filing and finishing thin nails. BLACK RING, very coarse crisscross flutes for filing very thick nails and thick calluses.

176 www.podiatech.com


HANDPIECE BURS



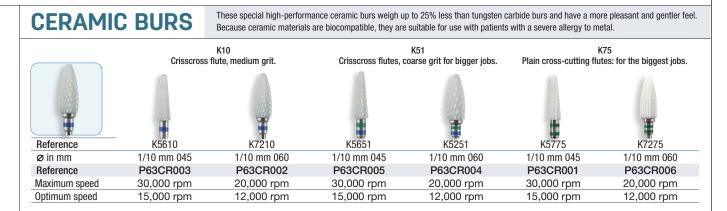

SPECIAL' TUNGSTEN CARBIDE BURS

SPIRAL FLUTES: the geometry of the longitudinal and cross-cutting blades optimises the abrasion action in line with the force and pressure exerted. Ergonomically shaped. Leave a smooth, unmarked surface. This special design also prevents the flutes from clogging.

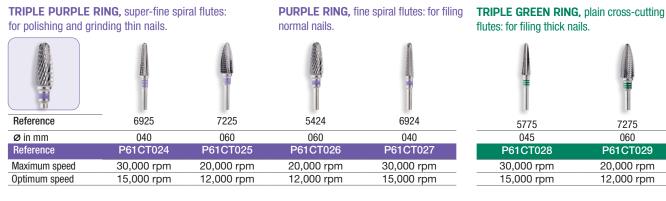

VOLCANO : carbide burs with a special DLC coating (twice as hard as a standard bur). Combines greater grinding power with a greater cutting ability. More stable when grinding and therefore a total absence of heating. Very long-lasting. Improved resistance to sterilisation.

BLUE RING DLC, standard crisscross flu	tes for filing normal nails.	GREEN RING DLC , coar for filing very thick nails a			
				Ø.	A
Reference	DLC 0110	5610	DLC 5510	5650	5450
ø in mm	023	045	060 HP	045	060
Reference	P61CT030	P61CT031	P61CT032	P61CT034	P61CT033
Maximum speed	40,000 rpm	30,000 rpm	20,000 rpm	30,000 rpm	20,000 rpm
Optimum speed	15,000 rpm	15,000 rpm	12,000 rpm	15,000 rpm	12,000 rpm

LEFT-HANDED BURS: enables a left-handed operator to work safely and accurately and provides a clear view of the work surface. **DOUBLE RED RING:** fluted for left-handed operator.

		ŧ	-	
Reference	0110	5620	5610	5410
ø in mm	023	045	045	060
Reference	P63CT012	P63CT013	P63CT014	P63CT015
Maximum speed	40,000 rpm	30,000 rpm	30,000 rpm	20,000 rpm
Optimum speed	15,000 rpm	15,000 rpm	15,000 rpm	12,000 rpm

HANDPIECE BURS



'STANDARD' DIAMOND BURS Only evenly concentrated high-quality diamond pieces ensure top-quality, long-lasting abrasion. There are several shapes available to suit the treatments you offer. ROUND CONE, STANDARD GRIT, WIDE FLAME, COARSE for polishing and working on the side edges of thick nails **GRIT,** for working on heels 850.104 856.104 856.104 854R.104 5860.104 Reference ø in mm 014 016 025 040 060 P63D012 P63D001 P63D002 P63D003 P63D004 Reference 30,000 rpm 30,000 rpm Maximum speed 40,000 rpm 30,000 rpm 20,000 rpm Optimum speed 25,000 rpm 20,000 rpm 20,000 rpm 20,000 rpm 15,000 rpm **BALL, STANDARD GRIT,** THIN CONE, STANDARD GRIT, for working on periungual grooves and side edges of thin nails for enucleation and nail piercing 858.104 801.104 801.104 Reference 858.104 859.104 801.104 ø in mm 012 016 018 010 014 021 P63D013 P63D014 P63D015 P63D009 P63D010 P63D01 Reference Maximum speed 40,000 rpm 40,000 rpm 40,000 rpm 40,000 rpm 40,000 rpm 40,000 rpm Optimum speed 20,000 rpm 20,000 rpm 20,000 rpm 25,000 rpm 25,000 rpm 20,000 rpm **'SPECIAL' DIAMOND BURS** Rubynit trimmer: special bur featuring a metallic binder and grains of natural ruby. Recognised for its abrasive action ideal for treating thin nails with fungal infections and for calluses.

	T	T		T	
Reference	3108.104	3110.104	3106.104	3101.104	3112.104
ø in mm	055	055	065	085	050
Reference	P63D005	P63D006	P63D007	P63D008	P63D021
Maximum speed	20,000 rpm	20,000 rpm	15,000 rpm	15,000 rpm	20,000 rpm
Optimum speed	15,000 rpm	15,000 rpm	10,000 rpm	10,000 rpm	15,000 rpm


Hollow burs: designed for effective abrasion of large areas while creating only a small amount of heat. Maximum rotation speed: 10,000 rpm. Apply a light pressure only.

Reference	490.104	490.104
ø in mm	090	110
Reference	P63D018	P63D019
Maximum speed	15,000 rpm	15,000 rpm
Optimum speed	10,000 rpm	10,000 rpm

Dia Duo: bur with a special diamond grain with grooves that creates less heat, for optimum abrasion of heel calluses. Its crisscross ridges optimise diffusion of the spray over the work area.

Reference	894.104	860.104
ø in mm	065	085
Reference	P63D017	P63D016
Maximum speed	15,000 rpm	15,000 rpm
Optimum speed	10,000 rpm	10,000 rpm

TURBINE BURS

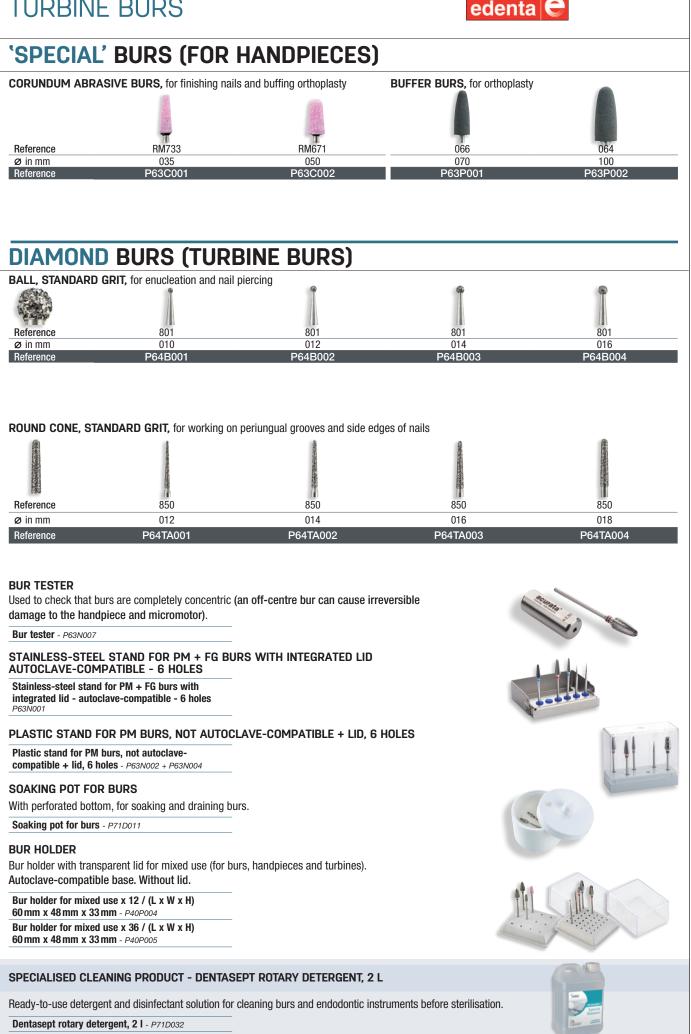

BANDAGES/BAND-AIDS AND PROTECTION

DIGITAL PAD

- Sold in boxes of 6.
- Elastic tube with circular cushion of polymer gel.
- Protects corns from pressure and friction.

Small (S) for the little toe Medium (M) for the middle toes - P45P004-05-06 Large (L) for the big toe

DIGITAL TUBE

- Length: 15 cm.
- Sold in boxes of 2 tubes.
- Elastic tube fully lined with polymer gel.
- Protects the whole toe. Cut to the required length.

Small (S) for the little toe Medium (M) for the middle toes - P45P007-08-09 Large (L) for the big toe

DIGITAL STRIP

- Available in 3 sizes: Small, Medium and Large. Length: 71 cm.
- Sold singly.
- Thick gel for additional protection. Relieves pressure on the top, underside or side of the toe. Cut to the required length.

Small (S) for the little toe Medium (M) for the middle toes - P45P010-11-12 Large (L) for the big toe

TUBOPROTE

- Sold in bags of 12 tubes.

Tuboprote No. 1: 15 mm diameter - P45P013
Tuboprote No. 2: 18 mm diameter - P45P014
Tuboprote No. 3: 21 mm diameter - P45P015
Tuboprote No. 4: 25 mm diameter - P45P016

FLEECY WEB

- Adhesive fabric with 2mm fleece lining.
- 100% cotton, hypoallergenic.
- Permeable.
- Pack of 4 sheets, 22.5 x 40 cm.
- Protection and cushioning. Cut to the required shape and size.

Fleecy web - P45P001

Fleecy web extra - P45P029

MOLESKIN

- Extra-thin 1mm adhesive fabric.
- 100% cotton.
- Permeable, prevents soaking.
- Roll, 23 cm x 0.70 m.
- Protects inflamed areas.

Moleskin - P45P002

FLEECY FOAM

- 5mm laminated polyester felt with polyurethane cellular foam.
- Resumes initial thickness immediately after compression.
- Permeable.
- Pack of 4 sheets, 22.5 x 40 cm.

Fleecy Foam - P45P003



BANDAGES/BAND-AIDS AND PROTECTION

BANDAGES/BAND-AIDS

COVERMED

Boxes sold singly. Roll, 5 m x 6 cm, wide - P44P002 Roll, 5 m x 8 cm, wide - P44P011

HYPAFIX

Boxes sold singly.

Roll, 5 m x 2.5 cm wide - P44P008	
Roll, 5 m x 5 cm wide - P44P003	
Roll, 2 m x 10 cm wide - P44P004	
Roll, 10 m x 5 cm wide - Transparent - P44P010	

TENSOPLAST

Boxes sold singly.

Roll, 2.5 m x 3 cm wide - P44P005 Roll, 2.5 m x 6 cm wide - P44P009

COTTON WOOL

Cotton wool balls - Sold in packs of 700 balls P44C001 Cotton wool squares - Sold in packs of 500

squares - P44C002

COTTON WOOL DISPENSER

- Diameter 10 cm.
- Height 10 cm
- Sold singly.

Cotton wool dispenser without spring - P40B002

COMPRESSES

- Sterile nonwoven compresses, 5 x 5 cm.

Sold in boxes of 50 packs of 5 compresses - P44P006

- Non-sterile nonwoven compresses.

Sold in packs of 100 compresses, 5 x 5 cm P44P012

Sold in packs of 100 compresses, 7.5 x 7.5 cm P44P013

GAUZE

- Tubular gauze No.1.5, particularly recommended for toe bandages.
- Diameter 19 mm.

- Boxes sold singly.

20m roll - P44P007













ACCESSORIES

RASPS

FOOT RASP

- 2 abrasive sides, coarse grit and fine grit.
- Sold singly.
- Colours: green white blue purple P42R005-4-6-7

AUTOCLAVE-COMPATIBLE RASP

- Autoclave-compatible and sterilisable stainless-steel rasp.
- Sold singly.
- Autoclave-compatible rasp P42R011

ABRASIVE STRIPS FOR AUTOCLAVE-COMPATIBLE RASP

- Sold in packs of 10.
- Coarse-grit P42R012

Fine-grit - P42R014

HANDLE FOR ABRASIVE STRIP

- Sold singly.

Handle for abrasive strip - P42R003

ABRASIVE STRIPS

- Sold in packs of 25 strips.

Coarse-grit / fine-grit abrasive strips - P42R001-02

RECTANGULAR ALUMINIUM BOXES

- For storing instruments for home visits.
- Sold singly.

Rectangular aluminium box 18 x 9 x 3 cm - P40B001 Perforated aluminium box, grey, 17 x 7 x 3 cm - P40B011 Aluminium blue box 18 x 9 x 3 cm - P40B012

ALUMINIUM TRAY

- 28.4 cm x 18.3 m. Height 1.7 cm.
- Designed to hold an instrument rack or podiatry tray, tray base of the same dimensions.
- Sold singly.

Aluminium tray - P40P001

ALUMINIUM INSTRUMENT RACK

- Designed to be placed inside the tray. Holds 17 instruments.
- Sold singly.

Aluminium instrument rack - P40P002

PODIATRY TRAY

- Single-use cover for aluminium tray, made of disposable plastic for faultless hygiene. Use a new cover for each patient.
- Sold in boxes of 400.

Podiatry tray - P40P003

BOTTLES

- 3 sizes of bottle: 60, 100, 250 cc, with red perforated cap.
- Sold singly.

60cc - 100cc - 250cc bottle - P40F003-04-05

WASH BOTTLES

- Made of flexible plastic for medical use, with curved nozzle.
- 2 sizes of wash bottle: 100 cc, 250 cc.
- Sold singly.

100cc wash bottle - P40F001

250cc wash bottle - P40F002



PROTECTION & HYGIENE

ALL TOGETHER, TOUS ENSEMBLE!

Because you are a healthcare professional and your patients are your main concern, our Protection & Hygiene range meets all your different needs.

MASKS

SURGICAL MASKS

- 3-ply healthcare mask, type II.
- Sold in boxes of 50 masks.

3-ply healthcare mask - SIMKALTOG_001

FFP2 MASK

- Sold in packs of 10 masks.

FFP2 mask - SIMKALTOG_024

AMBIDEXTROUS EXAMINATION GLOVES

- Meet European standards: CE.
- Sold in boxes of 100.

POWDERED LATEX EXAMINATION GLOVES

Examination gloves, sizes 5/6 - 6/7 - 7/8 - 8/9 - 9/10 - P43G000-01 to -04

NON-POWDERED LATEX EXAMINATION GLOVES, BLUE

Examination gloves, sizes XS = 5 / S = 6 / M = 7 / L = 8 - P43G009-10 to 12

NON-POWDERED NITRILE EXAMINATION GLOVES, BLUE

Examination gloves, sizes 5/6 - 6/7 - 7/8 - 8/9 - 9/10 - P43G004A-05 to 08

PAPER LINER

- PAPER LINER
- 2-ply paper liner for medical examinations, white tissue web, protects surfaces during clinical
- examination or to cover patient chair leg rests.
- Roll of 150 sheets.
- Rolls sold singly.

Paper liner / Width 50 cm, pre-cut every 35 cm - P43D001

LAMINATED PAPER LINER.

- Packaging: 6 rolls. 1 roll = 132 standard-size sheets, blue, laminated, embossed, glued.

Laminated paper liner / Width 50 cm, pre-cut every 38 cm \pm 5% - <code>P43D003</code>

SOCIAL DISTANCING FLOOR STICKERS

ROUND STICKER, DIAMETER 250 MM

Round social distancing sticker - SIMKALTOG_002

ADHESIVE STRIP

 $A dhesive \ strip \ 1000 \ x \ 80 \ mm \ - \ \text{SIMKALTOG_003}$

SOCIAL DISTANCING STICKER KIT

- 3 adhesive strips + 3 round stickers

Sticker kit - SIMKALTOG_010

ZAPP PURIFICATEUR D'AIR UV-C

Compact (can be mounted on wall or ceiling) Supplied with 8,000-hour bulb. Requires no maintenance. Active disinfection area: $15\ m^2$

ZAPP 18 - PPUR006















HYGIENE AND DISINFECTANTS



FOR HANDS AND SURFACES

HYDROALCOHOLIC GEL

Alcohol-based gel for hygienic treatment and disinfection of the hands by rubbing.

Aniosgel 85 NPC, 500 ml + pump - P71D017 Aniosgel 85 NPC, 1 I, airless bottle - P71D013 Wall-mounted holder for 1 I Aniosgel airless bottle - P71D014

FREQUENT-USE GENTLE SOAP

- Contains vitamin F and glycerine. For frequent hand washing.
- Maintains the skin's natural balance.

Frequent-use gentle soap with pump, 500 ml - P71D004 Frequent-use gentle soap. 11 bottle + pump - P71D005

MANUCLEAR SOAP Effectively combats hand-borne infection. Kind to hands, protects the skin.

Manuclear soap with pump, 500 ml - P71D036

PAPER HAND TOWELS

PAPER HAND TOWEL DISPENSER

Dispenser holds 200 x 2-ply paper hand towels.
42.5 x 29 cm.

Paper hand towel dispenser - P43E001

PACK OF PAPER HAND TOWELS

- 2-ply paper hand towels.
- Pack of 150.

Pack of paper hand towels - P43E002

INSTRUMENTS

PRE-DISINFECTION - ANIOSYME X3

- For cleaning and pre-disinfecting instruments in an ultrasonic cleaner or soaking tank.
- Bactericidal, acts on BK virus, Candida albicans, HIV.1 virus, hepatitis B and C.
- Concentrated: dilute to use (soaking time: 5 minutes). Rinse with water.
- Replace daily.

Aniosyme X3 11 with 25 ml integrated doser - P70P1005

1 | of concentrated product = 200 | of diluted product Aniosyme X3 51 with 25 ml doser - P70PI006

1 | of concentrated product = 200 | of diluted product

DISINFECTION - STERANIOS 2%

- For disinfecting instruments in a soaking tank. Ready to use.
- Bactericidal, acts on BK virus, Candida albicans, HIV.1 virus, hepatitis B and C.
- Soaking time: 20 minutes. Rinse with water.
- Use for up to 7 days

Sterianos 2%, 21 bottle - P71D001 Sterianos 2%, 51 bottle - P71D012










HYGIENE AND DISINFECTANTS



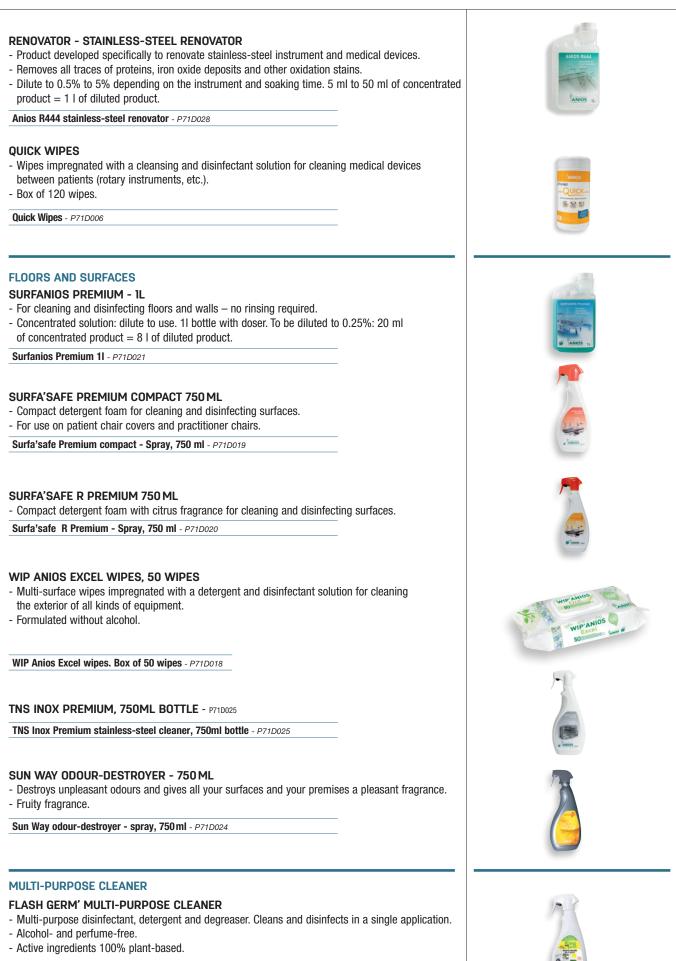

Flash Germ' Multi-purpose cleaner - Spray, 750 ml - P71D035







ORTHONYXIA - ONYCHOPLASTY

CENTRAL CLIP SYSTEM

Ovoid clips of different lengths for corrective nail treatments. - Box of 10 clips of the same size.

No.16: 17 mm - P79C018	
No.18: 19 mm - P79C019	
No.20: 21 mm - P79C020	
No.22: 23 mm - P79C021	

LATERAL CLIP SYSTEM

Clips for providing lateral traction on one side of the nail. - Box of 10 clips of the same size. No.18: 19 mm - P79C022 No.20: 21 mm - P79C023 No.22: 23 mm - P79C024

GLUE FOR CLIP SYSTEM

Glue for affixing nail clips. - 5g bottle.

Glue for clip system - P79C025

BACQUART CLIP

Bacquart clip: 10 x 0.5 cm, th.: 0.30 mm, titanium wire 0.12 mm - *P79C026* **Bacquart clip: 10 x 0.5 cm, th.: 0.30 mm,** titanium wire 0.14 mm - *P79C029*

SEMI-HARD WIRE

Stainless-steel wire, 6/10ths, for making staples.

- Cut to the required length.
- 10m roll.

Semi-hard wire, 6/10ths - P79C014

TWEED PLIERS

Stainless-steel tweed pliers, 14 cm, for making u-shaped and omega-shaped clips. **Tweed pliers, 14 cm** - *P79/001*

WIRE-BENDING PLIERS

Stainless-steel pliers for making lateral retention hooks.

Wire-bending pliers - P79/002

MATHIEU NEEDLE HOLDER

Stainless-steel Mathieu needle-holder, 14 cm, for rounding off ends.

Mathieu needle-holder, 14 cm - P791003

CEMENT SPATULA

Flexible stainless-steel cement spatula, 16 cm, for applying retention cement and composite cement.
Cement spatula, XL, extra-flexible - P79/004

GLASS

Apen glass

Apen glass - P40G001



188 www.podiatech.com

ORTHONYXIA - ONYCHOPLASTY

PHOTOPOLYMERISATION LAMP - POD ULTRA LIGHT®

Wireless hands-free lamp for photopolymerising nail resins for fast, long-lasting affixing of devices used in orthonyxia and onychoplasty.

- High-density LED, ultra-lightweight: 900 g.
- Pedal-operated automatic timer (10, 25 or 60 sec.).
- Integrated battery and polycarbonate UV filter.
- Supplied in box for home treatments.
- Box dimensions: width 11 cm depth 20 cm height 4.5 cm.
- Made in France.

Pod Ultra Light - PECI007

ORTHONYXIA STARTER KIT

Comprising: 1 PODOPRIMER bonding agent + 1 PODOCOMP D hard resin + 1 bonding agent spatula with pack of 50 brushes + 1 titanium wire in 3 diameters + 1 instrument A + 1 storage box.

Orthonyxia starter kit - P79PACKORTHONYXIE

ORTHONYXIA EXPERT KIT

 $\begin{array}{l} \mbox{Comprising: 1 PODOPRIMER bonding agent + PODOCOMP D + PODOCOMP F + PODOCOMP L resins} \\ + 1 brush handle + bags of brushes + tunnellers + 1 pair Adson tweezers + 1 titanium wire in 3 \\ diameters + 1 instrument A + 1 instrument B + 1 resin toenail + 1 storage box. \end{array}$

Orthonyxia expert kit - P79PACKONEXPERT

PODOCOMP DRY

Degreaser. For use on the nail in orthonyxia and onychoplasty treatments.

Podocomp Dry - P79C028

PODOPRIMER (BONDING AGENT)

Bonding agent with non-allergenic interface essential for resin applications. Photopolymerisable liquid composite bonding agent for podiatric use, essential for application of all resins to the nail (PODOCOMP D – PODOCOMP F – PODOCOMP L)

Podoprimer, 5 ml - P79C001

PODOCOMP

PODOCOMP D: transparent (hard) resin paste for pins and onychoplasty treatments.

Photopolymerisable composite paste specially developed for all types of orthonyxia treatments. Transparent paste that enables the practitioner to see the titanium in the sliding pin tunnel. In onychoplasty treatments, it provides a view of the nail bed.

PODOCOMP F: transparent (soft) resin paste for packing and onychoplasty treatments (nail reconstruction).

Photopolymerisable soft transparent composite resin specially designed for slightly flexible onychoplasty treatments that are resistant to microtraumas (sports and others) and neutral and impermeable peri- and subungual packing and drainage.

PODOCOMP L: transparent (hard) liquid resin for pins and onychoplasty treatments.

Podocomp D 5 g - P79C002
Podocomp F 5 g - P79C003
Podocomp L 2 ml - P79C004

TITANIUM WIRE

- Titanium wire, 0.12 / 0.14 / 0.16, for use in orthonyxia treatments.
- Cut to the required length.
- 20cm length in bag.

Titanium wire, 0,12 - P79C007
Titanium wire, 0,14 - P79C008
Titonium wire 0.16 process





ORTHONYXIA - ONYCHOPLASTY

BONDING SPATULA / BRUSH HANDLE

Sterilisable brush handle for applying bonding agent. Used for applying bonding agent (PODOPRIMER) to the nail. This instrument comprises a very strong high-quality plastic handle that is sterilisable. - Bag of single-use brushes, sold separately.

Bonding agent spatula - P79C005

Bag of 50 brushes - P79C006

TUNNELLER

System for guickly and easily creating the tunnel for the sliding pin. - Sold in packs of 10 units.

Tunneller (bag of 10) - P79C010

INSTRUMENT REF A

Sterilisable stainless-steel instrument for holding the titanium wire on to the nail. The spatula end is used to take resin from the syringe and place it on the nail. The angled notched end is used to hold the titanium wire in place during photopolymerisation.

Instrument ref. A - P79C011

INSTRUMENT REF B

Sterilisable stainless-steel instrument for holding the pin or clips in place.

The spatula end is used to take resin from the syringe and place it on the nail. The end with the hollowed-out area is used to hold the previously made pin firmly in place during photopolymerisation.

Instrument ref. B - P79C012

DEMONSTRATION TOENAIL

Resin toe for demonstrating devices and for training.

Demonstration toenail - P79C013

ADSON TWEEZERS

- Available in stainless steel, 12 cm, grooved jaws, for use with titanium wire. - Sold singly.

Adson tweezers - P61PR004

TITANIUM WIRE CUTTERS

Titanium wire cutters - P79/005

Wire cutters with oblique blades, 14 cm.




BUR KIT FOR ORTHONYXIA TREATMENTS

Plastic bur stand with cover for 12 burs (not autoclave-compatible)

- + standard diamond bur, round cone, standard grit 850.104 Ø 0.14
- + standard diamond bur, round cone, standard grit 856.104 Ø 0.16
- + special tungsten carbide bur, 5610 Ø 0.45 (blue ring)
- + special buffing bur, Ø 055 0669HP, grey colour, grinding for onychoplasty treatments
- + special buffing bur, EDENTA Ø 055 0659HP, yellow colour, fine smoothing for onychoplasty treatments
- + special buffing bur, EDENTA Ø RM 66 1025 HP, pink colour, for smoothing pins in orthonyxia treatments.

Bur kit for orthonyxia treatments - P63KITORTHONYXIE

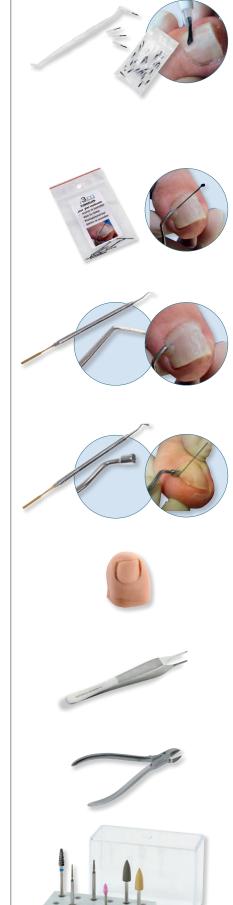

(photograph not legally binding)













NIGHT SPLINTS FOR HALLUX VALGUS

NIGHT SPLINT KITS FOR HALLUX VALGUS

Kits are available for every foot type and size (from 35-36 to 45-46). High-quality pre-cut materials. Products and protocols developed with several state-qualified podiatrist consultants. Training available on request, depending on dates and availability (contact us).

Night splint kit for 1 foot, right or left - P45CCONTNOCT

PROCEDURE (Contact us for details of the splinting procedure)

Watch the splinting protocol video.

PODIATECH

HALLUX VALGUS CORRECTOR

Foam-covered progressive corrector. Particularly suitable for post-operative treatments. 1 pair: 1 right foot + 1 left foot, one size.

Hallux valgus corrector (one pair) - P45C001

DOKI TAPE - ROLLS

Specialised application of DOKI TAPE can help tissues adapt and provide a response to various changes:

1. Decompression.


- 2. Relaxation of tissue next to pain receptors.
- 3. Decongestion of the underlying spaces.
- 4. Relaxation and normalisation of tissue.
- 5. Muscle toning and dynamisation.

A permanent decompressant

In contrast to conventional correction techniques that compress the tissues, this tape uses traction to expand previously infiltrated subcutaneous spacing reduced by infiltrates and the tension caused by dysfunctions.

By relaxing the tissue and forming cutaneous and subcutaneous microdepressions, it promotes a return to normal fluid exchanges.

Roll, 5 m x 5 cm wide Black - Blue - Pink P45P020 21 22

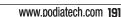











DOKI



ORTHOPLASTY

The leading products for use in all types of orthoplasty treatments. Suitable for use with infants and at-risk patients. • PROTECTIVE • CORRECTIVE • THIN • HARDWEARING • COMFORTABLE

J3M SILIPOD - 3 hardnesses: Shores 5 - 20 - 35, beige. SILIPOD; mix together the two parts in equal quantities: 13m - 50% Base paste (B) 5 (extra-soft) or 20 (soft) or 35 (firm) - 50% Catalysing paste (C). The Shore 3 hardnesses of the SILIPOD range can be used together to make inclusions, excess load protection and repairs. CINIER B PODORSIL - 3 hardnesses: Shores 20 - 30 - 45, beige. PODORSIL 3P; mix together the two parts in equal quantities: - 50% paste: B20 (soft) or B45 (firm) or B20 + B45 (medium) CINIER ****** B - 50% Catalysing paste (A). PODORSIL 30; mix together the two parts in equal quantities: - 50% B30 paste (medium). - 50% Catalysing paste (A). **PODORSIL 3P** 3P A catalyst, 500 g - PO8PODOR3PA-50 3P A catalyst, 250 g - P08PODOR3PA-25 3P A catalyst, 125 g - P08PODOR3PA-12 3P B20 soft part, 250 g - PO8PODOR3PB2-25 3P B20 soft part, 125 g - P08PODOR3PB2-12 3P B45 firm part, 250 g - PO8PODOR3PB4-25 3P B45 firm part, 125 g - P08PODOR3PB4-12 HARDNESS 30 developed **EXCLUSIVEL PODORSIL 30** for PODIATECH (MEDIUM HARDNESS, SOLD AS A KIT, A+B) Kit: Podorsil 30 A+B, 2 x 300g pots - PO8PODOR30AB SILIPOD 5 PROTECSIL Extra-soft viscosity, Shore 5. Suitable for orthoplastic devices designed for protection and comfort and recommended for at-risk patients (diabetes, arthritis, etc.) and inclusions. Silipod 5 Protecsil - 2 x 250 g - PO8SILIPOD5 SILIPOD 20 Soft viscosity, Shore 20. Extremely useful for protective and comfort orthoplastic devices and minor correction. Silipod 20 - 2 x 250 g - PO8SILIPOD20 SILIPOD 35 Firm viscosity, Shore 35. For orthoplastic devices designed for correction and light protection. Silipod 35 - 2 x 250 g - PO8SILIPOD35

ULTRA-FINE MESH ORTHOPLASTY INJECTION AND REPAIRS

SILITRAM

Nonwoven synthetic fabric mesh with low elasticity for ultra-fine and very resistant mesh orthoplasty treatments.

- Length 1 m.
- Width 8 cm.

SILITRAM - PO8SILIPODTRAM

SILIROLL

Ergonomic roll designed especially for mesh orthoplasty treatments. - Width 10 cm.

SILIROLL - PO8SILIPODROL

SILIPOD REP - DOUBLE-CARTRIDGE SYRINGE, 40 G (A+B) FOR REPAIRING ORTHOPLASTY TREATMENTS

Fluid silicone that acts like a glue, very useful for repairs, adding material to and finishing all types of orthoplasty treatments.

Grinding is possible after catalysing (reticulation) for 5 min.

Can be mixed with and sticks to all the other silicones in the range as well as other brands.

SILIPOD REP - PO8SILIPODREP

SILINJECT

Syringe for injecting all types of silicone in injected orthoplasty treatments under load on a podoscope (cold sterilisable).

SILINJECT - PO8SILIPODINJ

PACK OF ORTHOPLASTY BURS

Plastic bur stand with cover for 12 burs (not autoclave-compatible)

- + special tungsten carbide bur, 5510 Ø 0.60 For grinding hollowed-out area
- + special buffing bur, EDENTA Ø 055.0666HP for finishing buffing of silicones
- + special bur in the shape of an abrasive disc for finishing grinding: fine
- + support bur Diameter 0.5 for finishing grinding: rough work
- + support bur Diameter 0.7 for finishing grinding: rough work
- + support bur Diameter 10 for finishing grinding: rough work
- + abrasive tip Diameter 0.5
- + abrasive tip Diameter 0.7
- + abrasive tip Diameter 10.

Pack of orthoplasty burs - P63KITORTHOPLASTIE

CAF ADHESIVE

CAF adhesive, 100 g - PO8COLLECAF

ORTHOPLASTY BOXES

Pack of 10 green boxes - P40B006
Pack of 10 white boxes - P40B004
Pack of 10 blue boxes - P40B005







(photograph not legally binding)









WARTS

PODIATECH, EXCLUSIVE DISTRIBUTOR OF CRYOPEN® TO PODIATRISTS. Cryogenic treatment solution for verrucas and granulations. - Precise. - Effective.



- Fast.
- Pain-free.

Ask your technical sales advisor for advice.

Characteristics:

- -89°C
- Depth: 6 mm.
- Pressure: 55 bar.
- Assured elimination of necrotised cell: -27°C + fast descent.



ryoPen

Before

CRYOPEN B

+ 3 applicators (1 GOLD pinpoint applicator Ø 1 to 4 mm, 1 SILVER applicator 2 to 8 mm, 1 BRONZE applicator Ø 4 to 12 mm).

- For use with 23.5g cartridges.

CryoPen/B + 3 applicators - P77C001

APPLICATORS

CryoPen/B GOLD micro-applicator - 1 to 4 mm - P77C005 CryoPen/B SILVER micro-applicator - 2 to 8mm - P77C006 CryoPen/B BRONZE micro-applicator - 4 to 12 m - P77C007

CARTRIDGES

- 23.5g cartridges: 300 seconds of constant cryogenic application. Pack of 6.
- 16g cartridges: 200 seconds of constant cryogenic application. Pack of 6.
- 8g cartridges: 100 seconds of constant cryogenic application. Pack of 24.

CryoPen 23.5g cartridges (CryoPen B). Pack of 6 P77C002 CryoPen 16g cartridges (CryoPen M, O, X+). Pack of 6 P77C003

CryoPen 8g cartridges (CryoPen 0, X+). Pack of 24 P77C004



ONYCHOMYCOSIS - HYPERHIDROSIS

LUNULA LASER

A unique way to treat onychomycosis

An effective and pain-free non-thermal laser treatment for fungal nail infections. Treats the fungi responsible for nail detachment and discolouration. A small number of 12-minute sessions stimulates anti-fungal responses.

Lunula laser - P76L001

Contact us for marketing tools:

- "Waiting room" communication pack to inform your patients
- (video poster 100 flyers).



-



A MULTIFACETED TREATMENT FOR FUNGAL NAIL INFECTIONS

IN YOUR PODIATRY PRACTICE



FREE INSTALLATION AND FIRST USE INTRODUCTION

FREE "WAITING ROOM" COMMUNICATION PACK

IONOMAT®

Device for treating hyperhidrosis and plantar warts.

Set includes:

- 1 control unit (IONO cycle / WART cycle / AC current intensity adjustment).
- 2 electrodes + 2 podiatry-sized trays.
- 2 grids + 2 wires with crocodile clips.

Session duration: 20 minutes.

Treatment duration: 4 or 6 weeks or 8 to 10 sessions.

IONOMAT - PBAC007

100g bag of sodium salicylate - PBAC009







• 1-year warranty

QENOA[®] PRODUCTS







COMFORT CREAM 500 ml FOOT CARE LOTION

Natural post-treatment lotion made in France. *Rehydrates and intensely nourishes the skin.*



Shea butter: soothes and calms irritated skin.



Avocado oil: rich in vitamin E, nourishes deep down.



Organic witch hazel floral water: *revitalises and acts on skin scarring.*

Organic silicon: limits skin cell ageing.

Qenoa Comfort cream, 500 ml - P78POD001







Sudalun®, 180 ml - P78T003

Contact us at contact@qenoa.fr for details of our retail prices and to request samples. Essential oils have properties that are not recommended in certain situations. Please ensure you know how to use them safely.

www.podiatech.com 199

MARKETING TOOLS



Photographs not legally binding

30 cm x 30 cm











PRESCRIBER'S GUIDE TO PLANTAR ORTHOTICS P80GUIDEPRES_05 Size: 15 cm x 21 cm



PRACTICAL GUIDE TO STRETCHING P80LIVREKINE 2015 edition



CAR PHONE STAND P81SUPPORTTEL

RUNNING SHOES SIMKLIVRCHAUS

BOOK: HOW TO CHOOSE THE RIGHT



WAITING ROOM POSTERS

FEET: P80POSTFR22_VPIED

PATHOLOGIES: P80POSTFR22_PATHO

PHONE & TABLET DESK STAND P81SUPPORTORSO



PODIATECH USB STICK WITH BLUE LED P81CLEUSB17

BLOCKS OF PODIATECH STICKY NOTES

P81POSTITPODIA



BOOK: THE FOOT, HOW DOES IT WORK? SIMKLIVRPIED

Qenoa



QENOA TOTE BAG

Hardwearing bag with long handles,

100% cotton. Size: 38 x 8 x 42 cm

P78MK004

PODIATECH RIBBON KEY FOB P81PORTCLEFS



PODIATECH TOTE BAG P78MK004 Hardwearing bag with long handles, 100% cotton. Size: 38 x 8 x 42 cm



PHOTOGRAPHS NOT LEGALLY BINDING 200 www.podiatech.com









AN INTERNATIONAL PRESENCE

Podiatech provides its customers with the benefits of the world's leading podiatry network. Our distribution partners and frequent presence at international trade fairs remain a core part of the company's culture and development. This is how we help you to succeed.



For more information contact your local agent or go to www.podiatech.com

CONTACT HEAD OFFICE SIDAS SAS PA Champfeuillet Est 18 rue Láon Béridot - CS 7

18, rue Léon Béridot - CS 70353 38516 VOIRON - CEDEX France

International Sales Dpt. Ph.: +33 4 76 67 07 07 Fax: +33 4 76 67 35 12 exportmedical@sidas.com

CONTACT EXPORT AREA MANAGER DIRECT SALES Richard SERRE serre@sidas.com Ph.: +33 7 86 08 76 36



Always striving to perfect the quality of its products, PODIATECH reserves the right to make changes to the catalogue descriptions without notice. REF. :P80CATPODIAEXP22 Photograph credits: SIDAS SAS - No Comment - Photographs not legally binding.

SIDAS SAS - Capital €600,000 - APE code 3250 A - CORP. & TRADE REG. No. 304 577 661 00073 - Grenoble - INTERCOM. VAT No. FR 76 304 577 661